

SAMPLE COLLECTING DEVICE AND MASS
SPECTROMETRY OF DEVICE

Technical Field

The present invention is concerned with methods and devices for sample collection and simultaneous detection and/or quantitation of multiple trace elements in fluid samples.

Background Art

A wide range of trace metals and other elements is necessary for good health and physical well being in humans and other animals; deficiencies in essential elements have been shown to cause general malaise and lead to the induction of specific disease, commonly resulting in death. For many essential trace elements, it is not simply the absolute concentration, but also the inter-element balances that have a profound effect on health. For example, selenium deficiency is implicated in the aetiology of Iodine Deficiency Disorders amongst humans, whilst copper deficiency, associated with high levels of manganese, may be implicated as a predisposing or causative factor in induction of Bovine Spongiform Encephalopathy (BSE) in cattle and, by association, New Variant Creutzfeldt-Jakob Disease (nvCJD) in humans.

Dietary forages, vegetables, grains and fruits, which fix available trace elements as metal colloids within their tissue, have long been regarded as sources of essential trace elements. Such plant-based metal colloids are about ninety-eight percent absorbed and communities and animals that have a balanced range of plant products as essential components of diet may reasonably be expected to display markedly reduced incidence of specific trace element deficiency-related disease when compared with other groups lacking quality forage or a regular vegetable, fruit and grain intake.

The trace element content of vegetative material is directly related to the bioavailability of essential nutrients in soils supporting the vegetation. Soils vary in their trace element content from enriched to impoverished, according to local geology, soil degradation and nutrient impoverishment and as a function of inappropriate cropping practice, which is widespread throughout the world. In addition, soils throughout the world are sustaining increasing anthropogenic chemical damage threatening the existence of many plants and animals. Consequently, human health is being threatened through the food chain.

While the productivity of the soils may be maintained through the application of N-P-K fertilisers, food crops growing on these soils becomes, without the regular application of biologically-available 'balanced' trace elements, progressively impoverished in essential trace elements and minerals. If not corrected, this may result in sharply increased incidences of mineral deficiency-related disease.

BEST AVAILABLE COPY

Elements may be classified as being essential or toxic to human and animal health. In the case of animals, trace metal deficiency and/or toxicity is due largely to concentration levels controlled by environmental factors, whereas for humans, both environmental and occupational factors may be important; toxic response may a function of both natural and/or anthropogenic influences.

Ignoring carbon, hydrogen and oxygen, the biologically essential major elements are calcium, chlorine, magnesium, phosphorous, potassium, sodium, nitrogen and sulphur. Essential trace elements include bromine, chromium, cobalt, copper, fluorine, iodine, iron, manganese, molybdenum, selenium, silicon and zinc. If bio-available, many of these essential trace elements induce toxic responses, at elevated levels, or if out of balance with synergistic and/or antagonistic elements. Several other elements (lithium, scandium, rubidium, lanthanum) are minor essential elements.

In addition to dietary trace metal deficiency-induced disease, other cohorts of individuals are occupationally or environmentally exposed to a range of toxic element pollutants, which similarly induce general malaise and/or specific clinical symptoms commonly resulting in complications and death. Notable amongst these are arsenic, lead and mercury, which constitute the top three most hazardous substances on the US Environmental Protection Agency's Toxic Substances and Disease Registry priority list.

The leaching of heavy metals into the aquatic environment, and uptake by wildlife in the food chain, may have a profound impact on human health. Cadmium and mercury, in particular, are strongly bio-accumulated in fish and shellfish.

Although it is not possible to quantify the hazards and deleterious effects associated with all trace elements, some elements clearly present a more serious problem than others. Respectively ranked 1, 2, 3 and 7 on the NPL, arsenic, lead, mercury and cadmium, as elemental pollutants, are considered extremely toxic and the health effects of these elements have received a great deal of attention from research workers. Other elements on the list, in alphabetical order, are aluminium, antimony, barium, beryllium, chromium, cobalt, copper, manganese, nickel, plutonium, radium, selenium, silver, thallium, thorium, tin, uranium, vanadium and zinc.

Unlike many essential trace elements, the concept of a therapeutic index cannot be applied to toxic elements such as lead, cadmium, mercury and arsenic. These toxic elements play no known role in metabolism, as no enzyme has been identified which specifically requires any of them as cofactors. They are extremely hazardous to life and, resulting from ingestion, have been involved in historic poisoning episodes of both human and animal populations. They are increasing in concentration in both aquatic

and terrestrial environments due to anthropogenic inputs, and thus will continue to be a concern to toxicologists and clinicians.

Hence, proactive intervention to identify trace metal and element aberrations within general populations, thereby enabling the early implementation of targeted remedial strategies with consequent minimization of the huge social impact of trace metal-induced disease, is essential. However, mass screening of general populations for trace metal deficiencies and/or toxic metal excesses, with reference to age, sex, socio-economic status and physical geography, while acknowledged as being highly desirable in terms of preventative medicine, is presently impractical. So too, is the mass screening of human food chain components, such as slaughter animals, prior to their entering the food chain.

Present test methodologies require relatively large volumes of fluid samples (for example, 5-10 ml of blood) and are commonly trace element specific, that is, simultaneous measurement of other trace elements potentially present is not possible. Because of this, other relevant trace metals are either overlooked or require further fluid samples for their determination. In the case of blood, this involves invasive, often traumatic extraction, particularly for young children, babies and the elderly, using hypodermic syringes. The derivative body fluid products require stabilisation and preservation, and having regard for transmissible disease such as HIV, appropriate biohazard handling and disposal. Further, the large volumes required give rise to handling and storage problems.

There is no current technology available that can conveniently be used for the collection and broad-spectrum analysis of the trace element content of large numbers of blood and other body fluid samples. Presently available testing methods are cumbersome and expensive, placing the service outside the reach of the general population, particularly in underdeveloped regions where problems are often greatest. Further, there are no convenient and sensitive mass spectrometric methods for detecting pollutants or contaminants in fluids such as water or lubricants.

There is therefore a need for improved methodologies which will enable more efficient and cost effective screening of trace elements in fluid samples.

It is an object of the present invention to alleviate at least some of the disadvantages of prior art methods, or to provide a useful alternative.

Summary of the Invention

According to a first aspect there is provided a sample collection device comprising an inert collection matrix capable of adsorbing or absorbing a fluid sample, and a solid support, wherein the inert matrix is affixed to an area of the solid support

Particularly useful matrices may be selected from aragonite, aluminium hydroxide, titania, glucose, Starch "A", Starch "B", glucodin, cellulose powder/granules, fibrous cellulose, hydroxy butyl methyl cellulose, vegetable flour and the like, or mixtures thereof. Particularly preferred is fibrous cellulose. The fibrous cellulose matrix
5 may be modified by oxidation and/or acid hydrolysis to improve its properties and thus provide enhanced reproducibility and sensitivity.

The vegetable flour may be selected from rice, maize, wheat, soy, rye or corn flour, or mixtures thereof. Particularly preferred is rice flour.

The inert matrix may also contain, on or within, one or more pre-calibrated
10 selected analytes as internal standard, to aid in the quantitation of trace elements in the sample applied to the collection device.

The device of the present invention may also comprise an integral lancing member, capable of piercing for example skin or tissue, to aid in the collection and application of a blood or body fluid sample to the inert matrix. The lancing member may
15 be mounted adjacent to, within or below the area of inert matrix. There may be included a guiding channel in the inert matrix, to guide the lance should it be disposed below the inert matrix area.

The device may also be equipped with a laser-scannable bar code which may contain patient information or other information concerning the sample, its nature and
20 source. The device may also include an antibiotic barrier, to prevent contamination of the sample to analytical equipment and personnel.

Preferably the inert matrix is applied to only one side of the support. It is also preferred that the area to which the matrix is applied is smaller than the area of the solid support and that it be in the shape of a small tablet-sized disc.

25 The inert matrix may include hydrophobic and/or hydrophilic components, depending on the nature of the sample and the analysis to be performed.

Preferably the solid support is made of flexible material having sufficient durability to withstand transport and handling. Of course it will be understood that the support can be made of rigid material, depending on the nature of application. It is also
30 preferred that the device is of sufficiently small size to allow transport of the device through mail and for ease of storage. The device may have an integral or separate cover sheath, to protect the inert matrix and prevent possible contamination after collection. The cover sheath also protects the device during transport and handling.

According to a second aspect there is provided a sample collection device having
35 multi-layer construction wherein the collection matrix layer is sandwiched between two

supporting layers, one of said supporting layers having an opening, which exposes an area of the collection matrix.

Alternatively, the sample collection device may encapsulate a collection matrix tablet within the body of the support wherein the matrix is exposed flush with one surface of the support.

The collection device and methods of the present invention may be used for analysis of any fluid sample, including body fluids, oils and other lubricants, water from drinking supplies as well as waste water, and the like. Body fluids such as whole blood are particularly preferred, however, separated blood (eg. plasma or serum) and other body fluids, such as urine or sweat, can also be used with the same device.

It will be understood that a sample of body fluid, particularly blood, can be collected for analysis by conventional means, or by using for example a sample collection kit comprising a resealable, sterile sample collection device, embodying a bar coded support in which is embedded, or to which is affixed, a tablet, wafer, wad, strip or the like, of sample absorption/adsorption matrix, a sealed alcohol-saturated wipe, and a separate retractable, single use, spring-loaded lance for penetrating the skin and drawing blood. Of course a lance can be omitted from the kit if the sample to be collected is for example urine or sweat.

As indicated above, the analytical sample need not be a body fluid. Thus, the devices and methods of the present invention are equally applicable to collection and analysis of water or oil samples without significant adaptation of collection devices or analytical procedures and equipment.

The matrix of the sample collection device can include one or more matrix-matched standards either adsorbed/absorbed onto/into sample collection matrix or, alternatively, supported on an impermeable substrate. Here, the matrix may be spiked with elements, for example, Be, In and Hf and these elements will serve as internal standards that will be released simultaneously with the sample during ablation; this will facilitate matrix matching.

According to a third aspect there is provided a method of detecting simultaneously a plurality of elements in a fluid sample adsorbed onto or into an inert collection matrix, comprising:

(i) exposing the sample to high energy radiation capable of ionising at least a portion of the sample, and

(ii) detecting plurality of elements in the ionised portion of the sample by mass spectrometry.

According to a fourth aspect there is provided a method of quantifying simultaneously a plurality of elements in a fluid sample adsorbed onto or into an inert collection matrix, comprising:

- 5 (i) exposing the sample to high energy radiation capable of ionising at least a portion of the sample;
- (ii) measuring quantity of a plurality of elements in the ionised portion of the sample by mass spectrometry;
- (iii) measuring quantity of ionised portion of sample, and
- (iv) determining quantity of the plurality of elements in the sample.

10 According to a fifth aspect there is provided a method of quantifying simultaneously a plurality of elements in a fluid sample adsorbed onto or into an inert collection matrix having an internal standard applied thereto, comprising:

- (i) exposing the sample to high energy radiation capable of ionising at least a portion of the sample and a portion of said internal standard;
- 15 (ii) measuring quantity of a plurality of elements in the ionised portion of the sample by mass spectrometry;
- (iii) measuring quantity of ionised internal standard in the ionised portion of the sample by mass spectrometry, and
- (iv) determining quantity of the plurality of elements in the sample with reference
- 20 to quantity of ionised internal standard.

According to a sixth aspect there is provided a method of quantifying simultaneously a plurality of elements in a fluid sample adsorbed onto an inert collection matrix, comprising:

- 25 (i) introducing into the fluid sample a known quantity of a measurable internal standard
- (ii) exposing the sample to high energy radiation capable of ionising at least a portion of the sample and the internal standard;
- (iii) measuring quantity of a plurality of elements in the ionised portion of the sample by mass spectrometry;
- 30 (iv) measuring quantity of ionised internal standard in the ionised portion of the sample by mass spectrometry, and
- (v) determining quantity of the plurality of elements in the sample with reference to quantity of ionised internal standard.

According to a seventh aspect there is provided a method of quantifying
35 simultaneously a plurality of elements in a fluid sample adsorbed/absorbed onto or into an inert collection matrix comprising:

(i) exposing the sample to high energy radiation capable of ionising at least a portion of the sample;

(ii) measuring quantity of a plurality of elements in the ionised portion of the sample by mass spectrometry;

5 (iii) exposing a matrix-matched Certified Reference Material (CRM) to high energy radiation capable of ionising at least a portion of the CRM;

(iv) measuring quantity of ionised CRM in the ionised portion of the sample by mass spectrometry, and

10 (v) determining quantity of the plurality of elements in the sample with reference to the CRM.

According to an eighth aspect there is provided a method of quantifying simultaneously a plurality of elements in a fluid sample supported on an impermeable substrate, comprising:

15 (i) exposing the sample to high energy radiation capable of ionising at least a portion of the sample;

(ii) measuring quantity of a plurality of elements in the ionised portion of the sample by mass spectrometry;

(iii) exposing a matrix-matched Certified Reference Material (CRM) to high energy radiation capable of ionising at least a portion of the CRM;

20 (iv) measuring quantity of ionised CRM in the ionised portion of the sample by mass spectrometry, and

(v) determining quantity of the plurality of elements in the sample with reference to the CRM.

25 Details of some useful CRM's, for example, SARM 1, 3 and 46 (South African Bureau of Standards), and SY-2 (Canadian Certified Reference Material Project (CCRMP)) are given in Table 1. Other standard element cocktails may include elements such as Be, In, Hf, Bi, Th to cover the mass calibration range, but may include any element as a standard, that is not being analysed.

30 Preferably, the sample is whole blood and sample size is approximately 50 μ l to 100 μ l and even more preferred size of sample is 50 μ l or less. Of course, separated blood may also be used, eg. plasma or serum.

Also preferred is that the high energy radiation is UV laser radiation and that the sample is exposed to such radiation for a period of approximately 30 seconds, , but may be between 10 and 120 seconds.. The devices and methods of the present invention
35 may be used in conjunction with any Inductively Coupled Plasma-Mass Spectrometer

(ICP-MS) system. Particularly preferred are quadrupole and Time-of-Flight (TOF) ICP-MS systems.

The preferred elements to be detected and/or quantified are dietary trace elements, toxic elements and markers of pollution or wear and tear. For blood and other body fluids, these elements can include Li, Na, Mg, Al, P, K, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Rb, Sr, Mo, Cd, Sn, Sb, Te, Ba, La, Ce, Eu, Dy, Yb, Hg, Tl, Pb, Th and U. For wear metals in lubricants such as oil, the element array may include Li, B, Mg, Al, Si, P, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Sr, Y, Zr, Mo, Ag, Cd, Sn, Sb, Ba, La, Ce, Hf, Hg, Pb, and U.

In a preferred embodiment the matrix or the support comprise one or more wells or indentations to accommodate the fluid sample.

According to a ninth aspect there is provided a method of collecting a fluid sample for mass spectrometry analysis of multiple element content comprising the application of the sample to an inert matrix having a low background element content, wherein the matrix is selected from the group consisting of aragonite, aluminium hydroxide, titania, glucose, Starch "A", Starch "B", glucodln, cellulose powder/granules, fibrous cellulose, hydroxy butyl methyl cellulose, vegetable flour or mixtures thereof.

Description of the Preferred Embodiment

The present invention is in part based on Laser Ablation-Inductively Coupled Plasma-Mass Spectrometry technique, which allows rapid, automated, cost effective mass screening of general populations, bloodstock, zoo animals, pets and slaughter animals to identify trace element aberrations in body fluids. This technology facilitates proactive remedial intervention to target and correct essential trace element imbalances and/or toxic heavy metal excesses and enables identification and rejection of heavy metal-contaminated slaughter animals designed for human consumption. The methods and devices of the present invention are also useful for detection and quantitation of trace elements, metals and the like in fluids such water and lubricants, as indicators of for example water pollution or mechanical wear and tear.

The present invention in its various embodiments allows the simultaneous analysis and/or quantitation of a broad spectrum of up to 50 trace elements during a primary analytical run. A secondary run, using a screened torch may include Ca, Mg, Na, K and Fe. The analytical cost of a sample is lower than that of a large number of single element analyses currently being performed, on a chemically unmodified 50-100 micro-litre volume of body fluid sample or other fluid sample (single drop) adsorbed onto an inert collection matrix. In case of blood, the sample collection device, and collection protocol, may be so configured to eliminate the use of hypodermic syringes, and hence

potential for stick injuries, is non-invasive and hence, non-traumatic, and does not involve the preservation, movement and storage of large volumes of blood and urine, or involve large biohazard disposal facilities. Indeed, in the case of humans, samples may generally be self-acquired at any geographic location through absorption/adsorption of a drop of biological fluid, such as blood from a pin prick, into/onto a lightweight collection device as described herein, and dispatched to the nearest analytical facility by post or courier. Because an approximately 8000°C argon plasma is involved in ionisation of the samples, the body fluid samples are expected to be largely sterilized during analysis.

Certain embodiments of the present invention have been developed using an ultraviolet laser and quadrupole inductively coupled plasma-mass spectrometer (LA-ICP-MS) with manual sample handling. However, the present methods are equally applicable to Time-of-Flight (ToF) and High Resolution mass spectrometry techniques. Further, the methods of the present invention, whether they make use of quadrupole, ToF or High Resolution mass spectrometry, can be automated to allow rapid, high volume throughput screening of samples.

The methods and devices of the present invention permit cost effective, simultaneous, automated mass screening of blood, and other body fluids, for a wide range of essential and toxic trace elements on micro-litre volumes of test fluid absorbed onto inert collection matrices. In certain preferred embodiments the core of the analytical system comprises a quadrupole Laser Ablation-Inductively Coupled Plasma-Mass Spectrometer. The spectrometer may be used in conjunction with an associated automated sample insertion system.

In preferred embodiments of the present invention the collection device, or kit of parts, is envisaged to consist of the following components:

- housing mount that forms the surround of the actual collection matrix and acts as the support of this matrix and also increases robustness of the entire device allowing for transport of the entire system;
- the collection matrix itself consisting of an absorptive pellet;
- a mechanism for puncturing skin and facilitating the collection of a single drop of blood; and
- a bar code or equivalent which ultimately will facilitate the recognition of both the sample and its association with the client.

However, the collection device, or kits of parts, may exclude certain features or include additional features.

The invention will now be described in more detail with reference to non-limiting examples.

Examples

Example 1: Sample collection and application

5 Samples may be collected and applied to a chosen collection matrix of the present invention in a conventional manner well known in the art.

For example, blood from a subject may be collected using a kit which comprises a shielded, retractable, spring loaded 'pricker', as part of the sample kit, which also includes a sealed, alcohol-saturated wipe, or swab, for pre-cleaning the skin area to be
10 pricked to avoid unnecessary sample contamination.

It will be understood however that collection of samples of other body fluids, such as urine and sweat, or other fluids such as water or oil and other lubricants, will not require most of the components stipulated above for blood collection, but it will nevertheless be important to exclude contaminants. Conventional techniques for this
15 will be known to those skilled in the art.

The fluid sample, whichever fluid may be of interest, can be applied to the collection matrix for analysis by any known means. For example, a particular quantity may be applied to the collection matrix by a pipette, a capillary tube, a dip-stick or similar device. Exact quantity applied is not important but may be controlled if desired.

20 Alternatively, particularly for blood sample collection, a collection device such as described in Example 2 below may be used.

Example 2: Sample Collection Device

An example of one type of sample collection device of the present invention, particularly suitable for collection of a blood sample, incorporates an inert fluid
25 absorption matrix, most preferably a fibrous cellulose matrix (Whatman 540, but also 541, 542 and other cellulose filter papers, Whatman International Ltd, Maldstone, England), typically shaped in the form of a small tablet-size disc. The matrix is affixed to or encased within a small, lightweight, disposable or re-cyclable holder (disc holder or solid support material). Ideally the holder is made of relatively rigid material (for example
30 plastic, cardboard or similar material). The device is designed so that a drop of blood or body fluid can be placed on the absorption matrix and the device sealed at the site of collection. Thus immobilized sample can be easily transported via post or courier to a sample analysis center and/or stored.

Of course the device may be used for other samples, which are not body fluids.
35 For example water or a lubricants.

A collection device of this embodiment of the present invention, incorporating a number of features described below, is depicted in Figure 1. In plan view (A) the device is typically rectangular in shape and has an area of absorbent collection matrix (1) disposed on the surface, and may also have a bar code (2) containing relevant information about the sample and/or the subject. The collection matrix is preferably fibrous cellulose but other matrices described hereafter may also be used. The collection area shown is circular in shape but may be any other suitable shape. A cover sheath (B) may be provided, to cover the collecting matrix area after the sample has been collected. Figures 2 and 3 show the collection device in cross section, in closed and open positions respectively. The carrier or backing (support) portion (A) of the device can be suitably made of plastic or some form of card (stiff paper, cardboard and the like) material. The cover sheath (B) may be made of similar materials. Both the backing portion and the cover sheath may include a locking ridge (3), for positive engagement between the backing and cover sheath, and also to prevent the cover sheath, if used, from sliding off entirely.

Figures 2 and 3 also show the area of collection matrix (1) and a stylus or lance (5) disposed below the collection matrix and within the carrier or backing material. The lance may be guided by a channel (4) in the collection matrix, so that when the device is pressed between the thumb and a finger, the lance will be forced through the channel and into the finger, thus piercing the finger and enabling a sample of blood to be collected onto the collecting matrix. Once the sample has been taken, the cover or sheath can be slid over the collecting matrix, thus protecting the sample as well as individuals handling the used device.

Figure 4 is an enlargement of a section of figures 2 and 3, showing in more detail the preferred arrangement of the lance, collection matrix and the guiding channel.

Typically, a collection device contemplated herein, in a particular preferred configuration, will have dimensions of approximately 40x20 mm and will be about 2 mm thick. However, larger or smaller collection devices may be useful in different applications and can be designed along equivalent parameters.

The collection device is primarily designed for the collection of blood and other body fluids prior to analysis of the trace element content. However, similar design principles can be used for sample collection of other fluids, omitting the integral lance. Of course, even for blood sample collection, the device described above may be provided with a separate lance, packaged together in a kit of separate components if desired.

The design of the sample collection device provides for low manufacturing costs, a robust configuration, ease of transportation, ease of storage, and can be used to collect a drop of test sample from a remote site by an inexperienced collector.

5 The matrix, which forms an integral part of the device, is typically an inert material with respect to fluid interaction prior to analysis and does not interfere with the subsequent sample analysis. The sample adsorbed onto or into the matrix can be stored indefinitely, without the addition of preservatives that may add contaminants to the sample.

10 The preferred material suitable for the matrix is cellulose, either granular or fibrous and may be either formed or preformed. Typically, the sample of blood transferred to the blood collection device does not have a specific volume. Hence the matrix may be encoded with an internal standard to normalize the analytical data on analysis.

15 The matrix may also be composed of inorganic materials suitable for a matrix of the ceramic-type, for example compounds of lithium, boron, carbon, magnesium, aluminium and silicon. Although this list is not exhaustive, it does encompass the main ingredients for an appropriate robust thermo-ceramic.

20 Typically, a sample of blood is transferred to the collection device that has a small lance or puncturing needle incorporated into the matrix, or into the backing/support material. The patient grips the device and causes a small pinprick to be administered. The collected blood does not have to have a specific volume as the matrix can be encoded with an internal standard, which normalizes the analytical data on analysis.

25 The device can have a laser-scannable bar code for recognition of the patient or to include any other additional information on the sample and its source. The amount of blood required is usually less than 50 μ L. The device can also have a sealing mechanism to ensure that the device plus sample can be transported and will not be contaminated.

30 The matrix may be affixed to, or encapsulated within, the support material or holder by any known means and may employ adhesives. Further, an antibiotic barrier may be applied to prevent contamination of the sample or the analytical equipment and personnel.

35 The present invention also makes use of collection devices which do not possess a collection matrix affixed thereto. The collection matrix may be simply omitted and the sample applied directly to the support material (backing). This may be particularly useful in certain body fluid collection devices. In such devices it may be advantageous to

introduce indentations (wells) into the support material, to allow for sample immobilization or the application of multiple samples and/or standards to the same support material (device) by application to multiple indentations (wells) in the support material.

5 Sample of fluids applied to any of the collection devices describe herein may be dried before analysis.

Example 3: Sample Analysis System

Traditionally, quantitation in LA-ICP-MS has been approached by controlling the power coupling the laser to the sample, to ensure uniform ablation characteristics and transfer of uniform amounts of solid to the analytical plasma. While this has much to
10 recommend it when the nature of the matrix can be assured (eg. glass or similar), there are significant problems associated with standardisation of the coupling and transfer efficiency when matrices are not uniform. Furthermore, when the surface characteristics of the sample also vary it is extremely difficult to ensure uniform ablation.

15 Until the present invention laser ablation ICP-MS technology has been at best a semi-quantitative technique and more usually a comparative technique for the determination of trace element levels in any solid material. In this embodiment of the invention quantitation in LA-ICP-MS has been approached by quantitation of the amount of debris (ablated or ionised material) that is actually transported from the laser cell to
20 the analytical plasma.

When using an Infrared laser, where the particle size of ablated material is relatively large, Ultra-violet spectral interference can be used to quantify the amount of particles (ablation efficiency) entering the plasma. However, in the majority of cases the techniques currently employ either UV or Excimer lasers. These lasers produce
25 particles that are too small to have sensible UV scattering and consequently relatively inexpensive particle quantitation is not possible. However, laser interferometry can be used, as an appropriate alternative technique, to quantitate the amount of ablated material and thus the efficiency of UV lasers. Once transport efficiency is quantified, it is then possible to quantify the amount of particles that are entering the analytical plasma and hence quantify the resulting signal (ie. amount of any one element).
30

The quantification process can be further enhanced by using internal standards in the support matrix of the collection/transportation device described above, or by adding one or more standards to the sample to be analysed. A suitable internal standard can be selected from elements which are not commonly present or are below
35 detectable levels in a particular sample. Thus, for blood samples, elements such as Hf, Ir, Ru, Rh, Ta and heavy rare earths can be used as internal standards, and

incorporated into the inert matrix by bonding to the surface of the particles used to produce the matrix, or may even be present as a natural constituent of the sample itself.

In case where the internal standard is incorporated into the matrix, when the sample is ablated, the particles of the matrix are carried into the analytical plasma along with the sample. Quantitation of the transport efficiency of all debris is achieved using
5 laser interferometry, or an appropriate alternative technique, and supported by normalisation to the signal from internal standards. Since the bonding characteristics of the internal standards and the efficiency of absorption of the matrix are known, as is the transport efficiency, it is possible to calculate the concentration of the element in the
10 sample adsorbed onto the matrix, in this case blood.

In another embodiment of the present invention, quantitation by LA-ICP-MS has been approached by quantitation against matrix-matched standards.

Quantitation is achieved by using internal standards in the collection matrix, or by adding one or more standards to the sample to be analysed. A suitable internal
15 standard can be selected from elements that are not commonly present or are below detectable levels in a particular sample. Thus, for blood samples, internal standards are incorporated into the inert matrix through solution doping, or may even be present as a natural constituent of the matrix itself. The collection matrix is doped with the relevant standards to act as mass calibration standards. These may be Be, In and Bi, or
20 other suitable combination depending upon the analysis required. In addition any other analyte can be spiked into the matrix pad and the pads analyzed. The spiking of calibration standards onto the matrix pad allows for its analysis as a "blank". To the standard-spiked matrix pads, blood, sweat, urine or any other fluid sample may subsequently be added. The sample is dried at 105°C for 2 hours, but may be any other
25 suitable temperature and time, and then ablated. The sample plus the 'under' matrix is ablated and carried into the plasma simultaneously. Ionization is achieved for both components and, in this way samples are calibrated. Hence, because of this, the nature of the sample is not important as the sample and the matrix containing the internal standards are introduced simultaneously to the plasma. This protocol removes the
30 necessity for a spike as the spike is already in the matrix pad on which the sample is collected. Therefore, it does not matter what the sample is, as it will be introduced into the plasma with the standards thereby overcoming any matrix interference. In this embodiment, it is not necessary to add a range of analytes to the matrix because the Be, In and Bi act as the calibrants and can be calibrated against all other elements with
35 respect to mass response before the samples are analyzed. Of course there are a series of matrices that are spiked (detailed in text already) with standards from which

calibration curves may be established thereby facilitating quantification of trace elements contained in the blood or other fluid.

Thus, fibrous cellulose matrix pads are prepared and doped with the set of mass calibration elements and dried. Blood, or other fluid is added, dried and ablated using a
5 10x10 matrix raster. The data are collected and read against results obtained from a concentration range (100, 200, 500ppb etc) of multi-element standards prepared and measured in the same way. Quantitation for any matrix may thus be achieved because the standard and sample are being introduced in the same way which therefore negates potential matrix problems. The data are cross-referenced to Be, In and Bi in the
10 standards and in the matrix with sample, and their relative values in each normalized.

The core components of the Sample Analysis System of this embodiment comprise a laser for producing an aerosol of the sample (Laser Ablation), an argon plasma, or 'electrical flame', operating at temperatures in excess of 7,000°C (Inductively Coupled Plasma) in which the aerosol is ionized, a mass filter (Mass Spectrometer) for
15 separating the ions into 'packets' according to their mass to charge ratio, and an ion detector (Multi-channel Analyzer or Ion Multiplier) for detecting the ions in each 'packet'. The system operates with a routine sensitivity capable of achieving parts per billion detection limits. All data can be electronically stored for future reference.

Suitable ICP-MS system utilizes a quadrupole mass filter, controlled by
20 alternating RF and DC fields in the quadrupole, to allow transmission of ions of one selected mass to charge ratio at any specific time. Cycling of the quadrupole allows passage of any selected ion with a mass to charge ratio of <250amu at specific times during the cycling program. Each naturally occurring element has a unique and simple pattern of nearly integer mass to charge ratio, corresponding to its stable isotopes,
25 thereby facilitating identification of the elemental composition of the sample being analyzed. The number of registered element ions from a specific sample is proportional to the concentration of the element isotope in the sample.

For multi-element analysis, the quadrupole is generally configured to scan at 1Hz (once per second). Under this circumstance, if, for example, 100 isotopic masses are
30 being analyzed, each isotopic mass will be collected only one hundredth of the entire scan time.

It will be understood that other configurations and types of instrumentation can be used with the devices and methods of the present invention without undue modification of protocols presented herein.

35 In one exemplary operation, the sample is introduced into a laser ablation cell and ablated, using either an Excimer or Frequency Quadrupled Nd-YAG laser, for a

period typically not exceeding 30 seconds. Debris from the ablated sample passes down an interface tube, made from Nalgene as a suitable plastic material but other material could also be used, attached to the torch of an inductively coupled plasma (ICP). The sample debris passes through a zone in this tube, adjacent to the torch, into
5 which independent laser radiation is being passed. A concentric series of dynode detectors measures the photon flux, reflected from the sample debris particles, which facilitates quantitation of particle scattering. Knowing the amount of scattering allows linear correlation to the amount of particles doing the scattering. The Laser scattering device is calibrated using conventional smoke cells.

10 The level of scattering is a quantitative indication of the amount of debris passing down the tube. This debris contains the sample material (blood) in addition to particles of a pre-coded (with internal standard) carrier matrix. The particles now pass on into the Inductively Coupled Plasma (ICP) where they are ionised and separated using Time of Flight (ToF) segregation. The elemental composition for the sample is
15 established and quantified with reference to the signal obtained from each of the analyte isotopes. Quantitation of the concentration of elements present in the sample and hence the blood, is calculated with reference to the scattering signal from the Laser Interferometer. The amount of sample being analysed is normalized to the signal generation by ionisation of the components in the pre-coded matrix. In this way the
20 amount of material ablated is used to obtain the mass component of the transported material and the elemental signature of the pre-coded matrix facilitates normalization of the response with reference to an ionisation efficiency cross comparison.

Quantitation of elements in the sample may also be achieved by incorporating standards into the sample or into/onto the collection matrix/support, or both. The pre-
25 coded collection matrix may contain a cocktail of elements that are not naturally present in the sample such as blood or other fluid, at levels above the detection limit of the technique. These elements typically include one or more (ie. mixture of) Beryllium, Scandium, Zirconium, Niobium, Rhodium, Ruthenium, Indium, Hafnium, Tantalum, Rhenium, Osmium and Iridium. This requires doping of appropriate analytes at levels
30 between 1 and 10,000 ng/mL to the matrix or support. The elements are chosen to cover both mass and ionisation potential ranges present in the analytically significant analytes.

In another exemplary operation, the sample is introduced into a laser ablation cell and ablated, using a Frequency Quadrupled Nd-YAG laser operating at 266 nm, for
35 a pre-determined time interval typically dictated by the number of analytes being acquired. Debris from the ablated sample passes down an interface tube, made from

Nalgene or suitable other plastic, attached to the torch of an inductively coupled plasma (ICP). The pre-coded matrix may contain a cocktail of elements that are not naturally present in blood, at levels above the detection limit of the technique. These elements typically include one or more (ie. mixture of) Beryllium, Scandium, Zirconium, Niobium, Rhodium, Ruthenium, Indium, Hafnium, Tantalum, Rhenium, Osmium and Iridium. This requires doping of appropriate analytes at levels between 1 and 10,000 ng/mL to the matrix. The elements are chosen to cover both mass and ionisation potential ranges present in the analytically significant analytes.

Readout from the spectrometer, for reporting purposes, is expressed in concentration units appropriate to clinically accepted protocols. In addition, the readout contains information on the acceptable ranges of analytes in normal healthy individuals and indicate whether the sample under investigation is below, above or in the accepted range.

The methods and devices of the present invention enable the mass screening of a variety of blood or other body fluid samples for a wide range of essential and toxic trace elements, or of samples of other fluids such as water or lubricants, for contaminants indicative of pollution or wear. Only a small volume of sample liquid (one or two drops) is required for multiple element analysis. Sample collection of body fluids does not require the use of a hypodermic needle and consequently is essentially non-invasive and considerably safer than existing methods. The sample is collected and stored in an inert matrix without need for addition of preservatives. The sample can be handled and transported safely and easily. The preferred method of analysis, quadrupole Laser Ablation-Inductively Coupled Plasma-Mass Spectrometry, is very sensitive and can detect and measure trace/ultra trace amounts of an element. The methods described herein are suited to full automation and high throughput screening and analysis of samples. Further, the methods and devices of the present invention enable multi-element testing at a significantly lower cost than many current single element tests, thus making the economical mass-screening of target populations possible.

Examples of suitable internal standards which may be used for quantitation of elements, in conjunction with the devices and methods of the present invention, are detailed in Table 1 below.

Table 1:

Sample Name	SARM 1	SARM 3	SARM 46	SY-2
Alt. Name	NIM-G	NIM-L	S14	
Sample Type	Granite	Lujavrite	Stream Sediment	Syenite Rock

	ppm	ppm	ppm	ppm
Si	353848	244936		280975
Ti		2878		899
Al	63933	72190		63722
Fe 3+	4197	61410		16988
Fe 2+	10105	8784		27672
Mn	155	5963		2478
Mg	362	1689		18222
Ca	5575	23013		66889
Na	24926	62093		31974
K	41424	45741		36942
P	44	262		1877
Ag				0.029
As	19.3	1.92		17.3
Au	0.0011	0.00084		0.00052
B				88
Ba	120	450		460
Be	7.75	29.5		22
Bi	0.275	0.468		0.111
Br				
Cd	0.113	0.91		0.21
Ce	195	240		175
Cl	263	1200		140
Co	0.36	2.44	54	8.6
Cr	12	10	593	9.5
Cs	1.08	2.78		2.4
Cu	12	13	563	5.2
Dy	17	3.1		18
Er	10.5	2.6		12.4
Eu	0.35	1.2		2.42
F	4200	4400		5030
Ga	27	54		28
Gd	14	3.8		17
Ge		0.89		1.3
Hf	12.4	231		7.7
Hg	0.0189	0.0445		0.0043
Ho	3.6	0.9		3.8
I				
In				
Ir	0.0005			0.0005
La	109	250		75
Li	12	48		95
Lu	2	0.4		2.7
Mo	2.84	1.21		0.53
N				
Nb	53	980	26	29
Nd	72	48		73
Ni	8	2.2	122	10
Os				

Pb	40	43	14000	85
Pd	0.007			0.015
Pr	19.5	18.4		18.8
Pt				
Ra				3.7
Rb	325	190	18	217
Re				
Rh				
Ru	0.01			0.002
S		650		160
Sb	1.19	0.13		0.26
Sc	0.9	0.5		7
Se	0.012	0.014		20
Sm	15.8	5		16.1
Sn	3.3	7.4		5.7
Sr	10	4600	28	271
Ta	4.9	25.2		2.01
Tb	3	0.7		2.6
Te	0.007	0.009		0.002
Th	51	66		379
Tl	0.93	0.325		1.5
Tm	2			2.1
U	15	14		284
V	2	81	195	60
W	1.45	8.28		0.76
Y	143	22		128
Yb	14.2	3		17
Zn	50	395	6200	248
Zr	300	11000	95	280

The collection matrix, if one is used, may be impregnated with a trace metal cocktail, of known concentration using purpose prepared aqueous solution standards. In certain preferred embodiments, the matrix may contain 2ppm of Be, In, Hf as internal standards to calibrate the mass response for the system in blood analysis. In other embodiments describing wear metal analysis of oil, 2ppm of Be, In and Th may be used. In yet other embodiments, different suites of elements may be used.

Separate standard matrix pads may be used to calibrate the sensitivity and these may be as follows for blood and body fluids: a single pad containing, but not restricted to, Li, Na, Mg, Al, P, K, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Rb, Sr, Mo, Cd, Sn, Sb, Te, Ba, La, Ce, Eu, Dy, Yb, Hg, Tl, Pb, Bi, Th and U at 1 ppb, a second pad with all these at 2 ppb. A third pad with all of these at 5ppb a fourth pad with all of these at 10ppb a fifth pad with all of these at 20 ppb a sixth pad with all of these at 50 ppb a seventh pad with all of these at 100ppb an eighth pad with all of these at 200ppb a ninth pad with all of these at 500 ppb a tenth pad with all of these at 1000ppb. An appropriate

concentration can then be used for the set of elements being determined in a particular fluid sample. In another embodiment, a suite of elements appropriate to wear metal analysis in oil, for example, Li, B, Mg, Al, Si, P, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Sr, Y, Zr, Mo, Ag, Cd, Sn, Sb, Ba, La, Ce, Hf, Hg, Pb and U may be doped into matrix pads at 1ppb through 1000ppb as above, so that when ablated, a range of elements across the mass spectrum may be used as internal standards to standardise the system. Thus, the collection matrix, when used, may contain a pre-calibrated concentration of selected analytes. Both a broad-spectrum general collection matrix/device and a test specific matrices/device/s may be employed for specific elements or suites of elements. Further, any one, or combination or range of internal standards analytes may be spiked into the collection device to ensure its broad spectrum or specific use. For example, for broad spectrum, the preferred combination is , Li, Na, Mg, Al, P, K, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Rb, Sr, Mo, Cd, Sn, Sb, Te, Ba, La, Ce, Eu, Dy, Yb, Hg, Tl, Pb, Bi, Th and U and for specific applications, for example analyzing oils preferred is , Li, B, Mg, Al, Si, P, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Sr, Y, Zr, Mo, Ag, Cd, Sn, Sb, Ba, La, Ce, Hf, Hg, Pb and U and for blood the preferred combination is , Li, Na, Mg, Al, P, K, Ca, Ti, V, Cr, Mn, Fe, Co, Ni, Cu, Zn, Ga, As, Se, Rb, Sr, Mo, Cd, Sn, Sb, Te, Ba, La, Ce, Eu, Dy, Yb, Hg, Tl, Pb, Bi, Th and U.

A typical procedure of collecting and analyzing a sample is summarized in Figure 5. Of course, manual procedures can also be adopted, as can variations of the proposed exemplary scheme.

Example 4: Analysis of collection matrices

The purpose of the experiments described below was the definition and/or refinement of chemically and mechanically robust fluid adsorption/absorption matrix/matrices to facilitate the collection and quantitative analysis of micro-litre fluid samples by Laser Ablation-Inductively Coupled Plasma Mass Spectrometry (LA-ICP-MS). For purposes of this example fluids under consideration are blood, urine and oil. However it will be understood that any other fluid, biological or otherwise, may be analysed using similar matrices and techniques.

Preferably the sample collection matrices should be suitable for incorporation into a robust, transportable sample collection device. The device should have specific attributes such as but not limited to:

- be cheap and capable of precision mass production;
- be small and easily accommodated in laser cells for ablation prior to analysis;

- be able to be coded for automatic pre-analysis reading and referral of the sample back to the data, and the data to the client;
- for blood collection, contain a mechanism for penetration of individual patient's skin thereby minimising potential 'stick injuries'. There would be some form of shielding device, or mechanism, that would "shield" the puncturing mechanism such that it would not be able to penetrate the skin of another person subsequent to initial collection of blood;
- produce minimum biohazard with material after analysis and prior to disposal. This implies a small collection device and a small blood sample (less than 100 μ L), and a very small amount of material comprising the sampling device itself that would ultimately have to be incinerated;
- easy transportability to and from the collection site and through conventional mailing procedures. The device should be such that conventional postal systems can be used without the possibility of contamination and release of potentially bio-hazardous material; and
- be capable of being used by non-medical personnel.

MATRIX MATERIALS

The original preferred matrix material used for process testing was fibrous cellulose. Using this material, it was possible to readily form backed cardboard 'punch-outs' containing the cellulose absorptive medium. Micro-litre samples of blood, added to this material, were qualitatively analysed by LA-ICP-MS. Qualitative spectra and raw count data were generated, much of which reflected trace metals in the absorbed blood. However, it was reasoned that the cellulose, being a natural organic product, might be contributing to the analyte signal of a range of elements recorded. Hence, it was determined that cellulose, together with an array of other potential matrix materials, be further investigated, both in terms of its chemical and physical characteristics.

Some attributes of suitable sample collection matrices include but are not limited to:

- must be chemically "clean", that is, have a low concentration of analytes of interest;
- robust, that is, capable of transportation, often over long distances without fragmentation;
- have significant wettability, both by aqueous and non-aqueous (blood and oil) samples while still retaining integrity;
- capable of withstanding laser ablation removal of samples; and
- not contribute to analyte segregation during analysis.

MATRIX CHOICE

The parameters detailed above govern the choice of matrix and, as such, preclude certain materials. A list of matrices investigated follows with indications as to their potential suitability, or otherwise, which resulted in a final short list of potentially useful material to be subsequently tested. The choice of white metal oxides as potential matrices is based on the fact that the two detailed herein are locally manufactured in bulk, are extremely cheap and, using the modern generation of UV lasers (unlike IR lasers), are customarily considered not to have variable coupling efficiencies between light and dark matrices.

Potential organic and inorganic matrix materials investigated are:

- Pig-toe mussel shell (aragonite) – sourced from the WA pearl industry
- Aluminium hydroxide – Alcoa (WA)
- Titania – New Millennium (WA)
- Bacterial grade glucose – sourced by Professor Watling
- Starch "A" - BDH Analar analytical reagent
- Starch "B" – Ajax Chemicals Univar analytical reagent
- Glucodln – Boots Healthcare Australia
- Cellulose – high purity powder – Sigma Chemicals Microgranular
- Cellulose – high purity fibrous cellulose – Sigma Chemicals Medium Fibrous
- Hydroxy Butyl Methyl Cellulose – Sigma Chemicals
- Flour – rice, maize, wheat, soy, rye and corn flour commercially available grocery lines

All of the above matrices can be used for lubricants where the levels of metals are much higher. However, the following are particularly useful choices of matrices for blood and other body fluid analysis, which can also be used for analysis of lubricants or water samples.

Aluminium hydroxide [Al(OH)₃]: A very high quality aluminium hydroxide is produced in Western Australia. It is analytically relatively clean and cheap, and is being considered as a matrix.

Cellulose: Cellulose is an excellent theoretical matrix choice in that it is typically low in heavy metal concentration. A variety of ultra-pure cellulose was tested for compactability, wettability and metal content. The physical characteristics of cellulose as such (it was the original matrix) make it important material as a potential matrix. Particularly useful is fibrous cellulose in the form of cellulose filter papers (Whatman

540, but also 541, 542 and other cellulose filter papers, Whatman International Ltd, Maldstone, England).

Flour: Newly acquired rice flour has proved exceptionally robust under wetting and drying conditions and may also be advantageously used as a matrix.

5 In addition to simply using the matrix material as supplied, relevant matrices were leached and the leached residue tested to see if significant metals could be leached, thereby reducing the metal content of the matrix and possibly rendering it more useful by lowering the level of contaminant metals, or actually reducing the level of metals in the sample to a level where previously unsuitable material would now be suitable.

10 EXPERIMENTAL

(i) Chemical Characterisation

Solution ICP-MS: In order to assess the 'purity' of the respective potential matrices, appropriate sub-samples of water-soluble materials were dissolved in Milli-Q (mQ) water and made to volume. Water-insoluble samples, (primarily the inorganic materials) were subjected to both cold and/or hot (or both) hydrochloric, nitric, aqua regia and nitric-hydrofluoric acid leaches. The leachates were recovered, made to volume, appropriately diluted and analysed by solution introduction ICP-MS. The leached residues were recovered and a selection of sub-samples subjected to total dissolution followed by solution ICP-MS analysis using a VG PlasmaQuad 3 ICP-MS made by VG Elemental, Ion Path Road 3, Winsford, Cheshire CW7 3BX, United Kingdom. Further selected residue sub-samples, along with unleached equivalents, were subjected to total acid dissolution, made to volume, diluted and again analysed by solution introduction ICP-MS.

The solution experiments facilitated elimination of several of the potential matrix candidates, having unacceptable concentrations of analytes of interest in the raw material and analytes little, or not adequately, reduced by acid leaching. The 'solution' assessment indicated that cellulose and aluminium hydroxide were the best candidates but that both of these may contain certain analytes of interest. Because of the need to dilute the solutions for ICP-MS analysis, very low apparent concentrations in solution frequently translated to significant concentrations in the sample when corrected for mass and dilution; in many cases, these analytes may not be present or, if present, present at very much lower concentrations. To test this thesis, 'raw' sub-samples, and corresponding leached residues where applicable, were pressed into 'briquettes' (see below) and subjected to comparative qualitative UV LA-ICP-MS analysis.

35 **Laser Ablation ICP-MS:** It is not necessary that the sample matrix will contribute an equivalent amount of material to the analytical sample as the blood or other fluid.

The incorporation of the matrix and its ionisation will not be equal to that for the blood contained in it. Because of this, the contribution of matrix to the analytical signal will not necessarily be in proportion to its relative matrix/blood ratio. Hence, it was necessary to determine what relevant contribution the matrix has to the analytical signal during a real analysis. Laser ablation analysis of the matrix was therefore also undertaken. Because the use of argon as a carrier gas is the traditional method of transport of ablation debris to the plasma this was the initial gas used for all experimental purposes. However, helium is finding an increased following in the scientific community as a transport gas as it often gives improved sensitivity and reduced isobaric interferences. Consequently this gas was also investigated.

(II) Physical Characterisation

Physical characterisation of potential matrix materials included assessment of compaction integrity, both at 500 and 1000 kg/sq in, wettability to blood and aqueous solutions, integrity after sample addition, contrasting behaviour of single and multi-component matrices, and internal standard introduction. Results from some of these investigations are detailed below.

The use of an internal standard is necessitated because of the variability in ablation efficiency between samples. There is no way of controlling the "fluence" variation (variation in the efficiency of coupling and hence power transfer of the laser energy to the sample) from sample to sample. Because of this, varying amounts of analyte will reach the plasma depending on the relative fluence between samples. Consequently, it is necessary to ensure that there is a mechanism for estimating the amount of material being transported to the plasma for each sample. The method used for an infrared laser was to measure the scattering of light by the transported particles. However, this mechanism is not possible when a UV laser is used (the laser used for these experiments was a frequency quadrupled Nd-YAG UV Microprobe Laser System operating at 266nm in pulsed Q-switched mode. The Laser System was manufactured by VG Elemental, Cheshire, United Kingdom.

However, spiking a simple element cocktail into the matrix, either prior to, or concurrent with, sampling provides a useful and inexpensive internal standard for quantification experiments.

RESULTS AND DISCUSSIONS

Details of eighteen experiments completed during the period October-December 2002 are set out below. Sixteen of the experiments relate specifically to physical and chemical characteristics of the matrix, and analysis of absorbed aqueous standard, mineral CRM and blood samples. The remaining two experiments, Experiments 13 and

15, deal with the analysis of oil samples – these are reported together at the end of this section.

The resulting analytical data is presented in a series of Appendices identified by experiment number, for example, 'Appendix Experiment 12'. These appendices should
5 be viewed in conjunction with the relevant commentary on the individual experiments as contained herein. Frequently, averages of data and % standard deviations (coefficient of variations) have been computed.

In most appendices, isotopic data has been computed to 100 per cent elemental concentration using natural isotopic abundance relations. In a small number of cases,
10 data is presented solely as isotopic concentrations at the measured isotopic mass. This is clearly indicated in the respective appendices.

In an attempt to optimise signal response, peak hopping instead of normal scanning acquisition was employed. Under this analytical regime, data acquisition at each isotopic mass occurred on three channels only. Not uncommonly, transient
15 electronic spikes may be recorded on one of the three channels. The on-board computer processes the data from all three channels and reports the results as raw count 'concentrations'. Where a measurement includes a transient spike, the resulting raw counts for that analyte may be considerably elevated relative to duplicate or replicate analyses of the equivalent analyte in the same sample. This leads to often-
20 marked concentration contrasts for specific analytes in these samples. The problem may be overcome by increasing, to say seven, the number of channels over which individual isotopic mass data is collected. Under these circumstances, a normal 'smoothing' algorithm may be automatically applied across the seven channels to produce precision results for duplicate or replicate analyses. Having established this as
25 being a major cause of analyte variability, analytical protocols have been appropriately modified to allow data collection over the increased number of channels.

Another cause of analyte variability may be due to possible surface 'contamination' of the collection matrices. To minimise contamination, the top pad of a matrix wad has been removed so that there is no airborne contamination on the surface
30 to be analysed. In an embodiment of this process, the matrix pads are prepared in a sterile, dust-free clean room, enclosed in a container which may only be breached immediately prior to sample collection. Improved analytical precisions, following implementation of this protocol, are attributed to the sample preparation

Correction of data for identified transient spikes had led to a marked
35 improvement in analyte reproducibility and, hence, 'precision' data.

Example 5: Matrix And Blood-Related Experiments**Experiment 1**

The aim of this experiment was to develop and test procedures to produce 3 mm diameter test tablets as a prelude to physical characterisation of sample matrices. For this purpose, an XRF pressed powder vacuum press was modified, and new dies manufactured, to facilitate pellet production. Matrix materials chosen for the inaugural production tests were glucose, cellulose and a 1:1 mixture of the two; initial compaction pressure was 500kg/sq in. Initial physical and chemical investigations were undertaken concurrently until preferred matrices were identified.

Pelletising of glucose required the use of weighing paper between sample and metal on the press die. Absorption of liquid appears good.

Cellulose pelletised quite well, with very good strength. However, fluid absorption was slow. A 1:1 mixture of glucose and cellulose powder pelletised well without the need for weighing paper between pellet and die. Pellet strength was improved over glucose alone and fluid absorption was intermediate between rates for glucose and cellulose powder pellets compacted at equivalent pressure.

Experiment 2

The principal objective in this experiment was to assess the chemical purity of a range of potential matrix materials. Sample preparation for analysis was undertaken concurrently with pelletising press modifications. Various matrices, including pig-toe mussel shell, glucodin, glucose, cellulose, hydroxy butyl methyl cellulose (HBM cellulose), TiO_2 and $\text{Al}(\text{OH})_3$ were leached, dissolved or digested in preparation for solution ICP-MS purity assessment.

Method

Pig toe mussel (Sample A, B, C and D) - ~1.5g pearl seed taken, dissolved in 20mL 1:1 HCl:mQ water, then taken to dryness. 4mL of HNO_3 :mQ 1:1 added, heated and made up to 100mL with mQ water. Diluted x20 with mQ (2ppb Ir, Rh) water for ICP-MS.

Glucodin (Sample E and F) + Glucose (Sample G) - ~1.5g Dissolved in 100mL of mQ water. Diluted x5 for ICP-MS.

Cellulose (Sample H) + HBM Cellulose (Sample I) - ~0.5g digested in 20mL cHNO_3 for 36 hours, reduced to 10mL and made up to 100mL with mQ water. Diluted x5 for ICP-MS.

TiO_2 (Sample 001) + $\text{Al}(\text{OH})_3$ (Sample 003) - Leached with 1:1 HCl:mQ water for 36 hours, decanted and washed 3 times with mQ water (~20mL). Decanted solution (leachate) made up to 100mL with mQ water. Diluted x10 for ICP-MS.

TiO_2 (Sample 002) + $\text{Al}(\text{OH})_3$ (Sample 004) - Leached with 1:1 HNO_3 :mQ water for 36 hours, decanted and washed 3 times with mQ water (~20mL). Decanted solution (leachate) made up to 100mL with mQ water. Diluted x10 for ICP-MS.

Residues were dried and saved for LA-ICP-MS.

5

This experiment was concerned with the determination of the trace element concentrations in prospective matrices for blood (and other fluid) collection, together with looking at some of the results of leachates of titanium dioxide and aluminium hydroxide.

10

The results for the leachates are detailed (Appendix Experiment 2). It may be possible to indicate that aluminium is obviously leached from the aluminium hydroxide matrix, but also from the titanium dioxide matrix, and conversely titanium is leached from the titanium dioxide matrix and there is also some indication of leaching of titanium from the aluminium hydroxide matrix. In the case of titanium dioxide, HCl appears to be more aggressive than HNO_3 , whereas the reverse is the case for the aluminium hydroxide. Concentrations of manganese, copper, strontium, zirconium are found from the leachates of both matrices while zinc, rubidium, barium and lead appear to be quite concentrated in leachates from the titanium dioxide matrix. In the aluminium hydroxide matrix, tin, gallium, zirconium, hafnium and uranium appear to be present in leachates from this matrix.

20

Total digest and/or solubilization data of pig-toe mussel, glucodin, glucose, cellulose and HBM cellulose are also presented in Appendix Experiment 2. The pig-toe mussel contains significant concentrations of lithium, aluminium, titanium, manganese, copper, zinc, rubidium, strontium and barium. While this would imply that the matrix is not suitable as a blood collection matrix, because of the concentration of these elements, it is also necessary to analyse the pig-toe mussel material with sample attached under laser ablation conditions rather than solution conditions to make sure that these elements are also carried over by laser ablation and not just present in total digests. In the case of glucodin, glucose, cellulose and HBM cellulose all contain significant amounts of aluminium, titanium, chromium, manganese, nickel, copper, zinc, rubidium, strontium and barium while cellulose matrix alone, in addition to containing these elements, also contains significant concentrations of lead and bismuth; both cellulose and HBM cellulose also contain concentrations of zirconium, tin, thallium and thorium not found in the glucodin and glucose.

25

30

35

Although these matrices all contain significant amounts of trace elements in the ppb range, this does not necessarily preclude them from use as a sample collection

matrix as conventional blank correction can be used to overcome problems associated with blank content. This can be further emphasised by the fact that inter-element ratios could be used to determine, and to augment, blank corrections by looking at relationships between metals and tracing these through to the final analytical protocols

5 **Experiment 3**

The purpose of this experiment was to further test, the pelletising and adsorption characteristics of cellulose powder, glucose, and starch, and mixtures thereof, and to check the dissolution/absorption characteristics of the pellets by SY-2 (mineral CRM, , Canadian Certified Reference Material Project (CCRMP), Table 1 solution. The results of Experiment 3 are set out in Appendix Experiment 3

10 Cellulose powder alone works well. The glucose undergoes surface dissolution leaving holes on the surface. The starch absorbed water and expanded, causing the surface to bulge. Under the pelletising pressure of 500 kg/sq in, the cellulose powder is tightly compressed and it takes some 10 to 15 seconds for fluid absorption. This suggests that a more fibrous cellulose with an 'open' structure may be preferable. To this end, further experimentation with fibrous cellulose is indicated. In addition, further experimentation with powdered cellulose at differing packing pressures is warranted.

15 **Experiment 4**

The aim of this experiment was to assess the absorptivity and mechanical stability of cellulose powder pellets compacted under differing pressures. In the first instance, powdered cellulose was suspended in mQ water and vacuum filtered. The collected filter cake was mechanically incoherent. This caused it to flake and fall apart. However the adsorption of solution was rapid.

25 Cellulose powder compacted under a pressure of 100kg/sq in, while mechanically robust, still absorbed slowly. At low compaction pressure, estimated to be about 50kg/sq in and achieved by turning the tightening screw on the press just until there was resistance, the resulting pellets illustrated rapid absorption. Furthermore, the pellet holds together well. The experiment appears to confirm that compaction destruction of porosity rises with increasing pressure thereby rendering the matrix progressively less absorptive.

30 **Experiment 5**

The aim of this experiment was to quantitated trace elements in a blood sample using internal standards. The experiment also tested the absorption of SY-2 (mineral CRM) and blood onto cellulose pellets, robustness of the doped pellets when subjected to LA-ICP-MS analysis, assess levels of possible contaminants, evaluate results arising from the doped matrices and assess the comparability between 'wet' and 'dry' matrices.

35

The following instrument settings were used: Lens voltages – Lens 1, 2, 3, and 4 respectively –10.8, -22.6, 0.7 and –13.3 Volts, Collector – 4.6 Volts and Extraction, -332 Volts; Gas Flows – Cool gas 13.6 L/min, Aux gas 0.81 L/min Neb gas 0.74 L/min and Oxygen gas 0.00 L/min; Torch box positions – X, Y and Z axes respectively 932, 165
5 and 250 steps; Multiplier voltages – H.T. pulse count –2634 Volts and H.T. analogue) Volts; Miscellaneous settings – Pole bias –2.2 Volts, R.F. power 1500 Watts, Perl speed 0%; PlasmaScreen is OUT, S-Option pump is OFF.

Samples of blood were obtained from a subject with the aid of a SoftTouch
10 lancet device (used for home blood glucose testing and manufactured by Boehringer Mannheim, Germany) applied to a pre-cleaned (absolute ethanol wiped) area of a fingertip. Successive drops of blood were encouraged to form through application of pressure. The drops were directly 'touch' applied to 3mm diameter by 2mm deep sample collection matrix tablets formed by pressing granular cellulose (Sigma
15 Chemicals Microgranular powder) under a load of 500 kg/sq. in. The matrix tablets were affixed to a Perspex disc, 37.5 mm in diameter and 6mm deep, fabricated from Perspex rod, using 3M Scotch Permanent Double Stick Tape. The volume of the drops was estimated to range between 30 and 70 microlitres. No preservatives or anticoagulants were used and there was no requirement to store the blood prior to
20 application to the collection matrix, or subsequent analysis. However, there is provision for loaded sample collection matrix tablets to be refrigerated and stored following oven drying at 60°C for one hour.

Four blood samples were prepared; two were oven dried and two were maintained "damp". Duplicate sets of equivalent SY-2 CRM-doped (Syenite, Canadian
25 Certified Reference Material Project) matrix pellets were prepared by pipetting 50 µL of the standard solution onto the respective matrix tablets and drying thereby generating matrix matched standards. The SY-2 CRM contains calcium, iron, magnesium, potassium and so forth and this provides a high ion flux that is possibly equivalent to the ion flux expected of blood. Hence, any ion effects that were taking place would be
30 comparable in the blood and SY-2, as compared with a straight aqueous standard solution.

The sample holder, with affixed blood- and CRM- doped matrices was placed into the laser ablation cell of the UV Microprobe Laser System attached to a VG PlasmaQuad 3 ICP-MS both manufactured by VG Elemental, United Kingdom. The
35 laser is a frequency quadrupled Nd-YAG operating at 266 nm; 10x10 matrix raster

ablation of the samples was undertaken in pulsed Q-switched mode at a fluence of 6.2 millijoule for 60 seconds.

The output data was acquired as raw counts from on-board software and exported into Excel and manipulated. No algorithms were used for computations. The raw count data for both blood and CRM samples were matrix blank corrected by subtracting the averaged matrix blank value from the individual blood and SY-2 values. From these corrected data % Standard Deviations were computed as a measure of precision. Finally, trace element compositions for the 11 analytes examined in the exemplary run were computed with reference to matrix matched SY-2 CRM values.

Data obtained is set out in Appendices Experiment 5A and 5B.

As indicated above, part of the experimental design was to determine whether it was necessary to fully 'dry' the sample prior to analysis. Collection of blood onto a matrix without the drying step as detailed above, may lead to a sample being slightly damp. Hence, it was necessary to determine whether variation in the moisture content of the matrix would affect the readout of concentration of elements in the matrix. Consequently two sets of samples of cellulose were set up and, in addition to 'wet' and 'dry' blood, SY-2 certified reference material doped samples were also prepared in an attempt to quantify the concentration of metals in the blood. Blood samples and SY-2 were spiked onto cellulose in duplicate and one set of blood samples was analysed 'wet'. A second subset was taken and dried (as above) and the samples were analysed dry. Data from these experiments is also presented in Appendix Experiment 5A

Following analysis, results for the wet samples were blank corrected and data produced. Simple inspection of the data for the 'wet' blood samples indicates relatively high variability in analyte concentrations particularly in the case of lead and zinc where a variation of $\pm 100\%$ is recorded. The analysis of SY-2 certified reference material is far more uniform.

For the dry sample, the results are better. Reproducibility is improved and results are more uniform. From the blank corrected values for the dried blood sample it can be seen that, with the exception of barium, the results are meaningful. Barium results go negative and this is probably due to the fact that the barium signal is small relative to the blank – the blank is quite high. However, both lead and zinc are much improved and, if these are used to calculate concentrations of these elements in the blood, based on SY-2 concentrations (calculated in Appendix Experiment 5B) the blood values and expected blood values from the literature are quite close for the analytes under consideration. SY-2, a certified reference material, has been used for a number of reasons. First, use of simple aqueous solution on the collection matrix would not, on

ablation, have provided a significant ion flux. The SY-2 contains calcium, iron, magnesium, potassium etc (see Table 1) and this provides a high ion flux that is possibly equivalent to the ion flux of the blood. Hence, any ion effects that were taking place would be comparable in the blood and SY-2, as compared with a straight aqueous solution. Thus a normal CRM, that has a relatively high matrix concentration will suffice.

The above experiment, including instrument settings and internal standardisation as described, is equally applicable to simpler biological fluid samples such as components of whole blood (eg. serum or plasma), urine, sweat, tears, cerebrospinal fluid and the like. The sample collection, handling and analysis of such fluids is simpler and thus greater accuracy can be achieved.

Experiment 6

This experiment was conducted to analyse the titanium dioxide and aluminium hydroxide matrices, both before and after leaching (leached residues from Experiment 2). The data produced in this experiment ties in with the leachate data from Experiment 2. Upon total dissolution, solutions derived from titanium dioxide have very high concentrations of titanium, while those derive from digestion of aluminium hydroxide are similarly rich in aluminium. Accordingly, these two elements have not been measured.

The purpose of the experiment was to evaluate the efficacy of acid cleaning of the white oxide matrices. Hence, appropriate sub-samples of 'raw' titanium dioxide and aluminium hydroxide, together with their hydrochloric- and nitric acid-leached equivalents, were digested in a sulphuric/hydrofluoric acid, made up to volume, diluted and analysed by solution introduction ICP-MS. The leachates derive from HCl- and HNO₃-leaching of bulk titanium dioxide and aluminium hydroxide were analysed in Experiment 2 and the results reported in Appendix Experiment 2.

The comparison of the "raw" original material and the HCl- and HNO₃-leached residues show that, for titanium dioxide, its HCl-leached residue and associated leachate, weak to strong leaching of lithium, manganese, copper, zinc, gallium, rubidium, strontium, (zirconium), barium, lead, (thorium) and uranium has been achieved. Here, there is generally a good mass balance between concentration in the original versus the sum of concentrations in the leachate and leached residue. In contrast, concentrations of vanadium, chromium, nickel, germanium, yttrium, zirconium, niobium, tin, antimony, hafnium, tantalum and tungsten in the raw material are unaffected by HCl-leaching.

For titanium dioxide, its HNO₃-leached residue and associated leachate, weak to strong leaching of lithium, (chromium), manganese, copper, zinc, gallium, rubidium, strontium, (zirconium), barium, lead and (thorium) is evident. In contrast,

concentrations of vanadium, (chromium), nickel, germanium, yttrium, niobium, tin, antimony, hafnium, tantalum, tungsten, (thorium) and uranium are little or unaffected by HNO_3 -leaching.

Turning to the aluminium hydroxide matrix, HCl and HNO_3 both have a similar
5 leaching response with both acids weakly to strongly leaching all elements occurring in significant concentrations in the aluminium hydroxide matrix. The elements involved are lithium, beryllium, chromium, manganese, copper, gallium, strontium, zirconium, tin, hafnium, thorium and uranium. Hence, use of these acids to pre-clean the matrices is recommended. Both can be leached quite easily in both HCl and HNO_3 .

10 Of particular importance is the presence of gallium in the aluminium hydroxide matrix. A small amount is acid-leached but this does not impact its potential of being used as an internal standard; the same holds true for zirconium. Although not as high as zirconium in the titanium dioxide matrix, zirconium in aluminium hydroxide could still be used for a double internal standard based on gallium and zirconium. There is a
15 possible problem with the aluminium hydroxide matrix in that there is copper in it but the copper tends to be relatively uniform and if copper results in previous analyses are considered, reasonable results for copper are obtained by doing blank corrections. It should be remembered all the time that although these metals are present in the matrix, they may not contribute an equivalent amount to the determination of metals in blood
20 because they are not transported as much as the blood to the plasma. The blood tends to fill interstices and sit on top of the matrix; hence, these elements may not contribute a significant amount to the concentrations that are present in analysed, so-called blood.

This experiment demonstrates that it is possible to variably reduce and/or eliminate a range of trace elements from titanium dioxide and aluminium hydroxide
25 matrices. When combined with previous experiments, it would appear that possibly two matrices, aluminium hydroxide and cellulose, may constitute particularly suitable matrix materials.

Experiment 12

The purpose of this experiment was to examine the efficacy of a fibrous
30 cellulose mat (Whatman 540 filter paper, Whatman International Ltd) as a sample collection matrix. This material is an efficient absorber of fluids, but its 'coarse' fibrous texture may result in variable ablation characteristics. Six duplicate sub-samples of the cellulose mat were taken and pre-prepared as follows: Two duplicate sets were rinsed for 10 minutes with 50% aqua regia and dried; a further two duplicate sets were washed
35 overnight in aqua regia and dried while the remaining duplicate sets were left unwashed. One set each was doped with 2ppm multi-element standard and dried whilst

the second set of each was retained as blanks. It was observed that the fibrous cellulose mat, rinsed for 10 minutes with aqua regia, upon drying was rendered 'harder' than the other two (unwashed and overnight washed) mats.

The blanks and doped equivalents were analysed by LA-ICP-MS and the results of analysis are recorded in Appendix Experiment 12. Upon ablation, it was observed that for the 'hardened' rinsed matrix, the laser penetrated through the whole mat, whereas for the other two, the laser did not penetrate all the way through. This observation clearly implies that the contrasting physical characteristic of the fibrous cellulose mat impact upon laser penetration and, hence, lasing characteristics. With reference to the relevant Appendix, pages Experiment 12/3 and 12/4, it is clear that, for cerium-normalised data, data for the 'hardened' rinsed fibrous cellulose mat, which exhibited complete laser penetration, gives rise to the best overall precision data. Indeed, most analytes have precisions of less than 10% and frequently less than 5%. This outcome further emphasises the potential value of fibrous cellulose as a matrix material.

Experiment 16

The objective of this experiment was to evaluate potential sensitivity improvements for aqua regia and ammonium fluoride (NH_4F) doped 3:1 $\text{Al}(\text{OH})_3$:cellulose matrices.

From a 3:1 $\text{Al}(\text{OH})_3$:cellulose mixture, six triplicate sets of pressed pellets were prepared. These unwashed triplicate pellet sets were affixed to a Perspex disc. One set was left 'blank' and a further set was doped with 1ppm multi-element standard; both were oven baked. Two of the remaining four triplicate sets were doped with 5 μL of 50% aqua regia and oven at 105°C for 2 hours; the remaining two triplicate sets were doped with 5 μL of 1M ammonium fluoride (NH_4F) and oven baked. One set each of the aqua regia and ammonium fluoride treated pellets were further doped with 1ppm multi-element standard and dried.

A further sample of the 3:1 $\text{Al}(\text{OH})_3$:cellulose mixture was washed with aqua regia, rinsed and dried. This material is referred to as the washed matrix. From this washed matrix, equivalent triplicate sets of pellets were prepared as for the unwashed matrix described above. It was observed that the 50% aqua regia doped matrices were not as mechanically robust as other matrices prepared in this experiment. All triplicate sets were analysed by LA-ICP-MS. The results for the unwashed matrices are presented in Appendix Experiment 16A while those for the washed matrices comprise Appendix Experiment 16B.

When results for unwashed material, that is, no aqua regia wash, are considered, it is apparent that the results are significantly better for unwashed, than for the washed, material. For blank corrected matrices, normalised to cerium, precisions for the unwashed material are better than those of the washed matrix. This outcome suggests that there is no fundamental need to wash 3:1 $\text{Al}(\text{OH})_3$:cellulose matrix.

Disregarding, the blank corrected, cerium normalised data for the present, and considering only the 'raw' 1ppm doped matrix data, the recorded precision measurements for both unwashed and washed matrices show a general improvement in the NH_4F doped matrices. This apparent improvement in sensitivity may result from improved ablation of the matrix possibly through production of a more volatile atmosphere in the presence of NH_4F .

Experiment 18

The several previous experiments have sought to identify appropriate clean matrix materials together with preferred compaction, absorption, ablation and pre-treatment characteristics. Particularly preferred matrix and analytical conditions for most test samples, and particularly useful for blood and other body fluid samples, were identified as Whatman 540 filter paper, ablated at 10Hz at a fluence of between 4 and 9 Millijoule with a flow of argon between 900 and 1000mL per minute.

In the course of this work, consideration was given to the question as to whether it may be possible to prepare a blood sample in such a way that it was matrix supported, rather than matrix absorbed. If this could be achieved, then it may be possible to ablate blood samples free of matrix. In this way, analytes present in the analysis would be derived from the blood alone. Consideration of direct analysis of supported, rather than matrix-absorbed blood, arose from the observation that, during the experimental procedures segregation of blood serum and plasma appeared to occur. The observed probable segregation was not considered to be a significant problem; the laser ablation protocol was designed in such a way that the laser would penetrate through any dispersion front in the matrix, thereby sampling any segregated blood and consequently 're-assembling' or re-combining the analyte cocktail. Nonetheless this observation suggested that it might be possible to overcome any potential matrix interference by ablating only dried blood.

It was reasoned that if a shallow, 3mm diameter, 125 micron deep, depression was cast into the surface of the matrix pellet, then a drop of blood delivered to the depression would flow to fill the depression and present a flat surface away from the depression lip (meniscus) for subsequent lasing. A requirement would be that no chromatographic segregation of serum and plasma occurred. To this end, it was further

reasoned that if the 3:1 Al(OH)₃:cellulose powder was compacted under high pressure (at least 1 tonne/sq in), then the matrix may be rendered effectively impervious and simply support blood as it coagulated and dried.

Consequently, a new die for the vacuum press was fabricated to produce a 6mm
5 diameter pellet into which was impressed a 3mm diameter by 125 micron deep, flat bottomed circular depression. An appropriate number of new pellets were pressed at 1 tonne/sq in pressure.

Micro-litre samples of blood were delivered to, and contained within, the surface depressions on the surfaces of ten matrix pellets; five of these pellets were air dried at
10 ambient temperature and the remaining five oven dried at 60°C. A further two blood drops were applied to the Perspex mounting disc and dried. Here, the surface of the dried blood drops was not flat, but rather, strongly undulating.

On application, it was clear that some plasma segregation and absorption occurred, causing a volume increase and expansion in the tightly compressed cellulose
15 powder. However, the pellets retained sufficient mechanical integrity to allow LA-ICP-MS analysis. When ablated, the 'serum' tended to fragment in 'chunks' giving rise to somewhat variable results. Notwithstanding, the counts obtained were reasonable for most elements.

For the matrix free blood drops, dried onto the Perspex support, the ablated
20 blood was far more coherent, with nice ablation. However, as noted above, the surface was strongly undulating leading to changed laser focal conditions and, hence, non-optimal results.

Given that the aluminium hydroxide:cellulose matrix was not impervious, the matrix free approach described above can be adopted, i.e. use impervious substrate,
25 such as Perspex, into which 3mm diameter by 125 micron deep circular impressions have been pressed, moulded or machined. Each sample collection device can contain two such depressions, one for a matrix-matched, trace metal-doped standard reference blood, and the second to contain and confine the unknown blood sample. Alternatively, a matrix-matched, trace metal-doped reference blood could be inserted into the
30 analytical run such that each unknown had a standard immediately adjacent to it. This would lead to 33% reference samples in the analytical run as opposed to 50% if standard and unknown were applied to the same collection device.

The results from this Experiment are presented in Appendix Experiment 18. This experiment examined heat and air-dried blood partially absorbed into an aluminium
35 hydroxide:cellulose powder matrix, and matrix-free blood dried onto an impervious Perspex substrate.

If the corrected and normalised "no-matrix" blood is examined, the numbers are reproducible. Indeed, values are commonly comparable to the dried material. In the 'no matrix' blood, both mercury and lead are recorded and the reproducibility of lead is with a precision of 14%. Good numbers are also recorded for uranium on the dried material, but in the blood matrix alone, the numbers are considered to be 'below detection limit', consistent with a matrix uranium background and anticipated absence in the blood.

Example 6: Wear Metal Analysis In Oils

Experiment 13

The objective of this experiment was to carry out pilot analysis of wear metals in engine oil. It is held that the technology being investigated is equally applicable to the analysis of wear metals in oils, and that wear metals analysis is a major global industry aimed at early detection and prevention of catastrophic plant failure. Such early detection is of particular importance to the military, airline, shipping and mining industries where component failure (automotive, heavy machinery, weaponry and the like) may lead to tragic loss of life and destruction of expensive plant.

Oil from the engine of a 'new' Ford Fairlane was sampled hot, with the engine still running, via the dip-stick. Oil from a single dip of the dip-stick was transferred to both an unwashed and washed 3:1 $\text{Al}(\text{OH})_3$:cellulose powder matrix pellet pressed at 500kg/sq in. Duplicate pellets (without oil) were prepared as blanks and all four pellets analysed by UV LA-ICP-MS. Instrument settings as for Experiment 5 were used, with minor adjustments for day-to-day variations. The results of analysis are presented in Appendix Experiment 13.

When blank corrected, there is very little difference between results obtained on the unwashed and washed matrices. If the two matrices are treated as a single matrix, then precisions, with the exception of Iron, are excellent, commonly being <1 for the restricted range of analytes expected in oil. Reproducibility of the data, are thus excellent and this is graphically illustrated in the X-Y log plot of 'concentration' versus elements comprising Chart Experiment 13/1. Here, consistent with the precision/reproducibility data, Iron excepted, the two profiles are effectively superimposed upon each other.

The experiment clearly indicates the general reproducibility of the analysis and indicates considerable promise for the technique.

Experiment 15

This experiment had as its main objective, the analysis of oil from the engines of five different cars, collected under the same conditions as described above, that is hot

with the engines running, on three consecutive days, to assess whether contrasts in wear metal content in oil from cars of contrasting age, engine capacity and, presumably oil used, could be established. For one 'old' car, which required frequent oil top-ups between services, a sample of the new top-up oil was available for comparison. The oil
5 was collected as for Experiment 13, but in duplicate on unwashed 3:1 $\text{Al}(\text{OH})_3$:cellulose powder pellets pressed at 100kg/sq in pressure; new reference oil was dipped with a glass rod and applied, in duplicate, to equivalent pellets. All samples were analysed by UV LA-ICP-MS; the results of the expanded range of analytes are presented as Appendix Experiment 15.

10 During the course of the analysis, eleven glass standard measurements were made. The precisions on the raw glass data are generally in the range 10 to 20%. However, when the raw data are normalised to average cerium, precisions are generally excellent and, with the exception of selenium, cadmium and mercury, are <10; selenium and cadmium are just marginally higher and mercury sits at 24%. The cerium
15 normalised glass standard data have been plotted in a log X-Y line chart plot which comprises Chart Experiment 15/1. Here, it is clear that the several profiles essentially superimpose, consistent with the very good precisions and reproducibility. In addition to the glass standard, 10 air blank measurements were made throughout the analytical run. These have been drift corrected and the average drift corrected air blank has been
20 used to correct the reported data.

Assessment of the data clearly demonstrates significant, and often marked differences, in specific analytes between the engine oils from the different vehicles. Oil from two cars, 'John' and 'Scott', were selected to demonstrate these contrasts. 'John' engine oil is plotted as a log X-Y line chart in Chart Experiment 15/2 while 'Scott' oil
25 comprises Chart Experiment 15/3. Examination of the respective Charts illustrates that while, there is general profile superimposition for the respective replicate oil analyses, there are some clear difference in the shapes of the respective profiles as well as peak height contrasts between equivalent analytes. Chart Experiment 15/4 graphs the averaged composition of 'John' and 'Scott' oil (n=6). This latter Chart clearly
30 emphasises the marked compositional contrast between the two oils. Hence, from this experiment, it may reasonably be concluded that the technique can readily identify and measure analyte contrasts in the examined engine oils. It is clear from the pilot experiments that wear metal analysis of oils of plant in service by LA-ICP-MS techniques is feasible and useful. The experimentation into the analysis of wear metals
35 in oils indicates considerable potential economic benefits of being able to, for example, regularly monitor potential component wear, through 'dip-stick' sampling, in plant in

service, that is without the need to plant take off-line, are large. In this way plant down-time can be carefully scheduled with minimal impact upon operations.

The use of a defocused laser to ablate sample matrices is a variation of the protocols described, which can be used to improve laser coupling to the sample. If a
5 laser is focused on the surface of a sample, the first crater it produces is a response to the laser focal point being on the surface of the sample. As soon as the surface material has been ablated and removed, the next ablation event (laser shot) is into the crater area from the first shot where there is no focus and, therefore, the laser coupling is diminished. If, however, the laser is focused below the surface, that is, it is
10 defocused at the surface, potentially it is now possible to generate a more active ablation because a large amount of material can be ejected from the middle of the sample because the focussing is below the surface. Hence, it might be expected that at least the first and second shots will produce a lot of ablation debris and therefore this may increase the sensitivity because, at this stage the ablation ejecta is a
15 powder/aerosol and this may be more efficiently transported to the plasma torch. For the existing equipment, laser defocusing can be fairly readily achieved manually. Modern lasers have automatic defocus capabilities where the depth for defocusing can be simply programmed.

As a further modification of the present protocols, triple shot ablation, as compared with
20 double shot, at each point in a 10 point by 10 point raster grid, may be used.

Example 7: Quantitation using solution doped matrices (further experiments)

In this example three fibrous cellulose matrices, being Whatman 541, high purity Whatman 541 and old Whatman 540 filter papers (Whatman International Ltd, Maidstone, England), were prepared as blank material by affixing to a support substrate
25 using a backing tape; a sample of the backing tape (3M Scotch Permanent Double Stick Tape) was also analysed. The raw count data was analysed firstly as isotopic concentrations for the designated elements and secondly as elemental abundance concentrations derived from the isotopic data using natural abundance relations. All elemental data has been air blank corrected. Air blank correction has produced
30 negative values for isolated analytes implying that the analyte concentrations in the average air blank are significantly higher than in the matrices for those analytes. Examination of the data illustrates generally high analyte air blank values.

All elements have been spike corrected (ie. normalised to an average value for the spike) and 'old' refers to fibrous cellulose substrates that have previously been
35 opened and exposed to the laboratory environment through 'open' long-term storage. 'New' refers to sealed fibrous cellulose substrates opened for this experiment. With

respect to the single versus multiple layer substrate data, it appears probable that analysis of single layer substrates may have involved laser penetration into the backing tape. Hence, data for single layer substrates may reflect composite data whereas for the multiple layers, where the top layer was peeled off immediately prior to analysis, the data reflect only the cellulose matrix substrate.

The data illustrated lower concentrations for a significant number of analytes in multiple, relative to single, layer matrices; other analytes are essentially equivalent while some are higher. For many analytes, for example Cu, Zn, Sn, concentrations in the backing tape is very much greater than in the both the single and multi layer matrices but, here, the single layer matrices are much higher in these elements than the equivalent multi layer material. This strongly suggests that laser penetration to the backing tape has occurred and that much of the difference between single and multi layers has little to do with handling contamination.

Furthermore, the corresponding data for 'new' versus 'old' clearly demonstrates significantly lower overall concentrations in the new matrices, both single and multiple. This latter observation strongly suggests that long-term exposure of matrices to the laboratory environment has led to variable, but significant ambient laboratory contamination of exposed matrices.

Further experiments examined white and black Whatman 540 filter paper cellulose matrices (Whatman International Ltd, Maidstone, England) doped with 1ppm multi-element standard (details are provided in the table) and with blood.

The data have been matrix blank corrected. For many of the analytes the air blank is high and similar to the concentrations measured in the white and black cellulose blanks (matrices without samples applied).

The isotopic data, as obtained, was converted to elemental concentrations and the multi element standard and blood doped samples have effectively been doubly corrected. The respective white and black cellulose matrix blanks have first been air blank corrected using the average of two air blanks. Following this, the averaged data, for multi standard and blood doped white and black cellulose, have been corrected using the respective corrected air blank corrected white and black cellulose matrix blanks. There is good correlation between the averaged corrected values for white and black multi element standard doped matrix samples and white and black blood doped samples. Little difference exists between the multi element standard and the blood on white and black matrices. The data obtained in this experiment also illustrates excellent reproducibility for the vast majority of analyst across the mass spectrum in both multi element and blood doped matrices.

Comparison of the computed concentrations in the blood may now be compared with anticipated concentration ranges from the literature. Data for Fe, Cu Zn, Sn, Ba and Pb show very good agreement.

Hardware optimisation

5 This experiment was to evaluate hardware optimisation at low, medium and high mass, using respectively manganese, lanthanum and lead. The isotopic data (isotopic concentrations), as obtained, has been rearranged and treated in a manner analogous to that in Example 7. For the current data, air blank, 540 matrix blank, 1ppm multi element standard and blood doped matrices were examined during optimisation at the
10 relevant masses. Again, the respective 540 matrix blanks have been air blank corrected by subtracting the averaged values from the averaged matrix blank values. Using the corrected matrix blanks, both the 540 multi element and blood doped matrices have been matrix corrected. Again using the corrected data, concentrations in ppb in blood have been computed.

15 The current data appear to indicate that low mass optimisation may be preferable. When doubly corrected, the indications are that, both for the multi element and blood doped matrices, optimisation at the lower mass, that is manganese, appears preferable to the mid mass and to the high mass. Once again, it is clear, with respect to quantification of trace element in the blood, matrix matched standards are of particular
20 value.

Detection limits and precision

The experiment was designed to establish detection limits, precision and quantitation for solution doped cellulose matrices. A series of standards were used for these experiments. In addition a reagent blank was also used.

25 Deionised water samples were doped, using a 'stock' multi-element standard solution, to produce a series of aqueous multi-element standard solutions with element concentrations of 100, 200; 500; 1000; 2000; 5000 and 10000 ppb. 100 μ L of each of these aqueous standard solutions was transferred to fibrous cellulose matrix pads, prepared from Whatman 540 filter paper (Whatman International Ltd, Maldstone,
30 England), using a pipette; the pads were affixed to Perspex supports using 3M Scotch Permanent Double Stick Tape. Deionised water matrix blanks were also prepared by pipetting 100 μ L of deionised water onto the matrix pads. In addition, solutions of three Certified Reference Materials, SARM's 1, 3 and 46 (South African Bureau of Standards) were diluted 250 times, and 100 μ L aliquots of each were doped onto Whatman 540
35 matrix pads. In all, 10 matrix pads of each aqueous standard concentration and CRM were prepared along with deionised water matrix blanks. A 2ppm samarium internal

standard solution spike was added to the respective matrix pads to facilitate internal normalisation; the spike was added using a pipette. All doped matrix pads were dried at 105°C for two hours prior to ablation.

Five of each set of ten prepared matrices were analysed on successive days.

5 The sample holders, with affixed matrix pads, were placed in the laser ablation cell of a UP 266 UV Laser System connected to an X Series ICP-MS with Xi Cone System (Thermo Optek (Australia) Pty Ltd, Rydalmere, Australia) and ablated on a 10x10 matrix raster using a UV laser operating at 266 nm, 10Hz at a fluence of 6 Millijoule and an argon flow between 900 and 1000 mL per minute for 60 seconds.

10 Samples were analysed manually and results have been corrected for air blanks, facilitating cross comparison between CRM and standard matrix matched samples. The output data was acquired as raw counts from on-board software and exported into Excel and manipulated. No algorithms were used for computations. From these corrected data, Standard Deviations and Coefficients of Variation have been
15 computed as measures of reproducibility and precision. Finally, quantitative trace element compositions for the 44 analytes examined in the exemplary run were computed for the CRM's; sub-20ppb detection limits for most analytes were achieved.

Data obtained data is set out in Appendix Experiment M1. It is also quite apparent that data for the standards, when plotted, indicate excellent calibration can be
20 achieved. Quantitation of data for the CRM's indicated extremely good agreement for elemental concentrations for all elements with values (for samples once diluted) in the optimum analytical range of the technique.

There are a number of points that this data demonstrates.

- 1) It is possible to achieve sub 5% precision for a wide range of elements using the
25 analytical protocols developed in conjunction with ICP-MS.
- 2) It is possible to achieve sub 20ppb detection limits for a wide range of elements simultaneously.
- 3) It is possible to achieve accurate quantitative data, using matrix matched certified reference materials, or other equivalent CRM's.

30 Examples of useful areas of application of the methods and devices of the present invention are:

- screening occupationally exposed workers for anomalous levels of a range of toxic metals;
- monitoring environmental exposure of the general population to toxic metals;
- 35 • screening populations for trace/ultra trace element deficiencies for preventative medicine

- screening trace/ultra trace element deficiencies, and toxic heavy metal excesses, in bloodstock, general livestock, zoo animals (including animals in endangered species breeding programs), and domestic pets for veterinary medicine; and monitoring heavy metal pollutants in slaughter animals for meat product quality control in the human food chain.
- Monitoring/detecting wear of mechanical components of plant, machinery and the like by analysing lubricating oils.

Although the invention has been described with reference to certain preferred embodiments, variations in keeping with the broad principles and the spirit of the invention are also contemplated as being within its scope.

APPENDIX EXPERIMENT 2

Element - ppb* in original	Li	Be	Al	Ti	V	Cr	Mn	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Rb
TiO ₂ /HCl -001 leachate	7	<1	8,340	174,555	<1	<1	438	<1	<1	457	364	8	<1	<1	<1	76
TiO ₂ /HNO ₃ -002 leachate	11	<1	13,780	76,451	<1	14	638	<1	<1	527	438	13	1	<1	<1	108
Al(OH) ₃ /HCl -003 leachate	37	4	41,530	180	<1	118	48	<1	<1	14	<1	2,357	<1	<1	<1	<1
Al(OH) ₃ /HNO ₃ -004 leachate	45	4	48,312	1,456	<1	17	33	<1	<1	50	<1	2,523	<1	<1	<1	5
Pig Toe A digest	63	<1	11,800	1,779	<1	<1	761,988	<1	<1	113	817	<1	<1	<1	<1	23
Pig Toe B digest	84	<1	9,956	2,086	<1	<1	475,385	<1	<1	138	890	<1	<1	<1	<1	43
Pig Toe C digest	109	<1	10,314	2,165	<1	<1	760,369	<1	<1	126	922	<1	<1	<1	<1	72
Pig Toe D digest	67	<1	9,424	1,922	<1	<1	936,818	<1	<1	170	421	<1	<1	<1	<1	59
Guarodun E solute	8	<1	2,378	91	<1	359	265	<1	107	18	149	<1	<1	<1	<1	20
Guarodun F solute	4	1	2,218	92	<1	327	208	<1	103	29	181	<1	<1	<1	<1	31
Glucose G solute	9	2	1,886	89	<1	345	96	<1	110	21	131	<1	<1	<1	<1	18
Cellulose H digest	9	7	22,353	1,391	50	798	298	<1	953	523	982	<1	<1	<1	<1	62
HBM Cellulose I digest	71	3	25,313	1,278	50	2,382	1,538	<1	1,282	1,671	1,413	<1	<1	<1	<1	78
* ppb in solution for leachates																

APPENDIX EXPERIMENT 2

Element - ppb* In original	Sr	Y	Zr	Nb	Mo	Ag	Cd	Sn	Sb	Te	Cs	Ba	La	Ce	Pr	Nd
TiO2/HCl -001 leachate	134	<1	62	<1	69	<1	<1	<1	<1	<1	<1	2,808	8	8	<1	<1
TiO2/HNO3 -002 leachate	185	1	160	<1	<1	<1	<1	<1	<1	<1	<1	3,250	8	11	<1	<1
Al(OH)3/HCl -003 leachate	170	<1	1,289	<1	<1	<1	<1	168	<1	<1	<1	<1	<1	2	<1	<1
Al(OH)3/HNO3 -004 leachate	189	<1	818	<1	<1	<1	<1	174	<1	<1	<1	<1	<1	3	<1	<1
Pig Toe A digest	237,704	<1	<1	<1	10	<1	<1	<1	<1	<1	<1	66,117	4	9	<1	<1
Pig Toe B digest	233,800	<1	1	<1	34	<1	<1	<1	<1	<1	<1	40,257	4	15	<1	<1
Pig Toe C digest	332,025	<1	<1	<1	41	<1	<1	<1	<1	<1	<1	85,251	8	16	<1	<1
Pig Toe D digest	303,588	<1	<1	<1	61	<1	<1	<1	<1	<1	<1	101,341	10	28	<1	<1
Glucodin E solute	188	<1	<1	7	63	<1	<1	<1	<1	<1	<1	72	1	2	<1	<1
Glucodin F solute	229	<1	<1	6	61	<1	<1	<1	<1	<1	1	43	<1	<1	<1	<1
Glucose G solute	22	<1	<1	<1	12	<1	<1	<1	<1	<1	<1	8	<1	<1	<1	<1
CelMase H digest	357	<1	806	217	870	<1	<1	668	<1	<1	<1	166	8	12	<1	<1
NEB Cellulase I digest	13,800	<1	1,351	582	524	<1	<1	557	<1	<1	<1	480	6	11	<1	<1
* ppb in solution for leachates																

APPENDIX EXPERIMENT 2

Element - ppb* in original	Eu	Sm	Gd	Tb	Dy	Ho	Er	Tm	Yb	Lu	Hf	Ta	W	Hg	Tl	Pb	Bi	Th	U
TI02/HCl -001 leachate	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	19,014	<1	4	10
TI02/HNO3 -002 leachate	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	20,394	<1	3	<1
Al(OH)3/HCl -003 leachate	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	134	<1	<1	<1	1	<1	<1	3	135
Al(OH)3/HNO3 -004 leachate	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	131	<1	<1	<1	<1	<1	<1	2	152
Pig Toe A digest	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
Pig Toe B digest	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
Pig Toe C digest	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
Pig Toe D digest	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1
Glucodin E solids	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	1	<1	<1	<1	<1	<1	5	1	<1
Glucodin F solids	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	5	<1	<1
Glucose G solids	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	41	<1	<1	<1
Cellulose H digest	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	24	186	137	55	<1
HEBN Cellulose I digest	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	<1	25	<1	<1	32	<1
* ppb in solution for leachates																			

APPENDIX EXPERIMENT 3

Sample	Sample No.	Pelletize	Absorption Rate of SY-2	Dissolution	Comments
Glucose	1	POOR	Fast	Yes	Pellet dissolved, absorbed quickly
Cellulose	2	OK	10-15 sec	No	Solution absorbed slowly
AR Starch	3	OK	Slow	Partial	Pellet swells
UR Starch	4	OK	Slow	No	Pellet swells
Glucose + Cellulose 1:1	5	OK	Slow	Partial	Absorption OK, partial dissolution, holes on surface
Glucose + Cellulose 3:1	6	OK	Slow	Partial	Dissolution of pellet
Cellulose + Glucose 3:1	7	OK	V. Slow	Partial	Partial dissolution of pellet, holes left on surface
Glucose + AR Starch 1:1	8	OK	V. Slow	Partial	Dissolution and swelling
Glucose + UR Starch 1:1	9	OK	V. Slow	Partial	Dissolution and swelling
Cellulose + AR Starch 1:1	10	OK	Slow	No	Dissolution and swelling
Cellulose + AR Starch 3:1	11	OK	Slow	No	Dissolution and swelling
AR Starch + Cellulose 3:1	12	OK	Slow	No	Swelling of surface
Cellulose + UR Starch 1:1	13	OK	Slow	No	Swelling of surface
Cellulose + UR Starch 3:1	14	OK	Slow	No	Swelling of surface
UR Starch + Cellulose 3:1	15	OK	Slow	No	Swelling of surface
Glucose + Cellulose + AR Starch 1:1:1	16	OK	V. Slow	Partial	Dissolution and swelling
Glucose + Cellulose + UR Starch 1:1:1	17	OK	Slow	Partial	Dissolution and swelling

APPENDIX EXPERIMENT 5A

Isotope - Raw Counts	Mg 24	Ca 44	Mn 55	Fe 56	Cu 65	Zn 66	As 75	Se 77	Mo 98	Ba 138	Pb 208
WET											
"02/11/07 CELLULOSE AIRBL1"	38,010	14,080	2,719	25,180	2,686	377	660	432	138	111	73
"02/11/07 CELLULOSE AIRBL2"	35,740	13,480	2,578	24,210	2,592	309	626	443	108	38	58
"02/11/07 CELLULOSE BLANK1"	60,150	24,560	7,263	689,700	15,140	8,261	671	328	1,542	5,132	8,886
"02/11/07 CELLULOSE BLANK2"	58,520	20,620	10,250	701,400	10,720	5,452	704	393	2,254	3,988	6,359
"02/11/07 CELLULOSE ST21"	75,080	31,360	24,930	375,200	2,948	1,459	849	400	2,095	7,150	8,334
"02/11/07 CELLULOSE ST22"	73,650	28,060	22,240	337,700	3,588	1,065	714	426	1,663	5,975	5,185
"02/11/07 CELLULOSE BLOOD1"	128,300	28,240	4,941	2,803,000	6,377	15,490	688	447	735	3,213	10,030
"02/11/07 CELLULOSE BLOOD2"	101,900	26,030	5,736	2,218,000	6,518	7,604	714	468	817	4,711	2,713
"02/11/07 CELLULOSE GLSSTD1"	233,300	644,400	175,200	227,800	50,490	52,420	25,230	918	91,410	245,700	37,880
"02/11/07 CELLULOSE AIRBL3"	33,660	12,570	2,553	27,070	2,638	339	747	462	145	46	73
"02/11/07 CELLULOSE AIRBL4"	35,000	12,880	2,545	28,020	2,765	352	788	511	148	42	65
DRY											
"02/11/07 CELLULOSE AIRBL5"	25,660	10,520	2,391	23,630	2,197	327	860	511	145	95	74
"02/11/07 CELLULOSE AIRBL6"	25,490	10,700	2,465	24,380	2,211	336	831	532	128	41	73
"02/11/07 CELLULOSE BLANK5"	35,730	18,150	4,002	71,500	2,491	5,882	813	379	384	2,751	2,758
"02/11/07 CELLULOSE BLANK6"	38,820	19,450	4,194	76,720	2,500	5,450	882	398	348	2,147	2,319
"02/11/07 CELLULOSE ST23"	102,100	30,740	36,790	678,500	3,000	6,896	865	395	2,332	11,880	7,340
"02/11/07 CELLULOSE ST24"	117,400	35,750	43,580	791,600	3,104	5,782	948	465	2,869	14,010	8,050
"02/11/07 CELLULOSE BLOOD3"	107,400	32,000	4,320	2,898,000	6,533	8,471	829	539	382	1,066	3,126
"02/11/07 CELLULOSE BLOOD4"	108,200	33,000	4,300	2,785,000	6,308	7,468	957	540	392	1,179	3,369
"02/11/07 CELLULOSE GLSSTD7"	145,100	571,300	188,600	212,500	41,860	35,320	25,530	927	102,000	288,800	61,500
"02/11/07 CELLULOSE AIRBL7"	28,040	12,350	2,988	30,210	2,224	350	962	505	172	39	79
"02/11/07 CELLULOSE AIRBL8"	28,620	12,380	2,962	30,540	2,255	364	971	566	162	33	70
Ave ST2	71,975	14,940	36,137	660,940	557	673	59	62	2,246	10,496	5,167
Ave Blood	69,025	14,195	257	2,757,880	3,925	2,303	96	172	37	-1,332	709
Blank corrected											
"02/11/07 CELLULOSE ST23"	64,325	12,435	32,737	604,390	505	1,230	17	27	1,977	9,431	4,802
"02/11/07 CELLULOSE ST24"	79,625	17,445	39,537	717,490	609	116	100	97	2,514	11,561	5,512
% Std Dev	15	24	13	12	13	117	100	78	17	14	10
"02/11/07 CELLULOSE BLOOD5"	69,625	13,895	267	2,823,880	4,038	2,805	81	171	37	-1,393	588
"02/11/07 CELLULOSE BLOOD6"	69,425	14,885	247	2,697,890	3,813	1,800	110	173	37	-1,270	831
% Std Dev	1	5	6	3	4	31	21	1	0	-7	24

Experiment 5A 1/1

APPENDIX EXPERIMENT 5B

Isotope - Raw Counts	Mg 24	Ca 44	Mn 55	Fe 58	Cu 66	Zn 66	As 76	Se 77	Mo 98	Ba 138	Pb 208
¹⁰² Hf107 CELLULOSE AIRBL5*	25,660	10,520	2,391	21,630	2,197	327	860	511	145	95	74
¹⁰² Hf107 CELLULOSE AIRBL6*	28,490	10,700	2,465	24,380	2,211	338	831	511	128	41	73
¹⁰² Hf107 CELLULOSE AIRBL5*	25,660	10,520	2,391	21,630	2,197	327	860	511	145	95	74
¹⁰² Hf107 CELLULOSE AIRBL6*	26,460	10,700	2,465	24,380	2,211	338	831	511	128	41	73
¹⁰² Hf107 CELLULOSE BLANK5*	35,730	16,160	4,002	71,500	2,491	5,862	813	378	384	2,751	2,758
¹⁰² Hf107 CELLULOSE BLANK6*	39,820	18,460	4,104	76,720	2,500	5,450	882	358	346	2,147	2,318
¹⁰² Hf107 CELLULOSE SY23*	102,100	30,740	36,790	678,500	3,080	6,886	865	396	2,332	11,880	7,340
¹⁰² Hf107 CELLULOSE SY24*	117,400	35,750	43,580	781,600	3,104	5,782	948	465	2,869	14,010	8,050
¹⁰² Hf107 CELLULOSE BLOOD3*	107,400	32,000	4,320	2,889,000	8,533	8,471	928	539	382	1,058	3,126
¹⁰² Hf107 CELLULOSE BLOOD4*	109,200	33,000	4,300	2,768,000	6,306	7,466	967	540	392	1,179	3,358
¹⁰² Hf107 CELLULOSE GLSS100*	145,100	571,300	146,800	212,500	41,640	35,320	25,530	827	102,000	288,800	81,500
¹⁰² Hf107 CELLULOSE AIRBL7*	28,040	12,350	2,966	30,210	2,224	350	862	505	172	39	78
¹⁰² Hf107 CELLULOSE AIRBL6*	28,620	12,380	2,982	30,540	2,255	364	871	555	182	33	70
Blank Corrected											
¹⁰² Hf107 CELLULOSE SY23*	64,325	12,435	32,731	604,390	505	1,230	17	27	1,977	0,431	4,802
¹⁰² Hf107 CELLULOSE SY24*	79,625	17,445	39,537	717,490	608	116	100	97	2,514	11,581	5,512
¹⁰² Hf107 CELLULOSE BLOOD3*	68,625	13,695	267	2,823,890	4,038	2,805	81	171	37	-1,393	588
¹⁰² Hf107 CELLULOSE BLOOD4*	68,425	14,685	247	2,891,890	3,813	1,800	110	173	37	-1,270	831
Conc in ppm in SY-2	2,69	7,86	0,32	2,43	5,20	248,00	17,30	20,00	0,53	460,00	85,00
(MgO)		(CaO)	(MnO)	(Fe2O3+FeO)				conc ratio			
% in sample								for SY-2		197,07	
% Metal in SY-2	0,60	0,71	0,77	0,70							
				0,78							
Conc in ppm in SY-2	16220	56857	2478	17010	5,20	248,00	17,30	20,00	0,53	460,00	85,00
				27688							
Conc in ppm for SY-2 in 50mL sample	82,31	280,51	12,58	86,31	0,03	1,28	0,09	0,10	0,00	2,33	0,43
				140,50							
Average counts for SY-2	71875	14840	38137	660940	557	673	59	82	2246	10498	5157
Conc in ppm for blood samples (exp)	78,9	274	0,089	360	0,166	4,31	0,143	0,280	<0,001	<0,001	0,059
Expected concentrations for blood values where found in literature	50,0	300		500-1800	08-18	6,00					0,06

Experiment 5B/1

APPENDIX EXPERIMENT 12

Isotope - Raw Counts	Li 7	Mg 24	Ca 44	V 51	Cr 52	Mn 55	Fe 56	Cu 65	Zn 66	Ga 69	As 75	Sr 88	Zr 90	Mo 98	Cd 114
"0201127 HHG GLS STD 1"	107,400	194,900	690,900	182,200	152,300	252,900	258,100	41,720	25,830	180,900	25,180	415,400	177,500	112,700	38,070
"0201127 HHG GLS STD 2"	105,400	187,600	634,200	180,100	149,000	245,500	244,400	41,450	28,150	180,000	25,580	403,100	177,400	112,900	38,610
"0201127 HHG AIR BL 1"	1,919	94,140	21,220	122	1,688	10,620	50,128	1,434	1,246	231	3,955	1,761	139	252	186
"0201127 HHG AIR BL 2"	2,614	108,100	23,090	165	1,759	3,167	50,620	1,495	1,428	254	3,671	1,182	84	232	214
"0201127 HHG CELL ON BL 1"	2,024	101,800	27,540	235	6,288	3,582	61,890	1,602	1,984	445	2,785	1,181	241	341	4,647
"0201127 HHG CELL ON BL 2"	2,032	107,500	28,350	205	6,311	3,586	62,860	1,556	1,868	708	2,768	1,057	210	333	1,924
"0201127 HHG CELL R BL 1"	1,986	82,690	24,830	233	5,007	2,837	54,740	1,353	1,381	235	3,267	1,026	97	289	4,55
"0201127 HHG CELL R BL 2"	1,978	101,400	26,040	159	6,008	3,250	60,640	1,444	1,491	387	3,480	987	104	308	528
"0201127 HHG CELL UW BL 1"	2,213	118,000	37,410	1,080	7,191	4,622	78,450	1,822	1,778	568	3,571	1,495	127	325	1,183
"0201127 HHG CELL UW BL 2"	2,361	142,800	38,910	217	7,387	3,651	87,650	1,453	1,958	703	3,674	1,345	177	343	1,878
"0201127 HHG CELL ON ME 1"	5,211	122,380	22,180	2,751	8,867	5,811	92,480	1,988	2,422	2,653	3,013	7,890	3,179	2,007	1,580
"0201127 HHG CELL ON ME 2"	4,343	127,000	33,410	4,217	10,980	18,040	77,020	2,631	2,310	4,486	3,968	8,900	5,203	2,582	1,405
"0201127 HHG CELL R ME 1"	4,963	122,000	28,700	4,724	10,510	9,195	75,280	2,882	3,457	4,268	3,691	10,630	5,188	3,552	2,382
"0201127 HHG CELL R ME 2"	4,965	130,600	30,680	5,087	11,290	10,120	88,430	2,788	3,623	4,783	4,319	13,340	5,819	3,817	2,714
"0201127 HHG CELL UW ME 1"	2,850	124,200	23,210	2,241	6,764	14,920	57,330	1,980	5,443	2,185	2,935	7,087	2,364	1,891	4,907
"0201127 HHG CELL UW ME 2"	3,703	191,200	38,780	3,760	10,320	13,670	73,610	4,235	6,735	3,644	4,100	8,288	4,262	2,345	4,865
"0201127 HHG GLS STD 3"	86,400	188,700	684,900	184,900	137,500	222,400	235,800	34,300	21,930	162,200	22,170	383,800	180,200	89,780	30,370
"0201127 HHG GLS STD 4"	82,880	188,000	646,500	177,800	147,800	243,100	237,900	39,890	26,360	192,200	25,620	442,700	192,600	114,900	39,260
"0201127 HHG AIR BL 3"	2,128	120,200	28,330	162	2,625	3,701	57,110	1,508	1,804	306	4,043	1,135	109	375	280
"0201127 HHG AIR BL 4"	2,051	123,100	24,870	184	3,245	3,691	57,590	1,508	1,748	302	3,952	6,448	98	376	238
Blank corrected															
"0201127 HHG CELL ON ME 1"	1,188	17,640	5,755	2,531	301	2,032	9,945	410	508	2,083	237	6,581	2,968	1,700	-1,588
"0201127 HHG CELL ON ME 2"	2,315	17,350	5,465	3,867	4,860	12,461	14,685	1,053	474	3,921	1,210	8,791	4,982	2,225	-1,831
"0201127 HHG CELL R ME 1"	3,077	26,458	3,255	4,528	4,803	6,167	17,570	1,454	2,001	3,583	529	8,623	3,354	1,890	1,890
"0201127 HHG CELL R ME 2"	3,019	30,955	4,635	4,881	5,673	7,882	11,740	1,390	2,487	4,472	951	12,333	5,719	3,519	2,222
"0201127 HHG CELL UW ME 1"	528	-6,700	-12,450	1,588	-1,535	10,804	-95,420	488	3,505	1,558	-688	5,645	2,251	1,357	3,372
"0201127 HHG CELL UW ME 2"	1,401	300	-1,880	3,107	3,031	9,794	60	2,783	4,877	3,007	468	8,897	4,179	2,011	3,330
Normalized to carian															
"0201127 HHG CELL ON ME 1"	1,188	17,640	5,755	2,531	301	2,032	9,945	410	508	2,083	237	6,581	2,968	1,700	-1,588
"0201127 HHG CELL ON ME 2"	1,880	10,944	3,447	2,621	2,959	7,880	9,269	604	288	2,473	783	5,545	3,148	1,404	-1,188
"0201127 HHG CELL R ME 1"	1,093	18,343	2,011	2,767	2,967	3,808	10,854	868	1,256	2,482	108	6,943	3,143	2,072	1,188
"0201127 HHG CELL R ME 2"	1,722	17,444	2,644	2,700	3,178	4,045	8,897	783	1,419	2,551	542	7,033	3,262	2,009	1,268
"0201127 HHG CELL UW ME 1"	760	-8,880	-16,501	2,104	-2,834	14,356	-21,763	848	4,011	2,063	-885	7,482	2,803	1,789	4,468
"0201127 HHG CELL UW ME 2"	1,062	227	-1,425	2,355	2,298	7,418	45	2,085	3,698	2,280	377	5,207	3,168	1,525	2,504
Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Cu	Zn	Ga	As	Sr	Zr	Mo	Cd
"0201127 HHG CELL ON ME 1"	1,279	22,328	-278,883	2,539	358	2,032	-10,738	1,330	2,100	3,465	237	7,887	3,775	7,065	-5,588
"0201127 HHG CELL ON ME 2"	1,678	13,853	185,727	2,529	3,508	7,880	10,188	2,155	1,072	4,115	763	6,713	6,127	5,624	-4,133
"0201127 HHG CELL R ME 1"	2,122	20,688	98,675	2,808	3,540	3,809	11,857	2,815	4,431	4,958	189	7,197	8,115	8,568	4,088
"0201127 HHG CELL R ME 2"	1,882	21,701	127,107	2,798	3,793	4,045	7,303	2,573	5,085	4,244	542	8,517	8,347	8,329	4,417
"0201127 HHG CELL UW ME 1"	757	-11,240	-788,589	2,111	-2,428	14,356	-21,732	2,089	17,080	3,437	-685	8,058	5,808	7,464	15,570
"0201127 HHG CELL UW ME 2"	1,148	288	-68,530	2,383	2,742	7,418	50	8,890	13,254	3,784	377	6,303	6,164	6,327	8,788
"0201127 HHG CELL ON ME 3"	1,278	22,328	-278,883	2,539	358	2,032	-10,738	1,330	2,100	3,465	237	7,887	3,775	7,065	-5,588

APPENDIX EXPERIMENT 12

Isotope - Raw Counts	Sn 120	Ba 138	La 139	Ca 140	Eu 151	Dy 162	Yb 174	Hf 178	Pb 206	U 238
¹⁰² 1127 HGH GLS STD 1"	182,100	399,900	450,200	517,100	270,700	112,100	128,100	91,780	84,590	115,200
¹⁰² 1127 HGH GLS STD 2"	188,400	398,000	438,100	507,500	283,000	109,500	123,400	88,580	65,130	119,100
¹⁰² 1127 HGH AIR BL 1"	141	1,144	36	13	18	13	9	4	312	21
¹⁰² 1127 HGH AIR BL 2"	152	183	25	20	20	9	14	14	26	8
¹⁰² 1127 HGH CELL ON BL 1"	675	1,180	182	164	112	45	53	32	4,450	38
¹⁰² 1127 HGH CELL ON BL 2"	565	1,673	142	138	52	21	23	24	4,789	83
¹⁰² 1127 HGH CELL R BL 1"	528	242	52	30	64	17	12	10	889	24
¹⁰² 1127 HGH CELL R BL 2"	508	264	44	28	38	14	11	33	771	18
¹⁰² 1127 HGH CELL UW BL 1"	355	635	58	83	50	24	29	14	2,580	45
¹⁰² 1127 HGH CELL UW BL 2"	474	947	63	119	103	22	14	10	2,788	167
¹⁰² 1127 HGH CELL ON ME 1"	3,088	6,283	7,892	7,442	4,328	1,932	2,202	1,768	4,944	1,805
¹⁰² 1127 HGH CELL ON ME 2"	4,887	8,724	12,560	11,710	6,768	3,288	3,531	2,948	5,081	2,346
¹⁰² 1127 HGH CELL R ME 1"	5,747	10,898	12,480	11,830	6,827	3,112	3,407	2,525	6,512	2,378
¹⁰² 1127 HGH CELL R ME 2"	6,991	11,829	13,830	12,810	7,819	3,587	3,857	2,859	6,501	2,718
¹⁰² 1127 HGH CELL UW ME 1"	3,495	5,403	5,998	5,622	3,238	1,482	1,577	1,167	9,940	1,200
¹⁰² 1127 HGH CELL UW ME 2"	5,174	10,089	9,993	8,717	5,474	2,845	2,812	2,111	7,653	1,833
¹⁰² 1127 HGH GLS STD 3"	180,000	374,500	487,100	473,100	258,400	105,700	118,500	85,230	47,700	98,150
¹⁰² 1127 HGH GLS STD 4"	203,100	433,000	497,200	537,800	294,800	120,200	138,200	100,200	64,190	123,100
¹⁰² 1127 HGH AIR BL 3"	716	287	41	22	34	18	9	10	44	9
¹⁰² 1127 HGH AIR BL 4"	736	463	88	17	32	13	10	12	833	9
Blank connected										
¹⁰² 1127 HGH CELL ON ME 1"	2,488	4,877	7,830	7,291	4,244	1,919	2,164	1,880	340	1,546
¹⁰² 1127 HGH CELL ON ME 2"	4,277	8,399	12,398	11,559	6,708	3,238	3,493	2,818	457	2,267
¹⁰² 1127 HGH CELL R ME 1"	5,228	10,797	12,432	11,802	6,776	3,087	3,365	2,503	5,882	2,366
¹⁰² 1127 HGH CELL R ME 2"	6,073	11,567	13,882	12,782	7,583	3,562	3,878	2,837	5,881	2,686
¹⁰² 1127 HGH CELL UW ME 1"	3,081	4,912	5,981	5,301	3,152	1,483	1,568	1,155	7,168	1,094
¹⁰² 1127 HGH CELL UW ME 2"	4,760	9,899	9,940	8,918	5,368	2,622	2,781	2,088	4,979	1,727
Normalized to cadmium										
¹⁰² 1127 HGH CELL ON ME 1"	2,488	4,877	7,830	7,291	4,244	1,919	2,164	1,880	340	1,546
¹⁰² 1127 HGH CELL ON ME 2"	2,688	5,240	7,820	7,291	4,230	2,041	2,203	1,851	288	1,442
¹⁰² 1127 HGH CELL R ME 1"	3,230	6,533	7,880	7,291	4,188	1,913	2,088	1,547	3,518	1,455
¹⁰² 1127 HGH CELL R ME 2"	3,882	6,586	7,818	7,291	4,317	2,028	2,211	1,818	3,240	1,538
¹⁰² 1127 HGH CELL UW ME 1"	4,083	6,113	7,874	7,291	4,177	1,947	2,082	1,531	3,497	1,450
¹⁰² 1127 HGH CELL UW ME 2"	3,809	7,354	7,538	7,291	4,070	1,988	2,118	1,582	3,889	1,509
Element - Raw Counts	Sn	Ba	La	Ca	Eu	Dy	Yb	Hf	Pb	U
¹⁰² 1127 HGH CELL ON ME 1"	7,572	6,801	7,838	8,238	8,378	7,525	8,804	6,154	848	1,566
¹⁰² 1127 HGH CELL ON ME 2"	8,278	7,598	7,828	8,238	8,850	8,004	8,828	6,049	550	1,482
¹⁰² 1127 HGH CELL R ME 1"	9,808	9,251	7,888	8,238	9,037	7,502	8,588	5,885	6,711	1,488
¹⁰² 1127 HGH CELL R ME 2"	11,328	9,202	7,868	8,238	9,031	7,944	8,962	5,828	6,184	1,548
¹⁰² 1127 HGH CELL UW ME 1"	12,624	8,828	7,881	8,238	8,738	7,834	8,484	5,688	18,124	1,460
¹⁰² 1127 HGH CELL UW ME 2"	11,070	10,257	7,544	8,238	8,514	7,786	8,654	5,830	7,069	1,319
¹⁰² 1127 HGH CELL ON ME 1"	7,572	6,801	7,838	8,238	8,378	7,525	8,804	6,154	848	1,566

APPENDIX EXPERIMENT 12

Isotope - Raw Counts	Li 7	Mg 24	Cs 44	V 51	Cr 52	Mn 55	Fe 56	Cu 65	Zn 66	Ga 69	As 76	Sr 86	Zr 90	Mo 98	Cd 114
T0211Z7 HKH CELL ON ME 2"	1,579	13,653	165,727	2,529	3,508	7,860	10,108	2,168	1,072	4,115	763	6,713	6,127	5,824	-4,133
Std dev	212	5,094	312,851	7	2,228	4,121	14,739	584	727	459	372	887	248	870	1,019
% Std dev	15	33	-644	0	115	83	-4,693	34	46	12	76	12	4	14	-21
T0211Z7 HKH CELL R ME 1"	2,122	20,688	96,875	2,808	3,540	3,808	11,637	2,815	4,431	4,080	188	7,187	6,115	8,598	4,089
T0211Z7 HKH CELL R ME 2"	1,862	21,701	127,107	2,788	3,783	4,045	7,303	2,573	5,085	4,244	642	8,517	6,347	8,229	4,417
Std dev	894	716	21,519	5	179	187	3,208	242	482	105	242	833	164	180	248
% Std dev	9	3	19	0	5	4	34	9	10	3	85	12	3	2	6
T0211Z7 HKH CELL UWM ME 1"	767	-11,240	-789,309	2,111	-2,428	14,358	-33,792	2,088	17,090	3,437	-885	9,658	5,803	7,484	15,870
T0211Z7 HKH CELL UWM ME 2"	1,148	288	-68,930	2,383	2,742	7,418	50	6,880	13,254	3,784	377	8,303	8,184	8,327	8,786
Std dev	277	8,162	512,498	178	3,656	4,908	18,818	3,325	2,670	283	882	1,948	255	804	4,780
% Std dev	29	-168	-119	8	2,324	45	-142	76	18	7	-362	25	4	12	39

Experiment 12/3

APPENDIX EXPERIMENT 12

Isotope - Raw Counts	Sn 120	Ba 138	La 139	Ce 140	Eu 154	Dy 162	Yb 174	Hf 178	Pb 208	U 238
¹⁰² 11/27 HRH CELL ON ME 2'	8,276	7,308	7,828	8,238	8,680	8,004	8,028	9,049	550	1,482
Std dev	488	368	7	0	21	339	88	74	70	74
% Std dev.	6	5	0	0	0	4	1	1	12	5
¹⁰² 11/27 HRH CELL R ME 1'	9,909	9,251	7,638	8,238	8,767	7,502	6,586	5,885	8,711	1,488
¹⁰² 11/27 HRH CELL R ME 2'	11,328	9,202	7,828	8,238	9,031	7,944	6,852	5,028	6,184	1,548
Std dev	1,021	35	189	0	194	513	251	188	372	59
% Std dev.	9	0	2	0	2	4	4	3	6	4
¹⁰² 11/27 HRH CELL UW ME 1'	12,624	8,526	7,881	8,238	8,738	7,634	8,484	5,809	18,124	1,480
¹⁰² 11/27 HRH CELL UW ME 2'	11,070	10,257	7,544	8,238	8,514	7,788	6,854	5,830	7,859	1,519
Std dev	1,009	1,224	239	0	159	114	120	157	7,824	100
% Std dev.	9	13	3	0	2	1	2	3	82	7

Experiment 16A/1

"UNWASHED" MATRICES													
AR and RHF Beke													
Element - Row Counts													
LI	Hg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn	As	Se	Sr	Zr
Glass Standard													
180,379	178,035	54,275,269	282,339	275,236	352,091	373,770	400,083	202,827	157,619	34,725	28,845	848,428	423,431
213,282	186,388	55,148,258	288,275	283,518	380,856	350,116	408,883	221,517	131,888	36,200	28,882	688,743	440,172
Air Blank													
5,181	25,573	1,614,554	142	3,237	3,987	40,942	258,877	10,401	1,804	1,528	21,405	651	387
6,571	28,489	1,730,616	147	3,348	4,248	42,809	271,177	10,806	1,824	1,583	22,334	739	383
UW Blank													
8,410	65,725	1,901,033	689	9,175	5,856	252,296	283,005	14,520	5,483	1,987	23,820	3,593	9,488
8,882	67,228	1,920,185	725	7,906	5,876	241,588	284,748	14,208	5,764	1,915	23,887	3,245	10,630
8,743	71,788	1,907,465	677	8,108	5,882	274,062	284,882	13,898	5,828	1,910	24,723	2,946	10,694
4	5	1	4	10	2	8	1	3	2	1	2	10	7
% Std Dev													
UW AR W Blank													
5,206	34,313	1,513,148	1,273	8,851	4,898	226,219	258,142	10,877	3,317	1,519	19,707	3,838	10,438
5,477	38,142	1,589,014	1,383	7,821	5,481	249,889	275,535	11,818	3,242	1,684	19,477	3,679	11,138
5,191	34,517	1,554,742	1,340	7,288	5,448	240,561	282,780	11,828	3,445	1,873	19,588	3,678	12,978
3	6	2	4	7	5	5	5	4	3	2	1	2	7
% Std Dev													
UW NHAF W Blank													
4,803	36,102	1,557,277	1,037	7,811	5,174	141,978	263,891	10,428	3,580	1,754	20,657	3,485	9,284
4,873	38,117	1,583,758	1,123	7,288	5,820	151,987	288,086	11,956	3,886	1,899	20,088	3,828	10,510
4,881	38,089	1,541,887	1,141	7,817	5,791	157,858	278,838	10,110	3,448	1,861	19,109	3,411	9,505
2	3	1	2	4	6	5	3	5	6	6	4	3	7
% Std Dev													
UW ME 100m													
7,581	72,125	1,911,289	9,278	14,004	11,730	255,463	298,047	15,623	6,448	3,865	25,438	19,822	20,543
7,361	77,252	1,945,540	9,959	14,315	12,580	266,171	305,178	16,114	6,203	3,201	25,868	21,572	23,119
7,588	78,018	1,930,141	5,047	14,589	11,290	324,656	294,918	15,455	6,483	4,163	28,758	17,828	18,321
1	4	1	6	2	6	13	2	2	2	14	3	9	14
% Std Dev													
UW AR W ME 100m													
5,988	38,877	1,612,729	4,321	8,712	8,003	310,424	288,059	12,882	4,776	2,840	20,254	11,389	19,219
5,757	40,589	1,623,377	4,503	8,887	8,084	285,127	285,876	12,336	4,581	2,864	20,038	11,568	18,488
5,818	41,374	1,608,659	4,047	8,887	8,185	285,239	288,857	11,948	4,228	2,742	20,428	10,982	18,088
2	3	1	5	1	1	5	1	2	6	2	1	3	3
% Std Dev													
UW NHAF W ME 100m													
5,522	48,165	1,617,840	2,448	8,874	6,383	189,248	278,801	11,803	4,348	2,528	21,987	11,748	20,788
5,772	47,201	1,623,005	2,757	8,578	6,321	177,228	278,884	11,301	4,624	2,862	21,058	11,544	19,888
5,740	47,048	1,604,225	2,815	8,882	6,151	172,531	283,830	11,456	4,478	2,435	22,457	11,355	20,151
2	1	1	6	2	2	5	1	1	3	2	3	2	2
% Std Dev													
Matrix corrected													
UW ME matrix & UW Blank													

APPENDIX EXPERIMENT 16A

"UNWASHED" MATRICES													
AR and NH4F Balco													
Element - Raw Counts	Mo	Cd	Sr	Ba	La	Ce	Eu	Gd	Yb	Hf	Hg	Pb	U
Glass Standard													
"021209 HIGH GLS STD 5"	658,331	148,498	551,101	550,080	424,698	560,157	490,267	374,167	262,733	208,884	352	56,338	51,285
"021209 HIGH GLS STD 6"	683,278	154,948	580,025	685,925	442,521	585,830	487,238	333,957	272,971	213,268	349	56,243	55,018
AR Blank													
"021209 HIGH AR BL 5"	1,533	248	634	221	60	31	52	65	31	58	580	88	8
"021209 HIGH AR BL 6"	1,733	187	596	188	50	33	39	38	44	28	824	88	8
UW Blank													
"021209 HIGH 3:1UW BL1"	2,871	1,271	7,888	8,275	335	1,223	60	43	56	727	734	3,078	372
"021209 HIGH 3:1UW BL2"	2,871	750	8,285	11,414	232	384	89	41	52	841	734	1,878	285
"021209 HIGH 3:1UW BL3"	2,845	716	8,830	15,280	221	307	45	63	42	532	779	909	218
% Std Dev	2	34	10	41	24	80	21	24	14	28	3	58	26
UW AR W Blank													
"021209 HIGH 3:1UW AR W BL1"	2,889	483	10,783	1,589	138	172	66	51	25	537	532	637	428
"021209 HIGH 3:1UW AR W BL2"	3,288	339	11,269	1,040	74	138	37	48	33	575	532	488	508
"021209 HIGH 3:1UW AR W BL3"	3,855	388	11,560	1,153	74	112	35	24	41	804	527	584	548
% Std Dev	8	16	3	23	39	21	36	36	24	4	9	13	12
UW NH4F W Blank													
"021209 HIGH 3:1UW NH4F W BL1"	2,923	334	7,888	1,834	219	121	40	19	39	577	604	453	407
"021209 HIGH 3:1UW NH4F W BL2"	3,128	307	8,341	2,171	167	140	39	34	52	589	645	542	434
"021209 HIGH 3:1UW NH4F W BL3"	3,419	287	23,091	1,889	165	223	51	58	21	608	855	478	380
% Std Dev	8	8	67	9	17	31	16	33	41	2	19	8	5
UW ME Upgm													
"021209 HIGH 3:1UW ME1"	22,888	7,178	26,544	24,883	21,149	25,204	24,639	21,187	18,427	18,288	1,088	4,458	3,335
"021209 HIGH 3:1UW ME2"	21,542	6,185	27,931	28,278	20,589	24,311	25,369	21,242	18,161	18,883	858	3,824	3,288
"021209 HIGH 3:1UW ME3"	23,392	6,368	38,788	25,370	20,288	25,615	21,624	17,618	15,157	17,404	1,177	6,019	3,041
% Std Dev	4	8	5	3	2	3	8	10	11	6	8	18	5
UW AR W ME Upgm													
"021209 HIGH 3:1UW AR W ME1"	21,354	5,138	28,715	13,515	10,910	12,185	11,556	9,432	8,103	8,549	1,482	3,983	2,882
"021209 HIGH 3:1UW AR W ME2"	22,787	5,950	28,257	14,683	10,570	12,818	12,813	10,239	8,898	8,711	1,560	4,049	2,220
"021209 HIGH 3:1UW AR W ME3"	21,829	6,086	28,762	13,877	10,746	12,381	12,218	9,283	8,028	8,975	1,488	4,148	2,200
% Std Dev	3	9	1	4	2	3	5	6	6	1	3	2	3
UW NH4F W ME Upgm													
"021209 HIGH 3:1UW NH4F W ME1"	11,983	3,488	16,238	13,505	7,970	14,278	13,788	11,148	9,190	8,619	1,014	2,928	1,556
"021209 HIGH 3:1UW NH4F W ME2"	11,567	3,014	15,863	14,725	8,941	14,843	13,257	10,159	8,944	8,481	1,051	3,085	1,714
"021209 HIGH 3:1UW NH4F W ME3"	11,856	3,315	16,285	13,713	8,101	14,354	13,812	11,101	9,168	9,769	1,043	3,274	1,828
% Std Dev	2	7	1	5	7	2	2	5	5	6	2	6	5
Highly corrected													
UW ME initial Ar. UW Blank													

APPENDIX EXPERIMENT 16A

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn	As	Se	Sr	Zr
UW ME1	879	3,085	1,706	5,586	5,341	5,892	18,588	12,676	1,482	823	1,408	1,302	18,354	10,272
UW ME2	678	9,012	35,978	6,009	6,742	6,952	30,287	21,804	1,973	578	1,304	1,831	18,314	14,945
UW ME3	717	7,778	-18,421	5,257	5,026	5,452	88,082	11,544	1,313	858	2,266	2,619	14,889	9,050
% Std Dev	14	39	459	7	5	11	81	37	22	20	31	35	11	27
UW AR W ME minus UW AR W Blank														
UW AR W ME1	685	3,020	73,756	2,886	1,359	2,791	71,588	13,823	888	1,440	988	673	7,657	7,958
UW AR W ME2	456	4,712	90,810	3,177	1,535	2,762	98,238	13,730	895	1,248	1,012	457	7,938	7,249
UW AR W ME3	517	5,717	70,862	2,722	1,534	2,883	44,350	16,912	475	894	890	847	7,260	6,878
% Std Dev	21	30	14	8	7	2	24	11	31	23	7	39	4	8
UW ME-HF W ME minus UW ME-HF W Blank														
UW ME-HF W ME1	636	9,393	54,897	1,329	1,378	655	38,648	10,897	769	708	477	2,016	8,235	11,002
UW ME-HF W ME2	886	10,428	60,031	1,940	1,280	783	28,628	10,768	647	988	501	1,308	8,043	10,228
UW ME-HF W ME3	854	10,278	41,252	1,498	1,554	603	21,930	14,928	612	840	584	2,208	7,864	10,385
% Std Dev	17	6	19	10	10	18	30	19	17	16	11	33	2	4
Blank Corrected														
Normalised to Average Calcium														
UW ME1-UW BL1	673	3,959	1,695	5,550	5,308	5,853	19,431	12,588	1,472	818	1,458	1,293	18,257	10,205
UW ME2-UW BL2	700	8,281	37,091	6,195	5,827	6,951	31,203	22,479	2,834	588	1,344	1,888	18,881	15,305
UW ME3-UW BL3	701	7,800	-18,978	5,137	5,791	5,327	87,013	11,280	1,283	838	2,214	2,569	14,333	8,843
% Std Dev	13	40	429	9	5	14	78	40	24	18	28	33	14	30
UW AR W ME1-UW AR W BL1	702	3,086	75,631	3,072	1,384	2,882	73,265	14,277	810	1,477	1,013	691	7,852	8,208
UW AR W ME2-UW AR W BL2	441	4,558	87,554	3,074	1,484	2,873	54,404	13,242	837	1,205	878	442	7,877	7,013
UW AR W ME3-UW AR W BL3	522	5,768	71,530	2,746	1,548	2,888	44,750	16,880	478	902	888	855	7,326	6,538
% Std Dev	24	39	11	6	6	4	25	12	31	24	6	31	4	19
UW NHAF W ME1-UW NHAF W BL1	848	9,535	55,697	1,349	1,397	888	30,233	11,062	770	719	484	2,047	8,360	11,888
UW NHAF W ME2-UW NHAF W BL2	865	10,178	58,584	1,801	1,249	774	25,991	10,902	534	883	488	1,082	7,850	9,884
UW NHAF W ME3-UW NHAF W BL3	862	10,375	41,653	1,513	1,589	809	22,143	15,971	618	849	588	2,227	7,940	10,488
% Std Dev	16	4	17	9	11	17	31	20	19	14	11	35	3	6
Percent Standard Deviations														
Matrix Blank														
Av. UW BL %STDEV	4	6	1	4	10	2	8	1	3	2	1	2	10	7
Av. UW AR WASH BL %STDEV	3	6	2	4	7	6	5	5	4	3	2	1	2	7
Av. UW NHAF WASH BL %STDEV	2	3	1	2	4	6	5	3	6	6	8	4	3	7
1 ppm Multi-element Standard														
Av. UW ME %STDEV	1	4	1	8	2	8	13	2	2	2	14	3	9	14
Av. UW AR W ME %STDEV	2	3	1	5	1	1	3	1	2	6	2	1	3	3
Av. UW NHAF W ME %STDEV	2	1	1	6	2	2	5	1	1	3	2	3	2	2
Matrix Blank Corrected														
Av. UW ME-UW BL %STDEV	14	39	458	7	5	11	81	37	22	20	31	35	11	27
Av. UW AR W ME-UW AR W BL %STDEV	21	30	14	8	7	2	24	11	31	23	7	39	4	8
Av. UW NHAF W ME-UW NHAF W BL %STDEV	17	6	19	10	10	18	30	19	17	16	11	33	2	4

APPENDIX EXPERIMENT 16A

Element - Raw Counts	Mo	Cd	Sn	Ba	La	Co	Eu	Dy	Yb	Hf	Hg	Pb	U
UW ME1	20,067	6,266	20,846	14,003	20,887	24,568	24,561	21,118	18,278	18,528	339	2,569	3,043
UW ME2	18,812	5,273	20,234	15,288	20,326	23,873	25,250	21,993	18,141	17,827	249	3,736	3,008
UW ME3	20,482	5,454	23,088	14,380	20,024	24,977	21,688	17,570	15,107	16,687	428	4,131	2,769
% Std Dev	5	9	7	5	2	3	8	10	11	5	26	23	5
UW AR W ME minus UW AR W Blank													
UW AR W ME1	18,343	4,740	15,508	12,295	10,814	12,064	11,508	9,391	8,070	7,970	912	3,413	1,588
UW AR W ME2	19,558	5,581	15,050	13,322	10,474	12,777	12,787	10,188	8,865	8,132	980	3,480	1,727
UW AR W ME3	18,718	4,885	15,544	12,718	10,630	12,250	12,188	9,212	7,989	8,068	888	3,578	1,707
% Std Dev	3	10	2	4	2	3	5	6	5	1	5	2	4
UW NHAF W ME minus UW NHAF W Blank													
UW NHAF W ME1	8,808	3,179	3,200	11,544	7,787	14,116	13,745	11,112	8,153	8,042	313	2,437	1,198
UW NHAF W ME2	8,410	2,705	2,824	12,734	8,757	14,581	13,214	10,122	8,907	8,903	350	2,594	1,304
UW NHAF W ME3	8,738	3,008	3,245	11,722	8,917	14,182	13,789	11,084	9,071	8,191	341	2,783	1,215
% Std Dev	2	8	7	5	7	2	2	5	5	7	5	7	6
Blank Corrected													
Normalised to Average Curium													
UW ME1-UW BL1	19,035	6,225	20,710	13,912	20,760	24,405	24,420	20,980	18,256	18,408	337	2,552	3,024
UW ME2-UW BL2	19,188	5,438	20,880	15,761	20,865	24,405	26,032	21,949	18,702	18,481	256	3,851	3,088
UW ME3-UW BL3	18,894	5,329	22,568	14,051	18,568	24,405	21,073	17,167	14,781	16,285	418	4,036	2,888
% Std Dev	2	9	5	7	4	0	11	12	13	7	24	23	7
UW AR W ME1-UW AR W BL1	18,810	4,880	15,802	12,588	11,089	12,380	11,802	9,630	8,275	8,173	835	3,500	1,839
UW AR W ME2-UW AR W BL2	18,918	5,378	14,569	12,888	10,132	12,350	12,351	9,895	8,576	7,697	838	3,368	1,670
UW AR W ME3-UW AR W BL3	18,887	4,738	15,684	12,831	10,746	12,350	12,279	9,285	8,068	8,168	906	3,610	1,722
% Std Dev	6	7	6	8	5	0	2	3	3	2	3	4	2
UW NHAF W ME1-UW NHAF W BL1	8,839	3,227	3,248	11,688	7,804	14,330	13,953	11,280	8,282	8,183	317	2,474	1,173
UW NHAF W ME2-UW NHAF W BL2	8,208	2,640	2,766	12,428	8,548	14,330	12,898	9,880	9,870	8,880	341	2,532	1,273
UW NHAF W ME3-UW NHAF W BL3	8,834	3,035	3,277	11,238	9,004	14,330	13,802	11,171	9,159	8,281	345	2,819	1,227
% Std Dev	6	10	9	3	7	0	4	7	3	6	4	7	4
Percent Standard Deviations													
Matrix Blank													
Av. UW BL %STDEV	2	34	10	41	24	80	21	24	14	28	3	58	28
Av. UW AR WASH BL %STDEV	6	16	3	23	38	21	38	38	24	4	8	13	12
Av. UW NHAF WASH BL %STDEV	8	8	67	9	17	33	16	53	41	2	18	9	5
Typical Multi-element Standard													
Av. UW ME %STDEV	4	8	5	3	2	3	8	10	11	5	6	15	5
Av. UW AR W ME %STDEV	9	9	1	4	2	3	6	5	6	1	3	2	3
Av. UW NHAF W ME %STDEV	2	7	1	5	7	2	2	5	5	6	2	6	5
Matrix Blank Corrected													
Av. UW ME-UW BL %STDEV	6	9	7	5	2	3	8	10	11	5	28	23	5
Av. UW AR W ME-UW AR W BL %STDEV	3	10	2	4	2	3	5	5	6	1	5	2	4
Av. UW NHAF W ME-UW NHAF W BL %STDEV	2	8	7	5	7	2	2	5	5	7	6	7	6

APPENDIX EXPERIMENT 16A

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn	As	Se	Sr	Zr
Matrix Blank Corrected														
Normalized to Average Carbon														
Ax. UW ME-UW BL %STDEV	13	40	428	9	5	14	79	40	24	18	28	33	14	30
Ax. UW AR W ME-UW AR W BL %STDEV	24	30	11	8	5	4	25	12	31	24	6	31	4	10
Ax. UW NHF ME-UW NHF W BL %STDEV	16	4	17	9	11	17	31	20	19	14	11	35	3	8

APPENDIX EXPERIMENT 16A

Element - Row Counts	Mo	Cd	Sn	Ba	Li	Ce	Eu	Dy	Yb	Hf	Hg	Pb	U
Matrix Blank Corrected													
Normalised to Average Cerium													
Ax. UW ME-UW BL %STDEV	2	8	5	7	4	0	11	12	13	7	24	23	7
Ax. UW AR WME-UW AR WBL %STDEV	0	7	5	1	5	0	2	3	3	2	3	4	2
Ax. UW NHAF ME-UW NHAF WBL %STDEV	5	10	8	3	7	0	4	7	3	6	4	7	4

Experiment 16B.1

"WAS-ED" MATRICES															
AR and NH4F Beaks															
Element - Row Counts															
Glass Standard															
	220,284	194,784	59,436,620	314,958	206,820	401,800	408,763	421,254	231,058	156,843	38,110	28,424	704,094	474,787	629,583
	186,177	172,010	51,381,502	288,945	283,431	344,900	358,392	390,873	197,004	128,515	33,110	28,310	618,338	392,164	513,417
	202,475	179,353	54,340,745	289,957	278,958	370,025	381,598	394,180	210,348	141,065	35,673	27,894	839,008	438,212	614,209
	189,126	174,342	52,800,302	272,884	282,040	350,087	358,380	388,740	188,445	138,148	33,004	27,778	616,183	400,778	529,893
Air Blank															
	5,877	18,471	2,848,357	202	2,883	5,481	48,121	247,886	11,395	1,831	1,884	25,281	613	507	1,492
	5,538	18,628	2,739,808	213	2,883	5,482	48,139	248,151	11,478	1,877	1,887	25,227	649	543	1,364
	5,784	18,270	2,817,840	184	3,170	5,086	48,307	238,878	11,094	1,827	1,888	25,527	643	468	1,417
	5,628	18,208	2,701,878	181	3,380	5,089	48,143	237,853	11,344	1,800	1,892	25,125	672	563	1,361
W Blank															
	8,847	33,351	3,534,742	889	11,022	9,277	372,705	202,255	13,280	7,106	1,477	23,621	2,439	7,098	3,488
	8,498	32,813	3,773,709	793	10,525	9,020	392,898	208,088	13,336	6,717	1,585	23,158	2,398	7,071	3,107
	8,869	34,283	2,870,243	714	9,206	8,719	243,007	213,853	12,650	6,283	1,849	21,771	2,068	6,880	3,432
% Std Dev	4	3	14	9	9	3	24	3	3	6	6	6	9	3	8
W AR W Blank															
	6,247	28,687	3,495,305	635	11,746	8,721	637,849	183,118	11,922	6,802	1,828	22,404	6,302	14,024	2,229
	6,376	28,783	3,297,089	888	11,249	8,574	674,148	201,882	11,980	6,940	1,834	21,178	5,254	15,569	2,083
	8,202	28,738	3,140,376	872	10,750	8,780	682,041	228,748	12,088	6,018	1,817	22,903	5,604	13,370	2,303
% Std Dev	3	2	6	22	4	1	2	9	1	7	6	6	3	8	6
W NH4F W Blank															
	5,772	32,181	2,887,136	715	10,309	7,187	437,387	187,859	11,220	4,750	1,824	22,482	3,654	11,588	2,201
	6,388	31,884	2,365,389	784	9,820	8,783	423,314	212,844	11,330	5,488	1,829	23,828	3,228	9,950	1,835
	5,734	33,033	2,640,376	748	10,777	7,180	441,380	220,031	11,744	4,888	1,702	23,908	3,697	10,522	1,886
% Std Dev	6	2	7	5	5	3	2	5	2	8	6	3	6	12	8
W ME 1 ppm															
	7,407	33,219	3,415,882	618	11,687	10,339	423,785	224,732	14,364	7,091	3,857	23,948	11,845	28,487	15,350
	7,317	35,078	3,394,507	628	11,831	10,100	403,051	233,219	15,283	7,527	3,049	23,688	12,112	24,263	14,588
	7,156	30,751	3,539,886	639	11,309	10,878	413,829	231,651	14,562	6,848	3,273	24,078	11,878	28,687	13,708
% Std Dev	2	7	2	2	2	4	2	2	3	3	12	1	2	9	6
W AR W ME 1 ppm															
	9,887	30,216	3,165,728	923	11,523	9,738	700,848	215,158	13,858	6,753	3,655	22,198	10,713	21,451	17,728
	8,841	28,188	2,871,582	882	11,410	9,803	710,261	223,008	15,406	6,837	3,518	21,882	11,472	20,409	17,363
	8,770	29,831	3,155,888	921	11,885	10,125	707,124	228,478	14,908	7,057	4,058	21,870	12,065	22,401	17,808
% Std Dev	1	4	4	4	3	2	8	3	5	2	8	1	6	5	1
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840	813	11,887	9,924	484,047	218,882	15,183	5,783	2,678	24,134	14,158	21,839	19,617
	8,641	40,440	2,814,554	759	11,381	10,510	472,132	223,394	14,863	5,430	2,889	23,843	14,729	22,148	18,588
	8,882	32,443	2,889,531	833	12,348	10,380	508,540	232,492	16,503	5,861	2,853	24,007	14,971	23,507	20,785
% Std Dev	3	11	2	5	4	3	4	3	6	6	7	1	3	4	6
WV NH4F W ME 1 ppm															
	8,804	37,038	2,717,840												

APPENDIX EXPERIMENT 16B

WASHED MATRICES													
AR and NH4F Ballo													
Element - Raw Counts													
Glass - Standard													
	Cd	Sr	Ba	La	Ce	Eu	Dy	Yb	Hf	Hg	Pb	U	
"02/12/09 HHG GLS STD 1"	170,762	618,441	802,149	470,824	624,435	515,177	362,123	291,677	229,199	344	61,832	56,235	
"02/12/09 HHG GLS STD 2"	132,893	528,258	517,025	389,860	528,447	437,722	285,585	246,074	192,251	368	60,630	60,816	
"02/12/09 HHG GLS STD 3"	158,447	591,530	585,007	439,869	592,401	482,706	334,323	274,448	214,516	348	58,138	54,868	
"02/12/09 HHG GLS STD 4"	138,363	525,190	514,186	388,815	527,378	436,162	289,847	244,089	191,744	328	50,404	48,359	
AR Blank													
"02/12/09 HHG AR BL 1"	272	538	249	60	44	48	63	30	35	282	86	9	
"02/12/09 HHG AR BL 2"	220	542	241	43	36	63	38	26	32	224	81	13	
"02/12/09 HHG AR BL 3"	222	528	178	45	27	39	55	27	28	233	94	8	
"02/12/09 HHG AR BL 4"	244	637	183	58	28	39	55	37	34	303	97	6	
W Blank													
"02/12/09 HHG 3:1W BL 1"	1,951	7,868	3,111	231	283	73	79	58	260	394	1,048	148	
"02/12/09 HHG 3:1W BL 2"	1,243	8,119	3,205	189	282	71	84	54	285	394	1,182	182	
"02/12/09 HHG 3:1W BL 3"	1,117	8,946	3,684	188	228	68	81	80	289	382	1,023	143	
% Std Dev	8	9	12	9	83	4	3	5	7	0	6	14	
W AR W Blank													
"02/12/09 HHG 3:1W AR W BL 1"	2,183	15,584	1,824	74	189	68	80	35	577	448	2,268	604	
"02/12/09 HHG 3:1W AR W BL 2"	1,887	15,058	1,888	82	214	84	83	44	530	469	1,888	617	
"02/12/09 HHG 3:1W AR W BL 3"	1,807	16,187	2,084	88	224	54	80	48	598	408	1,800	705	
% Std Dev	10	4	13	7	8	10	7	15	4	4	12	8	
W NH4F W Blank													
"02/12/09 HHG 3:1W NH4F W BL 1"	726	9,280	2,189	111	183	49	82	42	474	394	1,721	390	
"02/12/09 HHG 3:1W NH4F W BL 2"	858	9,211	2,173	98	174	42	81	41	451	484	1,582	358	
"02/12/09 HHG 3:1W NH4F W BL 3"	885	8,855	1,781	90	175	49	70	42	451	431	1,581	378	
% Std Dev	10	3	11	11	3	8	9	3	5	8	12	4	
W ME 1ppm													
"02/12/09 HHG 3:1W ME 1"	5,072	19,811	12,882	8,709	12,082	10,880	8,693	7,582	7,588	781	3,457	1,887	
"02/12/09 HHG 3:1W ME 2"	8,629	17,774	12,870	10,448	11,751	11,507	9,168	7,980	7,914	701	4,898	1,564	
"02/12/09 HHG 3:1W ME 3"	4,088	17,828	12,706	10,068	11,888	10,868	8,433	7,217	7,802	883	2,809	1,536	
% Std Dev	23	7	1	3	2	3	4	5	3	10	28	6	
W AR W ME 1ppm													
"02/12/09 HHG 3:1W AR W ME 1"	6,058	30,308	10,181	8,888	8,817	7,889	6,286	4,912	5,002	1,408	2,810	1,843	
"02/12/09 HHG 3:1W AR W ME 2"	4,823	27,877	11,256	8,313	8,908	9,110	7,481	5,983	6,500	1,283	2,843	1,851	
"02/12/09 HHG 3:1W AR W ME 3"	5,320	30,648	11,832	8,081	8,912	8,482	7,089	5,747	6,289	1,424	2,872	2,071	
% Std Dev	5	5	8	9	8	8	9	10	14	8	5	7	
W NH4F W ME 1ppm													
"02/12/09 HHG 3:1W NH4F W ME 1"	5,817	37,458	17,885	13,282	18,704	15,016	11,821	9,888	8,840	1,088	3,088	2,888	
"02/12/09 HHG 3:1W NH4F W ME 2"	8,628	40,582	17,848	14,303	17,680	15,656	11,708	9,884	8,645	989	3,725	2,417	
"02/12/09 HHG 3:1W NH4F W ME 3"	5,770	36,851	18,827	13,882	18,884	16,887	12,188	9,808	8,008	1,088	3,305	2,713	
% Std Dev	7	8	3	4	3	3	2	1	4	8	18	6	
Matrix corrected													

APPENDIX EXPERIMENT 16B

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn	As	Se	Sr	Zr	Mo
W ME minus Av. W Blank															
W ME1	641	-200	20,031	-171	1,415	1,325	67,082	16,641	1,276	350	2,287	1,102	9,653	21,537	12,005
W ME2	561	1,657	-11,434	-104	1,380	1,095	68,847	25,127	2,195	628	1,470	916	9,817	17,353	11,043
W ME3	590	-2,988	134,065	-150	1,058	1,865	77,717	23,560	1,464	249	1,703	1,235	9,384	21,747	10,384
% Std Dev	24	-437	161	-4	15	28	13	21	30	45	23	20	2	12	7
W AR W ME minus W AR W Blank															
W AR W ME1	365	817	-125,186	92	231	1,046	27,670	7,288	1,863	468	1,863	357	4,297	6,133	15,509
W AR W ME2	300	-1,232	-515,582	30	129	1,211	37,256	15,680	3,403	650	1,824	121	9,655	6,088	15,142
W AR W ME3	428	532	-354,066	86	719	1,433	34,178	20,829	2,805	771	2,386	109	6,658	7,083	16,388
% Std Dev	16	2,444	-57	49	97	16	16	48	28	34	14	68	11	18	1
W WHAF W ME minus W WHAF W Blank															
W WHAF W ME1	629	4,643	163,537	84	1,595	2,887	48,687	9,748	3,752	774	881	796	10,736	11,449	17,509
W WHAF W ME2	586	7,747	260,280	10	1,069	3,473	37,982	13,550	3,438	411	1,251	305	11,308	11,758	16,580
W WHAF W ME3	807	50	145,227	84	2,048	3,353	74,180	22,248	5,072	933	1,355	689	11,548	13,117	16,788
% Std Dev	26	83	33	73	30	10	34	43	21	38	17	43	4	7	8
Blank Corrected															
Normalised to Average Centum															
W ME1	627	-198	19,610	-167	1,388	1,287	65,258	16,281	1,248	381	2,239	1,079	9,449	21,084	11,753
W ME2	656	1,672	-11,525	-165	1,392	1,104	67,440	23,550	2,214	631	1,489	923	8,804	17,487	11,941
W ME3	395	-2,703	136,793	-152	1,071	1,888	76,717	23,883	1,482	249	1,725	1,251	9,504	22,027	10,487
% Std Dev	23	-637	162	-5	14	29	12	22	31	46	21	20	3	12	6
W AR W ME1	381	900	-158,025	101	255	1,154	30,508	8,035	2,043	514	2,184	372	5,839	6,782	17,689
W AR W ME2	288	-1,177	-303,439	29	123	1,158	35,511	14,889	3,253	621	1,743	116	7,268	4,861	14,474
W AR W ME3	408	508	-129,019	68	881	1,369	32,651	18,708	2,775	738	2,280	104	6,370	6,767	14,702
% Std Dev	18	1,433	-82	63	83	10	8	41	23	88	13	77	5	18	5
W WHAF W ME1	845	4,780	167,875	88	1,835	2,880	50,944	9,894	3,847	704	883	818	11,006	11,738	17,952
W WHAF W ME2	547	7,482	251,689	9	1,053	3,358	38,703	12,717	3,325	597	1,210	285	10,834	11,372	16,035
W WHAF W ME3	685	51	148,586	85	2,065	3,365	74,879	22,457	6,120	941	1,947	675	11,657	13,240	18,988
% Std Dev	26	82	29	74	32	7	36	44	21	48	16	45	4	6	8
Percent Standard Deviations															
Matrix Blank															
Av. W BL %STDEV	4	3	14	9	9	3	24	3	3	6	6	4	9	3	6
Av. W AR W BL %STDEV	3	2	6	22	4	1	2	9	1	7	6	3	3	8	6
Av. W WHAF W BL %STDEV	6	2	7	5	5	3	2	5	2	8	6	3	5	12	8
100% Matrix element Standard															
Av. W ME %STDEV	2	7	2	2	2	4	2	2	3	3	12	1	2	9	6
Av. W AR W ME %STDEV	1	4	4	4	3	2	1	3	5	2	8	1	6	5	1
Av. W WHAF W ME %STDEV	3	11	2	5	4	3	4	3	6	5	7	1	3	4	8
Matrix Blank Corrected															
Av. W ME-W BL %STDEV	24	-537	161	-6	15	28	13	21	30	48	23	20	2	12	7
Av. W AR W ME-W AR W BL %STDEV	18	2,844	-57	48	87	16	15	48	28	34	14	68	11	18	1

Experiment 16B/3

APPENDIX EXPERIMENT 16B

Element - Raw Counts	Cd	Sn	Ba	La	Ce	Eu	Gd	Yb	Hf	Hg	Pb	U
W ME minus Av. W Blank												
W ME1	3,835	12,300	9,292	9,494	11,814	10,819	8,814	7,286	7,253	388	2,389	1,540
W ME2	5,292	10,163	8,471	10,140	11,484	11,437	9,084	7,933	7,639	308	3,811	1,407
W ME3	2,852	10,217	8,366	9,680	11,419	10,828	8,352	7,181	7,327	470	1,721	1,369
% Std Dev	31	11	1	3	2	3	4	6	3	21	41	8
W AR W ME minus W AR W Blank												
W AR W ME1	3,105	14,851	8,319	8,687	9,408	7,638	6,237	4,870	4,447	867	623	1,201
W AR W ME2	2,871	12,251	9,392	8,233	9,607	9,049	7,424	5,941	5,945	794	865	1,209
W AR W ME3	3,388	15,221	9,870	9,703	9,703	8,421	6,350	5,708	5,744	885	885	1,429
% Std Dev	8	11	8	9	8	8	8	10	16	42	18	10
W NHAF W ME minus W NHAF W Blank												
W NHAF W ME1	5,002	28,518	15,854	13,102	16,527	14,868	11,744	9,787	8,388	658	1,281	2,283
W NHAF W ME2	5,713	31,410	15,607	14,203	17,822	15,918	11,630	9,892	8,053	538	1,971	2,042
W NHAF W ME3	4,854	28,708	14,888	13,783	16,787	15,840	12,088	9,881	8,578	638	1,560	2,387
% Std Dev	8	8	3	4	3	3	2	1	4	10	22	7
Blank Corrected												
Normalized to Average Carbon												
W ME1	3,765	12,041	9,093	9,294	11,565	10,581	8,628	7,142	7,130	380	2,320	1,508
W ME2	5,359	10,253	9,656	10,287	11,568	11,538	9,183	8,004	7,707	310	3,845	1,420
W ME3	2,868	10,349	9,428	10,077	11,565	11,088	8,459	7,253	7,421	477	1,743	1,387
% Std Dev	31	8	3	6	0	4	4	6	4	21	41	4
W AR W ME1	3,423	16,186	9,172	7,583	9,270	8,421	6,878	5,989	4,902	1,088	688	1,324
W AR W ME2	2,744	11,711	8,978	7,889	9,270	8,650	7,086	5,678	5,883	758	818	1,156
W AR W ME3	3,218	14,540	9,524	7,843	9,270	8,045	6,039	5,451	5,458	941	845	1,365
% Std Dev	11	16	3	2	0	4	3	3	8	17	11	9
W NHAF W ME1	5,128	28,033	16,060	13,454	16,845	15,348	12,041	10,014	8,800	675	1,313	2,361
W NHAF W ME2	5,525	30,576	15,287	13,758	16,945	15,068	11,248	9,625	8,784	822	1,808	1,975
W NHAF W ME3	6,001	28,880	15,027	13,913	16,945	15,088	12,203	9,854	8,680	642	1,585	2,358
% Std Dev	5	6	3	2	0	3	4	2	1	13	19	10
Percent Standard Deviations												
Blank												
Av. W BL %STDEV	8	6	12	8	13	4	3	5	7	0	8	14
Av. W AR W BL %STDEV	10	4	13	7	8	10	7	15	4	6	12	9
Av. W NHAF W BL %STDEV	10	3	11	11	3	8	9	3	5	8	12	4
1ppm Multi-element Standard												
Av. W ME %STDEV	23	7	1	3	2	3	4	5	3	10	29	6
Av. W AR W ME %STDEV	5	6	8	8	8	8	8	10	14	8	5	7
Av. W NHAF W ME %STDEV	7	6	3	4	3	3	2	1	4	6	10	6
Blank Corrected												
Av. W ME-W BL %STDEV	31	11	1	3	2	3	4	6	3	21	41	8
Av. W AR W ME-W AR W BL %STDEV	8	11	6	9	8	8	9	10	15	12	18	10

Experiment 16B/4

APPENDIX EXPERIMENT 16B

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn	As	Se	Br	Zr	Mo
Av. W NHF MEAN NHF W BL %STDEV	25	63	33	73	30	10	34	43	21	35	17	43	4	7	6
North Blank Corrected															
Normalised to Average Cerium															
Av. W MEAN BL %STDEV	23	637	182	-5	14	29	12	22	31	46	21	20	3	12	6
Av. W AR MEAN AR W BL %STDEV	18	1,432	-52	53	83	10	5	41	23	18	13	77	5	18	9
Av. W NHF MEAN NHF W BL %STDEV	28	82	29	74	92	7	36	44	23	40	15	45	4	8	8

APPENDIX EXPERIMENT 16B

Element - Raw Counts	Cd	Ba	La	Ce	Eu	Dy	Yb	Hf	Ag	Pb	U
Av. W NHAF ME-W NHAF W BL %STDEV	8	8	3	4	3	3	2	1	4	10	22
Matrix Blank Corrected											
Normalized to Average Carbon											
Av. W NHAF BL %STDEV	31	9	3	5	0	4	4	6	4	21	41
Av. W AR W ME-W AR W BL %STDEV	11	18	3	2	0	3	3	3	8	17	11
Av. W NHAF ME-W NHAF W BL %STDEV	5	6	3	2	0	3	4	2	1	13	19
											10

APPENDIX EXPERIMENT 18

Isotope - Raw Counts	Li 7	Mg 24	Ca 44	V 51	Cr 52	Mn 55	Fe 56	Co 59	Ni 60	Cu 65	Zn 66
T021213 HKH GLS STD 1"	49,170	86,700	499,600	142,700	128,200	204,100	268,300	158,700	83,060	44,890	31,300
T021213 HKH AIR BL 1"	6,097	43,050	30,380	200	9,281	4,086	92,610	4,713	57,810	2,143	1,103
T021213 HKH AIR BL 2"	6,266	43,580	29,020	211	10,420	4,539	96,000	4,908	57,100	2,142	1,063
T021213 HKH BLOOD HEAT 1"	6,158	93,280	41,530	419	14,550	11,976	3,454,000	5,171	58,330	4,807	7,888
T021213 HKH BLOOD HEAT 2"	5,708	96,200	42,130	474	17,160	12,250	3,905,000	5,229	58,860	5,313	8,333
T021213 HKH BLOOD HEAT 3"	5,975	94,600	40,930	478	19,270	14,080	3,556,000	5,234	58,080	5,950	8,350
T021213 HKH BLOOD HEAT 4"	5,460	92,710	39,130	490	18,800	11,810	3,926,000	5,336	58,250	5,183	8,481
T021213 HKH BLOOD HEAT 5"	5,611	98,080	41,370	508	17,030	9,439	3,894,000	5,374	58,300	4,306	9,230
T021213 HKH BLOOD AIR 1"	5,142	104,600	43,810	475	19,060	11,200	3,502,000	5,280	58,010	5,641	9,320
T021213 HKH BLOOD AIR 2"	5,101	100,500	38,060	502	14,740	8,533	3,981,000	5,313	59,220	4,264	8,920
T021213 HKH BLOOD AIR 3"	5,364	124,400	40,090	460	16,840	8,338	3,497,000	5,362	59,480	4,139	9,310
T021213 HKH BLOOD AIR 4"	5,342	108,700	38,770	551	18,900	9,867	4,211,000	5,224	58,260	5,377	9,181
T021213 HKH BLOOD AIR 5"	5,469	111,100	38,580	628	18,710	9,405	4,763,000	5,337	59,530	4,642	8,960
T021213 HKH MATRIX BL"	4,989	38,400	31,890	713	13,480	9,868	477,300	4,168	67,860	2,435	1,786
T021213 HKH BLOOD 1" no matrix	5,276	102,900	39,780	245	13,780	5,998	2,779,000	4,441	58,110	5,068	8,127
T021213 HKH BLOOD 2" no matrix	5,511	133,500	52,230	267	14,880	6,401	3,997,000	4,568	59,050	7,003	12,500
T021213 HKH AIR BL 3"	5,574	37,660	23,580	280	12,450	6,069	110,100	4,932	57,120	1,930	1,602
T021213 HKH AIR BL 4"	6,882	38,930	24,410	268	12,770	6,228	111,000	5,120	57,100	1,980	1,653
T021213 HKH GLS STD 2"	42,650	66,880	435,700	122,900	108,000	176,500	235,400	128,300	78,170	37,790	24,760
Air Blank corrected											
T021213 HKH BLOOD HEAT 1"	169	52,280	14,815	179	3,115	6,872	3,350,950	251	1,220	2,746	6,536
T021213 HKH BLOOD HEAT 2"	-282	55,210	15,415	234	5,725	6,946	3,801,950	300	1,750	3,252	6,981
T021213 HKH BLOOD HEAT 3"	-15	53,810	14,215	238	7,835	8,776	3,452,950	314	970	3,889	6,988
T021213 HKH BLOOD HEAT 4"	-530	51,720	11,415	250	7,365	6,508	3,822,950	416	1,140	3,102	7,129
T021213 HKH BLOOD HEAT 5"	-379	57,060	14,655	266	5,886	4,135	3,790,950	454	1,100	2,245	7,878
T021213 HKH BLOOD AIR 1"	-848	63,610	17,055	236	7,625	6,896	3,398,950	360	1,900	3,580	7,868
T021213 HKH BLOOD AIR 2"	-888	58,510	11,335	263	3,305	3,229	3,887,950	393	2,110	2,203	7,588
T021213 HKH BLOOD AIR 3"	-628	83,410	13,375	220	5,405	3,034	3,393,950	442	2,370	2,078	7,958
T021213 HKH BLOOD AIR 4"	-648	67,710	12,055	311	7,465	4,363	4,107,950	304	1,150	3,316	7,828
T021213 HKH BLOOD AIR 5"	-521	70,110	11,885	388	7,275	4,101	4,659,950	417	2,520	2,581	7,598
Normalized to Ba											

Experiment 18/1

APPENDIX EXPERIMENT 18

Isotope - Raw Counts	As 75	Se 78	Mo 98	Cd 114	Sn 120	Sb 121	Ba 138	La 139	Ce 140	Eu 151	Dy 162
¹⁰² 12/13 HKH GLS STD 1"	99,680	11,340	132,309	69,000	214,900	200,200	438,300	500,900	551,600	259,000	68,710
¹⁰² 12/13 HKH AIR BL 1"	4,160	12,590	813	517	750	92	167	88	60	28	14
¹⁰² 12/13 HKH AIR BL 2"	4,254	12,580	868	538	649	81	163	108	85	33	21
¹⁰² 12/13 HKH BLOOD HEAT 1"	15,560	13,380	1,520	844	2,127	321	821	228	210	41	18
¹⁰² 12/13 HKH BLOOD HEAT 2"	16,840	13,780	2,005	631	2,142	407	964	222	144	37	10
¹⁰² 12/13 HKH BLOOD HEAT 3"	41,150	13,920	1,801	571	2,202	261	938	217	259	31	13
¹⁰² 12/13 HKH BLOOD HEAT 4"	22,330	13,860	2,050	561	1,915	217	814	145	109	31	22
¹⁰² 12/13 HKH BLOOD HEAT 5"	20,780	14,360	2,160	684	2,051	341	853	182	128	47	12
¹⁰² 12/13 HKH BLOOD AIR 1"	19,110	13,590	1,824	641	2,201	261	876	176	119	45	16
¹⁰² 12/13 HKH BLOOD AIR 2"	19,860	13,770	1,464	616	2,032	338	808	168	157	34	14
¹⁰² 12/13 HKH BLOOD AIR 3"	29,070	14,830	1,588	614	2,003	448	874	170	173	46	18
¹⁰² 12/13 HKH BLOOD AIR 4"	27,000	14,470	1,585	673	2,381	335	988	242	258	37	17
¹⁰² 12/13 HKH BLOOD AIR 5"	24,150	14,730	1,854	672	2,290	227	939	178	179	39	25
¹⁰² 12/13 HKH MATRIX BL"	30,810	13,080	2,809	640	3,371	251	504	180	133	71	17
¹⁰² 12/13 HKH BLOOD 1" no matrix	12,770	8,787	988	752	974	270	1,672	189	74	32	18
¹⁰² 12/13 HKH BLOOD 2" no matrix	16,230	11,140	1,138	725	1,268	283	2,175	214	82	34	20
¹⁰² 12/13 HKH AIR BL 3"	5,313	12,780	902	540	725	79	191	130	69	30	16
¹⁰² 12/13 HKH AIR BL 4"	5,397	12,100	948	529	684	86	189	143	83	32	23
¹⁰² 12/13 HKH GLS STD 2"	54,920	14,780	111,700	37,550	191,300	169,800	424,000	471,100	519,600	259,000	105,900
Air Blank corrected											
¹⁰² 12/13 HKH BLOOD HEAT 1"	367	786	635	111	1,423	229	643	109	134	10	-2
¹⁰² 12/13 HKH BLOOD HEAT 2"	347	1,205	1,120	98	1,438	315	786	103	68	6	-9
¹⁰² 12/13 HKH BLOOD HEAT 3"	357	1,335	916	38	1,498	170	780	98	183	0	-6
¹⁰² 12/13 HKH BLOOD HEAT 4"	437	1,375	1,165	29	1,211	126	736	26	33	0	2
¹⁰² 12/13 HKH BLOOD HEAT 5"	567	1,775	1,275	151	1,347	250	676	43	53	16	-8
¹⁰² 12/13 HKH BLOOD AIR 1"	417	1,005	738	108	1,497	170	698	57	43	14	-3
¹⁰² 12/13 HKH BLOOD AIR 2"	467	1,185	579	83	1,328	247	630	49	81	3	-5
¹⁰² 12/13 HKH BLOOD AIR 3"	377	2,245	704	81	1,299	356	696	50	97	15	-1
¹⁰² 12/13 HKH BLOOD AIR 4"	407	1,885	610	140	1,677	243	808	123	180	6	-3
¹⁰² 12/13 HKH BLOOD AIR 5"	357	2,145	969	- 139	1,586	136	761	59	103	8	6
Normalized to Ba											

APPENDIX EXPERIMENT 18

Isotope - Raw Counts	Yb 174	Hf 178	Hg 202	Tl 205	Pb 208	Th 232	U 238
T0212/13 HKH GLS STD 1"	100,400	72,550	172	11,530	55,250	84,200	98,250
T0212/13 HKH AIR BL 1"	14	18	108	17	267	10	14
T0212/13 HKH AIR BL 2"	14	8	85	14	153	12	7
T0212/13 HKH BLOOD HEAT 1"	10	31	789	15	1,415	10	203
T0212/13 HKH BLOOD HEAT 2"	18	30	1,028	17	1,200	15	278
T0212/13 HKH BLOOD HEAT 3"	18	32	1,139	23	1,840	26	362
T0212/13 HKH BLOOD HEAT 4"	9	39	551	12	1,389	15	183
T0212/13 HKH BLOOD HEAT 5"	20	53	538	18	1,391	14	219
T0212/13 HKH BLOOD AIR 1"	11	30	864	14	1,397	15	125
T0212/13 HKH BLOOD AIR 2"	14	53	617	15	1,289	12	211
T0212/13 HKH BLOOD AIR 3"	19	50	832	12	1,755	18	134
T0212/13 HKH BLOOD AIR 4"	18	67	485	15	1,785	23	407
T0212/13 HKH BLOOD AIR 5"	22	68	483	15	1,367	18	188
T0212/13 HKH MATRIX BL	14	97	185	18	1,344	19	378
T0212/13 HKH BLOOD 1" no matrix	14	17	1,010	11	1,502	9	8
T0212/13 HKH BLOOD 2" no matrix	15	17	1,178	30	1,316	14	10
T0212/13 HKH AIR BL 3"	14	15	232	13	157	17	5
T0212/13 HKH AIR BL 4"	19	18	269	12	143	11	17
T0212/13 HKH GLS STD 2"	108,300	74,610	281	6,293	47,550	87,250	98,340
Air Blank corrected							
T0212/13 HKH BLOOD HEAT 1"	-3	15	640	2	1,280	-2	192
T0212/13 HKH BLOOD HEAT 2"	4	14	858	4	1,045	4	265
T0212/13 HKH BLOOD HEAT 3"	2	15	981	9	1,885	14	352
T0212/13 HKH BLOOD HEAT 4"	-4	23	402	-2	1,235	3	153
T0212/13 HKH BLOOD HEAT 5"	6	37	380	2	1,238	3	208
T0212/13 HKH BLOOD AIR 1"	-2	14	708	1	1,242	3	114
T0212/13 HKH BLOOD AIR 2"	1	37	459	2	1,114	1	200
T0212/13 HKH BLOOD AIR 3"	6	33	674	-1	1,600	6	123
T0212/13 HKH BLOOD AIR 4"	4	51	328	1	1,630	12	395
T0212/13 HKH BLOOD AIR 5"	8	51	324	2	1,212	7	187
Normalized to Ba							

APPENDIX EXPERIMENT 18

Isotope - Raw Counts	Li 7	Mg 24	Ca 44	V 51	Cr 52	Mn 55	Fe 56	Co 59	Ni 60	Cu 65	Zn 66
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD HEAT 1"	169	52,290	14,815	179	3,115	6,672	3,350,950	251	1,220	2,748	6,536
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD HEAT 2"	-230	45,206	12,622	192	4,598	5,687	3,113,017	246	1,433	2,863	5,716
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD HEAT 3"	-12	45,380	12,033	202	6,632	7,429	2,922,845	266	821	3,282	5,923
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD HEAT 4"	-463	45,213	9,979	218	6,438	5,688	3,341,997	364	997	2,712	6,232
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD HEAT 5"	-361	54,377	13,959	253	5,329	3,639	3,610,816	432	1,133	2,138	7,603
%Stdev	<det limit	9	16	14	27	22	8	27	21	15	11
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD AIR 1"	-761	58,628	15,766	217	7,028	5,434	3,132,643	332	1,751	3,300	7,343
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD AIR 2"	-907	60,737	11,568	268	3,373	3,296	3,968,120	401	2,154	2,248	7,724
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD AIR 3"	-578	77,062	12,357	203	4,894	2,803	3,135,670	408	2,190	1,920	7,352
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD AIR 4"	-516	53,911	9,598	248	5,944	3,474	3,270,755	242	916	2,640	6,233
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD AIR 5"	-440	59,269	10,030	328	6,150	3,467	3,939,361	353	2,130	2,182	6,423
%Stdev	<det limit	14	21	19	25	27	12	19	30	22	9
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD 1" no matrix	5,276	102,900	39,780	245	13,780	5,998	2,779,000	4,441	58,110	5,056	8,127
⁷⁰ Zn/ ⁷¹ Zn HIGH BLOOD 2" no matrix	5,511	133,500	52,230	267	14,880	6,401	3,997,000	4,588	58,050	7,003	12,500
(Median air blank)	5,990	40,990	26,715	240	11,435	5,304	103,050	4,920	57,110	2,051	1,353
Blank corrected	<dl	61,910	13,065	5	2,345	694	2,675,950	<dl	1,000	3,005	6,775
	<dl	92,510	25,515	27	3,445	1,097	3,883,950	<dl	940	4,942	11,148
Normalized to Ba	<det limit	61,910	13,065	5	2,345	694	2,675,950	<det limit	1,000	3,005	6,775
	<det limit	69,211	19,089	20	2,577	821	2,913,224	<det limit	703	3,697	8,340
%Stdev	<det limit	8	26	84	7	12	8	<det limit	25	15	15

APPENDIX EXPERIMENT 18

Isotope - Raw Counts	As 75	Se 76	Mo 98	Cd 114	Sn 120	Sb 121	Ba 138	La 139	Ce 140	Eu 151	Dy 162
¹⁰² 12/13 HGH BLOOD HEAT 1"	367	785	635	111	1,423	228	643	109	134	10	-2
¹⁰² 12/13 HGH BLOOD HEAT 2"	284	987	917	80	1,177	258	643	84	56	5	-7
¹⁰² 12/13 HGH BLOOD HEAT 3"	471	1,130	776	32	1,268	144	643	83	155	0	-5
¹⁰² 12/13 HGH BLOOD HEAT 4"	382	1,202	1,019	25	1,058	110	643	23	29	0	2
¹⁰² 12/13 HGH BLOOD HEAT 5"	540	1,591	1,215	144	1,283	238	643	41	51	15	-7
%Stdev	24	29	24	65	11	33	0	62	66	<det limit	<det limit
¹⁰² 12/13 HGH BLOOD AIR 1"	384	926	681	99	1,379	156	643	53	40	13	-3
¹⁰² 12/13 HGH BLOOD AIR 2"	476	1,209	591	85	1,355	252	643	50	83	3	-5
¹⁰² 12/13 HGH BLOOD AIR 3"	348	2,074	651	75	1,200	328	643	47	89	14	-1
¹⁰² 12/13 HGH BLOOD AIR 4"	324	1,501	645	112	1,335	184	643	88	143	5	-2
¹⁰² 12/13 HGH BLOOD AIR 5"	301	1,813	819	118	1,340	115	643	50	87	7	5
%Stdev	19	30	43	18	6	40	0	38	41	<det limit	<det limit
¹⁰² 12/13 HGH BLOOD 1" no matrix	12,770	8,787	999	752	974	270	1,572	190	74	32	18
¹⁰² 12/13 HGH BLOOD 2" no matrix	16,230	11,140	1,138	725	1,268	283	2,175	214	82	34	20
(Median air blank)	4,784	12,585	885	533	705	91	178	119	78	31	19
Blank corrected	<dl	<dl	115	219	270	178	1,494	71	<dl	<dl	<dl
	<dl	<dl	253	182	564	182	1,997	95	<dl	<dl	<dl
Normalized to Ba	<det limit	<det limit	115	219	270	178	1,494	71	<det limit	<det limit	<det limit
	<det limit	<det limit	189	144	422	144	1,494	71	<det limit	<det limit	<det limit
%Stdev	<det limit	<det limit	35	29	31	15	0	1	<det limit	<det limit	<det limit


APPENDIX EXPERIMENT 18

Isotope - Raw Counts	Yb 174	Hf 178	Hg 202	Tl 205	Pb 208	Th 232	U 238
²⁰² Hf178 HIGH BLOOD HEAT 1"	-3	15	840	2	1,260	-2	192
²⁰² Hf178 HIGH BLOOD HEAT 2"	4	11	710	3	856	3	217
²⁰² Hf178 HIGH BLOOD HEAT 3"	2	13	830	8	1,427	12	298
²⁰² Hf178 HIGH BLOOD HEAT 4"	-4	20	352	-2	1,079	3	133
²⁰² Hf178 HIGH BLOOD HEAT 5"	6	35	352	2	1,178	3	188
%Stddev	<det limit	59	37	<det limit	18	<det limit	29
²⁰² Hf178 HIGH BLOOD AIR 1"	-2	13	650	1	1,145	3	105
²⁰² Hf178 HIGH BLOOD AIR 2"	1	37	468	2	1,137	1	204
²⁰² Hf178 HIGH BLOOD AIR 3"	5	31	622	-1	1,478	6	114
²⁰² Hf178 HIGH BLOOD AIR 4"	3	40	260	1	1,288	10	315
²⁰² Hf178 HIGH BLOOD AIR 5"	7	43	274	2	1,025	6	158
%Stddev	<det limit	37	41	<det limit	14	<det limit	48
²⁰² Hf178 HIGH BLOOD 1" no matrix	14	17	1,010	11	1,602	9	9
²⁰² Hf178 HIGH BLOOD 2" no matrix	15	17	1,178	30	1,318	14	10
(Median air blank)	14	16	158	14	155	11	11
Blank corrected	<dl	<dl	852	<dl	1,447	<dl	<dl
	<dl	<dl	1,020	<dl	1,161	<dl	<dl
Normalized to Ba	<det limit	<det limit	852	<det limit	1,447	<det limit	<det limit
	<det limit	<det limit	783	<det limit	859	<det limit	<det limit
%Stddev	<det limit	<det limit	8	<det limit	35	<det limit	<det limit

APPENDIX EXPERIMENT 13

Isotope - Raw Counts	Mg 24	Ca 44	Cr 52	Mn 55	Fe 56	Cu 65	Zn 66	Sr 88	Mo 98	Sn 128	Ba 138	Pb 207
"02/11/29 HKH GLS STD 1"	94,550	631,500	134,200	203,300	210,500	36,830	21,900	378,700	88,200	145,300	302,700	12,200
"02/11/29 HKH GLS STD 2"	105,400	687,700	151,700	233,900	236,200	43,820	25,290	434,100	113,900	175,000	368,300	16,610
"02/11/29 HKH AIR BL 1"	37,290	48,350	2,361	4,460	38,320	2,936	381	555	276	315	87	23
"02/11/29 HKH AIR BL 2"	34,630	41,380	2,390	4,175	34,240	2,682	347	532	272	254	84	23
"02/11/29 HKH 3:1 UW BL"	62,890	49,770	4,238	5,866	159,200	3,022	5,574	1,775	559	1,589	2,326	2,748
"02/11/29 HKH 3:1 W BL"	54,710	48,510	4,833	5,177	168,200	3,338	6,135	1,899	561	1,749	1,684	2,678
"02/11/29 HKH 3:1 UW OIL"	1,717,000	199,000	23,600	45,040	195,800	49,360	1,253,000	7,618	3,083	22,850	4,233	14,150
"02/11/29 HKH 3:1 W OIL"	1,691,000	198,300	24,160	43,490	194,000	49,220	1,081,000	7,678	3,340	20,840	3,879	13,620
Netto blank corrected												
"02/11/29 HKH 3:1 UW OIL"	1,654,110	149,230	19,364	39,174	38,600	46,328	1,048,428	5,844	2,545	21,261	1,907	11,402
"02/11/29 HKH 3:1 W OIL"	1,636,290	149,790	19,327	38,313	25,800	44,881	1,074,665	5,777	2,779	19,091	2,166	10,942
Element - Raw Counts	Mg	Ca	Cr	Mn	Fe	Cu	Zn	Sr	Mo	Sn	Ba	Pb
"02/11/29 HKH 3:1 UW OIL"	2,093,810	7,006,103	23,107	39,174	38,913	150,416	3,757,799	7,075	10,558	65,218	2,660	54,012
"02/11/29 HKH 3:1 W OIL"	2,071,253	7,032,394	23,063	38,313	28,135	145,718	3,652,583	6,994	11,532	58,561	3,061	51,833
% Std dev.	0.8	0.3	0.1	1.6	24.5	2.2	1.8	0.8	6.2	7.6	9.9	2.9

Experiment 13/1

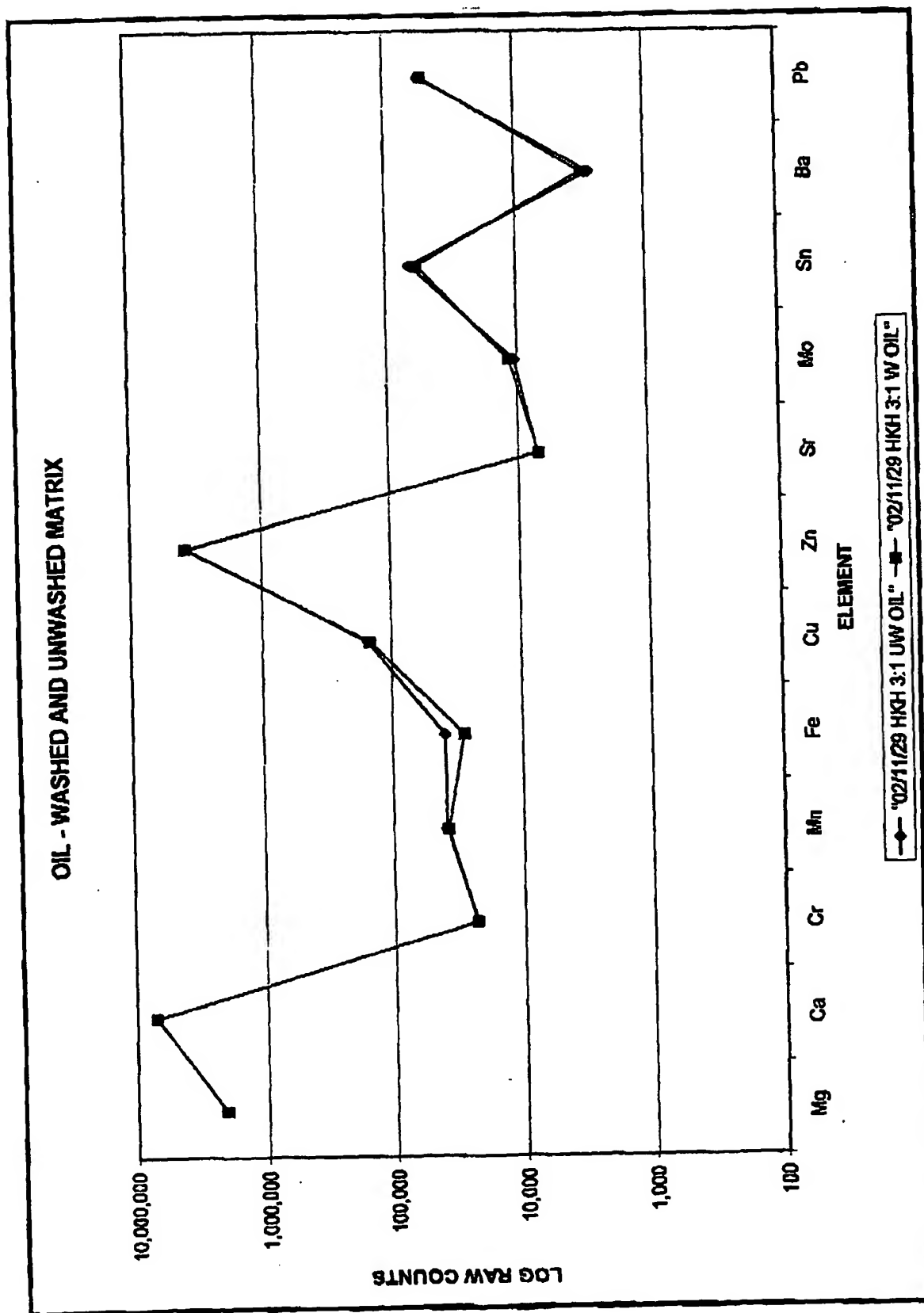

APPENDIX EXPERIMENT 15

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn
T021206 HGH GLS STD 1"	47,490	65,250	314,800	91,720	84,220	129,400	187,500	116,900	27,130	16,370
T021206 HGH GLS STD 2"	41,942	57,354	271,565	78,799	70,067	105,366	164,876	107,511	22,207	11,341
T021206 HGH GLS STD 3"	41,018	65,479	274,201	77,534	74,012	122,292	181,008	115,328	25,437	15,406
T021206 HGH GLS STD 4"	40,624	66,151	266,149	78,201	72,208	116,400	174,192	116,401	23,432	14,478
T021206 HGH GLS STD 5"	38,540	62,445	269,884	75,257	72,523	116,193	178,409	107,941	22,457	14,211
T021206 HGH GLS STD 6"	48,258	68,644	316,450	89,011	86,965	129,212	191,707	118,299	25,952	14,902
T021206 HGH GLS STD 7"	45,680	64,516	298,838	81,820	75,276	117,909	176,308	104,553	21,660	13,946
T021206 HGH GLS STD 8"	47,022	63,180	285,341	78,841	76,177	117,169	175,239	103,415	22,190	13,141
T021206 HGH GLS STD 9"	53,517	66,282	369,379	109,351	100,166	152,187	212,044	115,211	31,787	21,203
T021206 HGH GLS STD 10"	38,574	54,486	230,407	68,320	64,884	100,749	163,475	107,854	21,080	11,485
T021206 HGH GLS STD 11"	47,238	64,809	300,668	91,882	80,741	127,156	189,277	116,902	25,975	17,487
Average Glass Standard	44,827	63,414	290,791	83,704	77,931	121,276	181,276	111,891	24,474	14,906
% Std dev.	10	6	12	13	42	11	7	5	12	18
Cerium Normalized										
T021206 HGH GLS STD 1"	47,490	65,250	314,800	91,720	84,220	129,400	187,500	116,900	27,130	16,370
T021206 HGH GLS STD 2"	51,307	70,161	332,202	96,394	85,713	128,881	201,891	131,517	27,166	13,874
T021206 HGH GLS STD 3"	48,516	77,449	324,325	91,708	87,541	144,846	214,096	138,411	30,087	18,221
T021206 HGH GLS STD 4"	48,406	78,823	317,132	83,181	88,040	138,888	207,559	138,689	27,921	17,251
T021206 HGH GLS STD 5"	47,537	77,072	332,887	92,825	89,453	143,318	220,068	133,139	27,700	17,528
T021206 HGH GLS STD 6"	49,803	70,845	328,598	91,865	89,753	133,366	197,854	122,092	26,692	15,380
T021206 HGH GLS STD 7"	56,074	77,500	350,182	98,287	90,427	141,639	211,791	125,594	26,020	18,793
T021206 HGH GLS STD 8"	56,314	75,642	341,730	94,421	91,231	140,324	208,689	123,852	28,575	15,737
T021206 HGH GLS STD 9"	45,341	55,309	312,952	92,647	84,864	128,939	179,652	97,611	26,931	17,964
T021206 HGH GLS STD 10"	61,511	72,734	307,687	91,235	86,647	134,541	218,306	143,761	28,150	15,338
T021206 HGH GLS STD 11"	46,494	63,787	295,949	90,444	79,469	125,152	186,295	114,784	25,588	17,211
Average Glass Standard	49,890	71,320	324,222	93,157	88,851	135,354	203,152	125,849	27,267	16,512
% Std dev.	7	10	5	2	4	5	6	10	4	8
Drift corrected air blanks										
T021206 HGH AIR BL 1"	3,694	20,190	11,549	152	2,468	3,047	36,855	63,302	808	327
T021206 HGH AIR BL 2"	3,594	20,611	12,257	184	2,720	3,306	40,498	65,600	821	371
T021206 HGH AIR BL 3"	4,690	23,283	12,023	120	3,043	4,094	42,535	69,816	703	406
T021206 HGH AIR BL 4"	4,308	23,124	11,818	144	3,162	4,058	44,044	70,354	725	423
T021206 HGH AIR BL 5"	4,143	25,567	12,948	161	3,528	4,674	48,968	76,409	887	509
T021206 HGH AIR BL 6"	4,059	25,874	13,325	172	3,369	4,495	47,960	76,205	875	454
T021206 HGH AIR BL 7"	4,481	22,498	12,679	172	3,113	4,039	42,523	63,628	752	420
T021206 HGH AIR BL 8"	4,065	21,677	12,662	180	3,087	3,817	42,876	61,863	713	387
T021206 HGH AIR BL 9"	3,888	21,353	11,540	145	2,790	3,535	38,598	68,968	814	395
T021206 HGH AIR BL 10"	3,871	21,358	12,933	192	2,837	3,477	42,447	66,395	833	389
Average	4,083	22,651	12,372	162	3,010	3,854	42,790	68,034	799	406
Element - Raw Counts										

Experiment 15/1

APPENDIX EXPERIMENT 15

Element - Raw Counts	Ga	As	Se	Sr	Zr	Mo	Cd	Sn	Ba	La
"02/12/06 HGH GLS STD 1"	97,640	17,950	5,077	233,800	108,100	64,430	10,920	108,900	235,800	263,700
"02/12/06 HGH GLS STD 2"	80,014	14,411	4,775	203,320	92,902	51,228	8,370	83,587	198,445	226,040
"02/12/06 HGH GLS STD 3"	85,443	14,494	5,015	211,943	98,237	58,704	8,090	94,701	217,337	223,229
"02/12/06 HGH GLS STD 4"	81,681	15,295	5,583	207,032	95,489	54,328	7,412	91,701	195,728	214,195
"02/12/06 HGH GLS STD 5"	84,958	14,524	4,985	200,520	94,666	53,912	7,045	88,138	194,117	211,640
"02/12/06 HGH GLS STD 6"	98,635	16,941	5,688	220,573	105,481	64,205	8,313	100,315	229,030	249,680
"02/12/06 HGH GLS STD 7"	82,557	14,885	5,368	201,842	91,566	53,590	7,146	85,117	199,608	224,982
"02/12/06 HGH GLS STD 8"	83,899	15,447	5,247	193,725	87,890	53,525	7,710	90,454	199,979	212,689
"02/12/06 HGH GLS STD 9"	120,855	20,448	5,202	281,520	131,249	79,064	12,516	131,784	273,378	313,117
"02/12/06 HGH GLS STD 10"	70,750	13,024	4,770	187,271	72,202	48,183	6,450	78,170	170,847	180,676
"02/12/06 HGH GLS STD 11"	97,820	18,164	4,905	228,149	103,497	68,426	11,540	112,414	245,398	268,020
Average Glass Standard	88,479	15,882	5,201	213,609	98,121	59,062	8,682	98,753	214,333	235,271
% Std dev.	14	13	8	13	14	16	22	18	13	15
Cerium Normalized										
"02/12/06 HGH GLS STD 1"	97,640	17,950	5,077	233,800	108,100	64,430	10,920	108,900	235,800	263,700
"02/12/06 HGH GLS STD 2"	97,880	17,828	5,841	248,719	113,646	62,667	10,289	102,252	240,309	276,512
"02/12/06 HGH GLS STD 3"	101,082	17,144	6,641	250,885	116,184	69,435	9,569	112,012	267,065	284,034
"02/12/06 HGH GLS STD 4"	97,339	18,224	6,652	246,691	113,781	64,734	8,832	109,266	233,221	255,228
"02/12/06 HGH GLS STD 5"	104,791	17,915	6,149	247,330	116,765	66,497	8,690	108,714	239,433	281,046
"02/12/06 HGH GLS STD 6"	101,797	17,484	5,870	227,645	108,842	68,263	8,580	103,531	236,373	257,685
"02/12/06 HGH GLS STD 7"	99,171	17,881	6,448	242,484	109,994	64,376	8,566	102,247	239,781	270,273
"02/12/06 HGH GLS STD 8"	100,479	18,500	6,284	232,008	105,342	64,102	9,234	108,328	239,498	254,721
"02/12/06 HGH GLS STD 9"	102,402	17,323	4,407	238,515	111,199	66,986	10,804	111,652	231,816	265,285
"02/12/06 HGH GLS STD 10"	94,400	17,382	6,370	223,375	98,419	64,317	8,613	105,724	228,150	241,278
"02/12/06 HGH GLS STD 11"	96,278	17,877	4,828	224,564	101,868	67,348	11,457	110,643	241,531	283,797
Average Glass Standard	99,393	17,766	6,870	237,789	106,105	65,650	9,676	107,388	238,434	261,232
% Std dev.	3	2	12	4	5	3	11	3	3	3
Drift corrected air blanks										
"02/12/06 HGH AIR BL 1"	280	832	3,019	266	108	284	19	165	122	32
"02/12/06 HGH AIR BL 2"	345	971	3,304	275	128	326	28	182	152	44
"02/12/06 HGH AIR BL 3"	306	908	3,129	320	97	362	20	206	147	38
"02/12/06 HGH AIR BL 4"	315	928	3,241	293	103	353	19	186	153	38
"02/12/06 HGH AIR BL 5"	385	1,091	3,859	314	134	382	25	231	158	46
"02/12/06 HGH AIR BL 6"	388	1,057	4,001	309	122	380	23	223	170	41
"02/12/06 HGH AIR BL 7"	368	939	3,288	286	128	354	28	184	149	40
"02/12/06 HGH AIR BL 8"	368	947	3,228	277	132	350	22	193	156	41
"02/12/06 HGH AIR BL 9"	307	918	2,937	285	113	330	23	189	136	39
"02/12/06 HGH AIR BL 10"	359	954	3,432	278	133	333	23	182	141	41
Average	342	939	3,345	281	120	347	23	196	148	40
Element - Raw Counts										

APPENDIX EXPERIMENT 15

Element - Raw Counts	Ce	Eu	Dy	Yb	Hf	Hg	Pb	U
"02/12/06 HKH GLS STD 1"	305,900	145,300	57,670	61,330	42,160	367	36,940	54,670
"02/12/06 HKH GLS STD 2"	250,064	127,020	51,079	52,610	36,902	412	27,794	43,100
"02/12/06 HKH GLS STD 3"	258,624	121,397	47,081	47,634	32,567	525	25,563	43,145
"02/12/06 HKH GLS STD 4"	256,723	114,252	45,268	46,559	31,276	483	25,892	43,881
"02/12/06 HKH GLS STD 5"	248,005	111,211	45,510	45,148	30,859	416	22,469	38,761
"02/12/06 HKH GLS STD 6"	256,387	135,559	53,917	56,454	38,642	426	30,187	54,131
"02/12/06 HKH GLS STD 7"	254,661	121,501	47,787	51,349	35,756	251	26,924	42,254
"02/12/06 HKH GLS STD 8"	255,423	116,918	45,224	47,694	33,289	269	27,444	45,918
"02/12/06 HKH GLS STD 9"	361,055	165,458	65,438	69,503	47,354	338	34,320	54,089
"02/12/06 HKH GLS STD 10"	229,068	101,413	38,879	40,738	27,482	325	21,044	41,430
"02/12/06 HKH GLS STD 11"	310,788	147,527	56,844	61,549	42,538	421	32,514	60,233
Average Glass Standard	275,155	127,860	50,418	52,833	36,266	387	28,281	47,419
% Std dev.	13	14	14	16	16	20	16	14
Carlum Normalized								
"02/12/06 HKH GLS STD 1"	305,900	145,300	57,670	61,330	42,160	367	36,940	54,670
"02/12/06 HKH GLS STD 2"	305,900	155,382	62,485	64,502	45,142	504	34,001	52,724
"02/12/06 HKH GLS STD 3"	305,900	143,588	55,687	56,341	38,520	621	30,236	51,031
"02/12/06 HKH GLS STD 4"	305,900	136,137	53,939	55,478	37,268	576	30,852	52,287
"02/12/06 HKH GLS STD 5"	305,900	137,172	56,134	55,688	39,186	613	27,715	47,810
"02/12/06 HKH GLS STD 6"	305,900	139,905	55,846	58,264	39,881	439	31,155	55,866
"02/12/06 HKH GLS STD 7"	305,900	145,953	57,405	61,684	42,952	302	32,342	50,758
"02/12/06 HKH GLS STD 8"	305,900	140,023	54,161	57,119	39,668	347	32,868	54,992
"02/12/06 HKH GLS STD 9"	305,900	140,182	56,440	59,225	40,120	287	29,077	45,826
"02/12/06 HKH GLS STD 10"	305,900	135,428	52,053	54,401	36,699	434	28,103	55,325
"02/12/06 HKH GLS STD 11"	305,900	145,202	55,751	60,578	41,868	415	32,002	59,284
Average Glass Standard	305,900	143,207	56,034	58,610	40,242	437	31,390	52,779
% Std dev.	6	4	5	5	8	24	8	7
Drift corrected air blanks								
"02/12/06 HKH AIR BL 1"	11	21	6	9	9	262	65	8
"02/12/06 HKH AIR BL 2"	18	23	12	11	10	302	72	8
"02/12/06 HKH AIR BL 3"	14	23	8	10	9	319	74	10
"02/12/06 HKH AIR BL 4"	13	23	8	9	7	317	63	7
"02/12/06 HKH AIR BL 5"	22	29	12	12	7	453	69	11
"02/12/06 HKH AIR BL 6"	14	22	11	10	10	432	63	4
"02/12/06 HKH AIR BL 7"	11	20	8	9	9	228	62	8
"02/12/06 HKH AIR BL 8"	15	19	8	6	11	223	61	8
"02/12/06 HKH AIR BL 9"	16	25	10	11	9	312	74	10
"02/12/06 HKH AIR BL 10"	14	21	8	11	11	287	69	7
Average	15	23	9	10	9	317	67	8
Element - Raw Counts								

APPENDIX EXPERIMENT 15

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn
"02/12/06 HGH SVEN OIL BL 2"	3,821	235,018	41,480	687	10,950	5,483	150,553	73,186	1,189	167,148
"02/12/06 HGH SVEN OIL BL 3"	3,888	201,744	39,846	683	8,118	5,682	157,177	73,459	2,225	143,782
"02/12/06 HGH SVEN OIL WED 1"	3,742	190,075	33,354	584	8,487	9,138	161,819	71,368	4,519	137,949
"02/12/06 HGH SVEN OIL WED 2"	4,128	195,788	34,980	711	10,163	6,988	268,814	74,881	3,343	143,612
"02/12/06 HGH SVEN OIL THUR 1"	4,719	276,925	62,367	745	11,530	11,882	563,657	81,668	7,485	182,213
"02/12/06 HGH SVEN OIL THUR 2"	4,824	258,792	45,329	1,031	13,952	10,300	534,454	81,280	10,451	176,523
"02/12/06 HGH SVEN OIL FRI 1"	4,810	288,334	68,590	2,446	18,629	19,376	528,987	77,330	18,538	221,004
"02/12/06 HGH SVEN OIL FRI 2"	5,029	238,601	45,334	1,105	13,938	10,325	508,588	83,148	16,947	189,439
"02/12/06 HGH JOHN OIL WED 1"	5,385	580,487	55,967	346	13,776	19,858	234,195	82,858	20,828	304,144
"02/12/06 HGH JOHN OIL WED 2"	5,147	604,576	60,978	417	18,936	22,912	306,614	88,485	20,456	314,960
"02/12/06 HGH JOHN OIL THUR 1"	4,518	408,802	44,199	448	13,941	16,549	270,544	83,824	13,895	212,212
"02/12/06 HGH JOHN OIL THUR 2"	4,282	418,970	45,512	425	14,472	18,970	213,334	83,907	14,674	218,577
"02/12/06 HGH JOHN OIL FRI 1"	4,222	467,852	49,288	415	18,658	18,436	214,237	85,038	15,914	242,640
"02/12/06 HGH JOHN OIL FRI 2"	4,394	455,915	49,409	461	17,280	19,570	285,871	84,323	15,748	266,535
"02/12/06 HGH RYAN OIL WED 1"	5,532	409,850	50,572	619	23,880	10,525	470,647	82,108	5,760	358,710
"02/12/06 HGH RYAN OIL WED 2"	5,315	269,141	37,981	908	17,157	11,958	554,841	87,060	5,272	286,034
"02/12/06 HGH RYAN OIL THUR 1"	5,135	585,490	64,218	607	27,065	15,071	555,053	85,204	8,876	493,518
"02/12/06 HGH RYAN OIL THUR 2"	5,015	413,186	48,900	672	17,325	9,512	387,147	84,519	5,325	391,813
"02/12/06 HGH RYAN OIL FRI 1"	4,885	618,761	67,912	560	24,139	10,701	424,569	85,514	8,871	680,379
"02/12/06 HGH RYAN OIL FRI 2"	5,093	601,154	95,593	588	27,817	11,352	475,080	86,087	7,080	673,978
"02/12/06 HGH DAVE OIL WED 1"	6,284	54,719	49,159	583	14,019	18,012	485,381	82,729	4,151	188,777
"02/12/06 HGH DAVE OIL WED 2"	5,625	53,475	49,934	548	11,967	10,856	418,508	81,447	3,872	168,231
"02/12/06 HGH DAVE OIL THUR 1"	5,731	68,496	61,802	815	12,045	11,243	339,597	83,326	4,070	235,505
"02/12/06 HGH DAVE OIL THUR 2"	5,619	55,528	61,737	608	12,589	9,874	268,282	84,838	4,189	195,804
"02/12/06 HGH DAVE OIL FRI 1"	5,678	97,435	172,212	508	21,019	13,060	357,339	85,822	6,315	200,078
"02/12/06 HGH DAVE OIL FRI 2"	5,618	91,816	162,196	421	19,831	10,788	198,769	85,450	4,682	176,146
"02/12/06 HGH SCOTT OIL WED 1"	7,178	359,173	78,240	921	27,762	88,903	11,838,207	119,587	8,650	1,591,134
"02/12/06 HGH SCOTT OIL WED 2"	6,901	218,524	52,418	820	17,884	52,411	10,702,080	104,254	5,678	1,243,243
"02/12/06 HGH SCOTT OIL THUR 1"	6,355	197,533	50,333	900	18,788	72,574	9,738,842	98,617	6,188	843,194
"02/12/06 HGH SCOTT OIL THUR 2"	6,488	241,758	64,444	1,435	23,479	98,567	13,984,018	111,528	8,980	1,683,237
"02/12/06 HGH SCOTT OIL FRI 1"	6,356	168,149	48,849	1,069	18,013	86,219	8,987,866	101,870	5,466	1,090,938
"02/12/06 HGH SCOTT OIL FRI 2"	6,385	220,839	59,311	1,015	22,930	75,366	10,140,408	109,390	7,714	1,704,582
Average Air Blank Connected										
Sven Reference Oil										
"02/12/06 HGH SVEN OIL BL 2"	-281	212,487	28,117	525	7,980	1,828	107,823	5,161	398	166,742
"02/12/06 HGH SVEN OIL BL 3"	-195	179,194	27,474	521	5,108	1,828	114,448	5,425	1,433	143,376
Sven Engine Oil										
"02/12/06 HGH SVEN OIL WED 1"	-341	167,524	20,581	432	5,458	5,284	318,890	3,334	3,728	137,443

APPENDIX EXPERIMENT 15

Element - Raw Counts	Ga	As	Se	Sr	Zr	Mo	Cd	Sn	Ba	La
"02/12/06 HIGH SVEN OIL BL 2"	24,526	1,977	4,304	1,917	9,882	887	127	919	1,035	82
"02/12/06 HIGH SVEN OIL BL 3"	29,525	2,031	4,800	1,882	12,522	676	63	1,128	738	93
"02/12/06 HIGH SVEN OIL WED 1"	25,965	1,828	3,601	2,130	4,661	820	56	3,242	3,040	1,147
"02/12/06 HIGH SVEN OIL WED 2"	30,868	2,157	3,907	1,831	5,203	795	66	2,729	3,389	1,414
"02/12/06 HIGH SVEN OIL THUR 1"	32,120	2,818	5,043	4,285	9,881	1,349	98	1,527	5,424	1,164
"02/12/06 HIGH SVEN OIL THUR 2"	36,567	2,537	4,582	4,238	8,228	1,274	153	3,330	5,778	1,583
"02/12/06 HIGH SVEN OIL FRI 1"	37,388	2,604	4,194	7,829	10,880	1,986	84	9,445	4,276	1,476
"02/12/06 HIGH SVEN OIL FRI 2"	40,695	2,638	4,626	3,411	18,410	1,195	204	3,850	4,497	855
"02/12/06 HIGH JOHN OIL WED 1"	9,870	1,773	3,967	2,911	4,180	2,368	65	11,459	1,955	148
"02/12/06 HIGH JOHN OIL WED 2"	12,719	1,920	4,405	3,386	5,247	2,998	84	11,801	2,433	210
"02/12/06 HIGH JOHN OIL THUR 1"	20,970	1,731	3,924	2,411	8,571	1,631	60	8,203	1,795	269
"02/12/06 HIGH JOHN OIL THUR 2"	19,686	1,807	3,771	2,600	6,313	1,807	38	11,414	1,800	430
"02/12/06 HIGH JOHN OIL FRI 1"	19,641	1,859	4,148	2,595	4,743	1,683	49	7,343	1,379	85
"02/12/06 HIGH JOHN OIL FRI 2"	18,636	2,021	4,138	2,730	3,164	2,004	85	8,186	1,591	538
"02/12/06 HIGH RYAN OIL WED 1"	34,832	1,845	4,188	5,115	1,478	1,855	425	2,205	14,046	186
"02/12/06 HIGH RYAN OIL WED 2"	40,453	1,825	4,163	4,187	2,088	1,879	87	2,916	11,678	408
"02/12/06 HIGH RYAN OIL THUR 1"	30,594	2,186	5,082	4,212	1,571	1,458	135	3,713	9,008	326
"02/12/06 HIGH RYAN OIL THUR 2"	38,800	2,043	4,710	3,317	2,045	1,613	158	4,642	3,163	227
"02/12/06 HIGH RYAN OIL FRI 1"	26,133	2,506	4,666	5,494	826	2,030	191	2,728	9,848	205
"02/12/06 HIGH RYAN OIL FRI 2"	19,987	2,357	4,752	7,552	1,184	2,647	143	2,640	93,280	211
"02/12/06 HIGH DAVE OIL WED 1"	39,625	1,871	3,984	2,142	4,657	1,311	86	3,028	2,242	226
"02/12/06 HIGH DAVE OIL WED 2"	38,853	1,877	3,815	2,218	4,073	972	58	3,465	2,100	236
"02/12/06 HIGH DAVE OIL THUR 1"	64,661	2,107	4,433	3,038	5,477	2,575	138	2,625	2,087	193
"02/12/06 HIGH DAVE OIL THUR 2"	43,001	2,254	4,543	2,688	4,590	1,174	76	1,854	651	101
"02/12/06 HIGH DAVE OIL FRI 1"	32,320	2,838	4,719	5,484	3,744	1,266	156	1,603	1,583	108
"02/12/06 HIGH DAVE OIL FRI 2"	32,793	2,955	4,653	5,137	3,748	1,220	155	1,857	1,610	110
"02/12/06 HIGH SCOTT OIL WED 1"	31,712	2,523	4,503	4,233	8,295	3,284	147	4,314	12,088	116
"02/12/06 HIGH SCOTT OIL WED 2"	48,230	2,385	4,437	2,724	9,820	5,003	86	4,241	10,009	124
"02/12/06 HIGH SCOTT OIL THUR 1"	48,711	2,680	4,320	2,559	8,751	2,085	88	4,173	11,778	365
"02/12/06 HIGH SCOTT OIL THUR 2"	48,853	2,900	4,379	3,483	8,709	4,374	233	6,877	18,437	222
"02/12/06 HIGH SCOTT OIL FRI 1"	55,886	3,031	4,616	2,686	11,979	2,139	217	4,371	11,678	280
"02/12/06 HIGH SCOTT OIL FRI 2"	44,363	3,122	4,509	3,446	10,427	2,297	158	4,528	14,550	869
Average Air Blank Corrected										
Sven Reference Oil										
"02/12/06 HIGH SVEN OIL BL 2"	24,184	1,019	959	1,626	9,742	550	105	723	887	42
"02/12/06 HIGH SVEN OIL BL 3"	23,183	1,072	1,255	1,370	12,402	330	40	930	889	53
Sven Engine Oil										
"02/12/06 HIGH SVEN OIL WED 1"	25,623	970	258	1,838	4,542	473	34	3,046	2,881	1,107

APPENDIX EXPERIMENT 15

Element - Raw Counts	Ce	Eu	Dy	Yb	Hf	Hg	Pb	U
"02/12/06 HKH SVEN OIL BL 2"	120	32	13	19	103	604	440	82
"02/12/06 HKH SVEN OIL BL 3"	66	27	15	18	33	606	438	102
"02/12/06 HKH SVEN OIL WED 1"	314	28	15	14	97	502	41988	155
"02/12/06 HKH SVEN OIL WED 2"	108	28	22	19	94	498	43,195	113
"02/12/06 HKH SVEN OIL THUR 1"	197	32	22	14	107	749	66,643	19
"02/12/06 HKH SVEN OIL THUR 2"	673	42	24	22	234	685	65,056	136
"02/12/06 HKH SVEN OIL FRI 1"	528	44	23	28	109	488	77,539	171
"02/12/06 HKH SVEN OIL FRI 2"	845	46	27	25	181	508	59,094	165
"02/12/06 HKH JOHN OIL WED 1"	81	28	10	15	32	876	21,561	53
"02/12/06 HKH JOHN OIL WED 2"	191	30	13	16	74	833	21,248	68
"02/12/06 HKH JOHN OIL THUR 1"	96	24	11	17	110	687	11,754	86
"02/12/06 HKH JOHN OIL THUR 2"	139	26	16	19	122	889	13,188	80
"02/12/06 HKH JOHN OIL FRI 1"	72	25	12	14	24	736	12,871	70
"02/12/06 HKH JOHN OIL FRI 2"	112	23	12	10	107	801	15,171	60
"02/12/06 HKH RYAN OIL WED 1"	300	28	12	18	44	730	13,378	156
"02/12/06 HKH RYAN OIL WED 2"	770	31	19	21	60	778	10,142	190
"02/12/06 HKH RYAN OIL THUR 1"	248	29	18	15	148	1,023	15,181	118
"02/12/06 HKH RYAN OIL THUR 2"	502	40	14	23	58	1,018	10,079	155
"02/12/06 HKH RYAN OIL FRI 1"	396	35	18	42	28	721	9,711	115
"02/12/06 HKH RYAN OIL FRI 2"	233	26	18	21	34	742	11,567	142
"02/12/06 HKH DAVE OIL WED 1"	195	27	15	17	83	450	34,766	160
"02/12/06 HKH DAVE OIL WED 2"	128	25	13	28	82	460	41,522	145
"02/12/06 HKH DAVE OIL THUR 1"	574	25	14	27	78	568	37,894	273
"02/12/06 HKH DAVE OIL THUR 2"	96	78	14	19	33	586	35,358	144
"02/12/06 HKH DAVE OIL FRI 1"	85	27	17	22	17	487	40,138	102
"02/12/06 HKH DAVE OIL FRI 2"	59	27	16	21	18	465	43,944	107
"02/12/06 HKH SCOTT OIL WED 1"	281	29	19	28	181	630	7,987	164
"02/12/06 HKH SCOTT OIL WED 2"	130	29	17	18	44	625	8,630	164
"02/12/06 HKH SCOTT OIL THUR 1"	108	28	16	18	107	608	6,244	198
"02/12/06 HKH SCOTT OIL THUR 2"	96	37	26	22	64	744	7,980	173
"02/12/06 HKH SCOTT OIL FRI 1"	108	35	18	24	114	508	5,951	185
"02/12/06 HKH SCOTT OIL FRI 2"	152	33	18	19	114	639	6,900	151
Average Air Blank Corrected								
Sven Reference Oil								
"02/12/06 HKH SVEN OIL BL 2"	105	9	4	9	94	287	372	74
"02/12/06 HKH SVEN OIL BL 3"	50	5	6	6	24	289	371	94
Sven Engine Oil								
"02/12/06 HKH SVEN OIL WED 1"	300	4	6	4	88	186	41,920	147

APPENDIX EXPERIMENT 15

Element - Raw Counts	Li	Mg	Ca	V	Cr	Mn	Fe	Ni	Cu	Zn
"02/12/06 HGH SVEN OIL WED 2"	45	174,218	22,567	549	7,154	3,114	224,084	6,847	2,560	143,208
"02/12/06 HGH SVEN OIL THUR 1"	636	254,375	49,996	583	8,541	8,008	515,927	13,032	6,682	181,807
"02/12/06 HGH SVEN OIL THUR 2"	741	217,242	33,156	869	10,943	8,446	491,725	13,028	9,658	176,117
"02/12/06 HGH SVEN OIL FRI 1"	727	283,784	58,218	2,294	13,820	15,522	487,257	9,296	17,745	220,598
"02/12/06 HGH SVEN OIL FRI 2"	948	216,051	32,982	943	10,926	6,871	463,856	15,114	16,154	198,033
John Engine Oil										
"02/12/06 HGH JOHN OIL WED 1"	1,302	557,937	43,586	184	10,766	16,102	181,465	14,824	20,035	303,738
"02/12/06 HGH JOHN OIL WED 2"	1,065	581,825	48,603	256	13,826	18,058	253,884	18,451	19,663	314,554
"02/12/06 HGH JOHN OIL THUR 1"	435	387,252	31,826	286	10,931	12,695	227,815	15,780	13,102	211,806
"02/12/06 HGH JOHN OIL THUR 2"	189	396,420	33,139	263	11,462	13,116	170,605	15,873	13,881	218,171
"02/12/06 HGH JOHN OIL FRI 1"	139	446,371	36,915	283	15,648	14,581	171,508	18,004	15,121	242,234
"02/12/06 HGH JOHN OIL FRI 2"	311	433,384	37,037	298	14,281	16,715	243,141	16,288	14,955	265,129
Ryan Engine Oil										
"02/12/06 HGH RYAN OIL WED 1"	1,449	387,300	38,200	457	20,670	6,871	427,918	14,074	4,987	359,304
"02/12/06 HGH RYAN OIL WED 2"	1,232	246,591	25,609	743	14,147	8,105	512,112	19,026	4,479	285,828
"02/12/06 HGH RYAN OIL THUR 1"	1,052	562,939	51,846	445	24,055	11,217	522,323	17,170	8,083	493,112
"02/12/06 HGH RYAN OIL THUR 2"	932	399,515	38,528	510	14,315	5,658	344,417	16,485	4,532	391,207
"02/12/06 HGH RYAN OIL FRI 1"	903	597,211	56,539	397	21,129	6,947	381,839	17,480	8,078	658,972
"02/12/06 HGH RYAN OIL FRI 2"	980	578,604	83,221	425	24,607	7,498	432,381	18,053	8,287	673,571
Dave Engine Oil										
"02/12/06 HGH DAVE OIL WED 1"	2,211	32,168	38,786	420	11,009	14,158	442,652	14,694	3,358	166,371
"02/12/06 HGH DAVE OIL WED 2"	1,542	30,924	37,582	385	8,958	7,101	378,178	13,413	3,079	167,825
"02/12/06 HGH DAVE OIL THUR 1"	1,648	45,948	49,530	652	9,035	7,389	286,868	15,291	3,277	235,088
"02/12/06 HGH DAVE OIL THUR 2"	1,536	32,977	49,365	444	9,580	5,820	223,553	16,804	3,398	195,388
"02/12/06 HGH DAVE OIL FRI 1"	1,595	74,885	158,840	345	18,009	9,205	314,809	17,887	5,522	188,672
"02/12/06 HGH DAVE OIL FRI 2"	1,535	68,365	149,823	259	16,622	6,904	158,040	17,416	3,900	176,740
Scott Engine Oil										
"02/12/06 HGH SCOTT OIL WED 1"	3,056	338,623	65,868	759	24,753	95,049	11,786,478	51,553	8,858	1,590,728
"02/12/06 HGH SCOTT OIL WED 2"	2,818	195,974	40,044	658	14,854	48,557	10,659,351	38,270	4,885	1,242,837
"02/12/06 HGH SCOTT OIL THUR 1"	2,272	174,883	37,961	738	15,778	68,720	9,694,113	31,583	5,385	942,788
"02/12/06 HGH SCOTT OIL THUR 2"	2,405	218,208	52,072	1,333	20,469	92,713	13,941,288	43,484	8,187	1,682,831
"02/12/06 HGH SCOTT OIL FRI 1"	2,273	145,599	36,477	897	15,004	62,365	8,945,136	33,836	4,693	1,090,532
"02/12/06 HGH SCOTT OIL FRI 2"	2,303	186,288	46,838	853	18,921	71,511	10,097,676	41,358	6,921	1,704,156
Average Engine Oil - John										
Average Engine Oil - John	575	467,018	38,519	258	12,836	15,211	211,403	18,538	16,126	258,272
Average Engine Oil - Scott										
Average Engine Oil - Scott	2,528	211,446	48,660	873	18,463	73,152	10,855,674	39,674	6,490	1,375,645

Experiment 15/7

APPENDIX EXPERIMENT 15

Element - Raw Counts	Ga	As	Se	Sr	Zr	Mo	Cd	Sn	Ba	La
"02/12/06 HIGH SVEN OIL WED 2"	30,524	1,198	582	1,340	5,083	448	44	2,533	3,251	1,374
"02/12/06 HIGH SVEN OIL THUR 1"	31,778	1,859	1,698	3,994	9,761	1,002	76	1,331	5,275	1,124
"02/12/06 HIGH SVEN OIL THUR 2"	36,225	1,579	1,237	3,846	9,109	927	130	3,134	5,629	1,543
"02/12/06 HIGH SVEN OIL FRI 1"	37,048	1,846	849	7,638	10,581	1,640	62	9,249	4,127	1,436
"02/12/06 HIGH SVEN OIL FRI 2"	40,353	1,680	1,281	3,119	18,291	848	181	3,654	4,348	815
John Engine Oil										
"02/12/06 HIGH JOHN OIL WED 1"	9,528	815	622	2,620	4,060	2,022	42	11,263	1,808	108
"02/12/06 HIGH JOHN OIL WED 2"	12,377	882	1,061	3,085	5,127	2,651	42	11,608	2,285	170
"02/12/06 HIGH JOHN OIL THUR 1"	20,628	773	579	2,119	8,451	1,284	38	8,007	1,647	229
"02/12/06 HIGH JOHN OIL THUR 2"	19,344	849	426	2,309	6,194	1,460	14	11,218	1,651	390
"02/12/06 HIGH JOHN OIL FRI 1"	19,299	901	803	2,304	4,623	1,336	26	7,147	1,230	45
"02/12/06 HIGH JOHN OIL FRI 2"	18,294	1,082	794	2,438	3,044	1,658	63	7,980	1,543	488
Ryan Engine Oil										
"02/12/06 HIGH RYAN OIL WED 1"	34,480	886	843	4,823	1,358	1,508	402	2,009	13,898	146
"02/12/06 HIGH RYAN OIL WED 2"	43,111	866	818	3,895	1,979	1,533	85	2,720	11,529	359
"02/12/06 HIGH RYAN OIL THUR 1"	30,252	1,227	1,747	3,921	1,451	1,111	113	3,517	8,861	288
"02/12/06 HIGH RYAN OIL THUR 2"	36,558	1,084	1,366	3,019	1,925	1,267	134	4,446	3,014	187
"02/12/06 HIGH RYAN OIL FRI 1"	25,781	1,548	1,311	5,203	706	1,684	168	2,530	9,700	185
"02/12/06 HIGH RYAN OIL FRI 2"	19,845	1,338	1,407	7,260	1,085	2,300	121	2,444	93,131	171
Dave Engine Oil										
"02/12/06 HIGH DAVE OIL WED 1"	39,283	812	638	1,850	4,538	965	43	2,832	2,083	186
"02/12/06 HIGH DAVE OIL WED 2"	38,511	919	470	1,824	3,953	625	35	3,268	1,951	195
"02/12/06 HIGH DAVE OIL THUR 1"	64,319	1,148	1,086	2,747	5,367	2,228	117	2,428	1,938	153
"02/12/06 HIGH DAVE OIL THUR 2"	42,659	1,295	1,188	2,398	4,470	827	53	1,658	703	61
"02/12/06 HIGH DAVE OIL FRI 1"	31,978	1,880	1,374	5,173	3,624	919	134	1,407	1,414	68
"02/12/06 HIGH DAVE OIL FRI 2"	32,451	1,907	1,308	4,846	3,628	873	132	1,451	1,481	71
Scott Engine Oil										
"02/12/06 HIGH SCOTT OIL WED 1"	31,370	1,565	1,158	3,942	8,175	2,938	125	4,118	11,947	76
"02/12/06 HIGH SCOTT OIL WED 2"	47,888	1,438	1,092	2,432	9,700	4,656	64	4,045	9,861	84
"02/12/06 HIGH SCOTT OIL THUR 1"	48,369	1,701	978	2,267	8,632	1,719	65	3,977	11,629	326
"02/12/06 HIGH SCOTT OIL THUR 2"	48,521	1,942	1,034	3,182	8,589	4,027	210	6,581	16,289	182
"02/12/06 HIGH SCOTT OIL FRI 1"	55,344	2,072	1,271	2,385	11,859	1,793	195	4,175	11,527	220
"02/12/06 HIGH SCOTT OIL FRI 2"	44,011	2,184	1,184	3,154	10,307	1,951	136	4,333	14,401	819
Average Engine Oil - John										
Average Engine Oil - John	18,578	883	714	2,481	5,250	1,735	37	9,339	1,694	240
Average Engine Oil - Scott										
Average Engine Oil - Scott	45,917	1,813	1,116	2,897	9,544	2,847	132	4,555	12,608	284

APPENDIX EXPERIMENT 15

Element - Raw Counts	Cs	Eu	Dy	Yb	Hf	Hg	Pb	U
T021206 HKH SVEN OIL WED 2"	93	5	3	9	84	182	43,127	105
T021206 HKH SVEN OIL THUR 1"	182	10	13	4	98	433	68,576	11
T021206 HKH SVEN OIL THUR 2"	638	20	14	12	225	359	65,027	128
T021206 HKH SVEN OIL FRI 1"	511	22	14	19	100	172	77,492	163
T021206 HKH SVEN OIL FRI 2"	830	24	12	15	172	191	59,027	157
John Engine Oil								
T021206 HKH JOHN OIL WED 1"	66	5	0	5	23	359	21,483	45
T021206 HKH JOHN OIL WED 2"	178	7	4	7	65	516	21,181	60
T021206 HKH JOHN OIL THUR 1"	81	2	2	7	100	371	11,586	78
T021206 HKH JOHN OIL THUR 2"	124	3	7	9	112	372	13,121	72
T021206 HKH JOHN OIL FRI 1"	57	3	3	4	15	419	12,803	63
T021206 HKH JOHN OIL FRI 2"	97	1	2	0	98	284	15,103	52
Ryan Engine Oil								
T021206 HKH RYAN OIL WED 1"	285	5	3	9	35	414	13,311	148
T021206 HKH RYAN OIL WED 2"	756	9	10	11	51	483	10,075	182
T021206 HKH RYAN OIL THUR 1"	231	6	7	5	139	706	15,113	111
T021206 HKH RYAN OIL THUR 2"	487	17	5	13	48	701	10,011	147
T021206 HKH RYAN OIL FRI 1"	390	13	7	32	19	403	9,644	107
T021206 HKH RYAN OIL FRI 2"	218	4	8	11	25	428	11,499	134
Dave Engine Oil								
T021206 HKH DAVE OIL WED 1"	180	4	6	7	84	134	34,687	152
T021206 HKH DAVE OIL WED 2"	111	2	4	18	53	143	41,494	137
T021206 HKH DAVE OIL THUR 1"	559	3	5	17	69	252	37,827	205
T021206 HKH DAVE OIL THUR 2"	61	58	5	9	24	279	35,291	136
T021206 HKH DAVE OIL FRI 1"	50	5	7	12	8	170	40,070	94
T021206 HKH DAVE OIL FRI 2"	44	5	7	11	9	149	43,876	99
Scott Engine Oil								
T021206 HKH SCOTT OIL WED 1"	246	6	9	18	172	314	7,919	158
T021206 HKH SCOTT OIL WED 2"	115	6	8	8	35	208	6,563	158
T021206 HKH SCOTT OIL THUR 1"	63	6	7	8	97	292	6,177	190
T021206 HKH SCOTT OIL THUR 2"	80	14	17	12	53	427	7,912	185
T021206 HKH SCOTT OIL FRI 1"	94	12	9	14	105	191	5,894	177
T021206 HKH SCOTT OIL FRI 2"	137	11	9	9	104	322	8,832	143
Average Engine Oil - John								
Average Engine Oil - John	102	4	3	5	69	387	15,888	62
Average Engine Oil - Scott								
Average Engine Oil - Scott	128	9	10	11	95	292	6,863	164




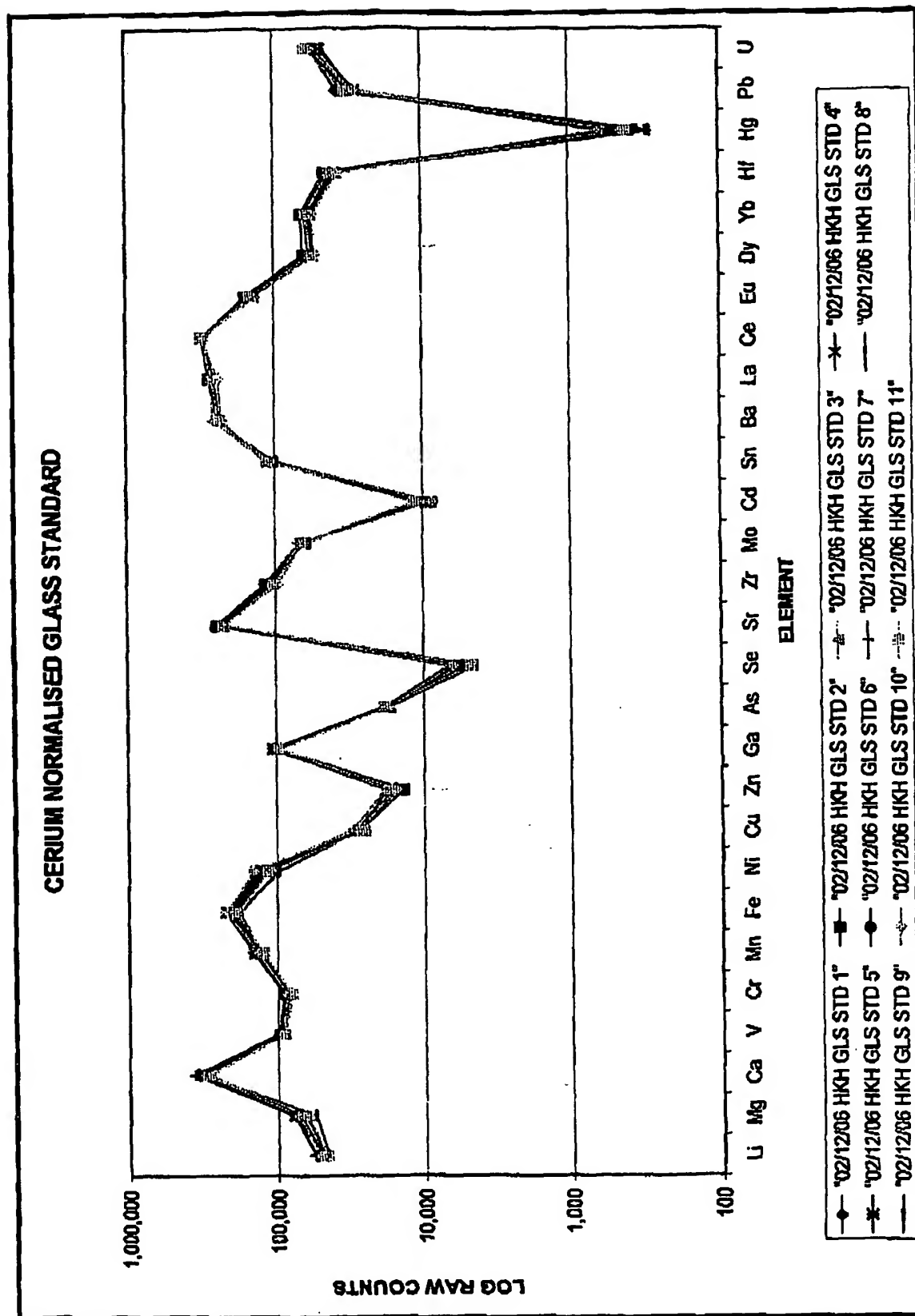


Chart Experiment15r1

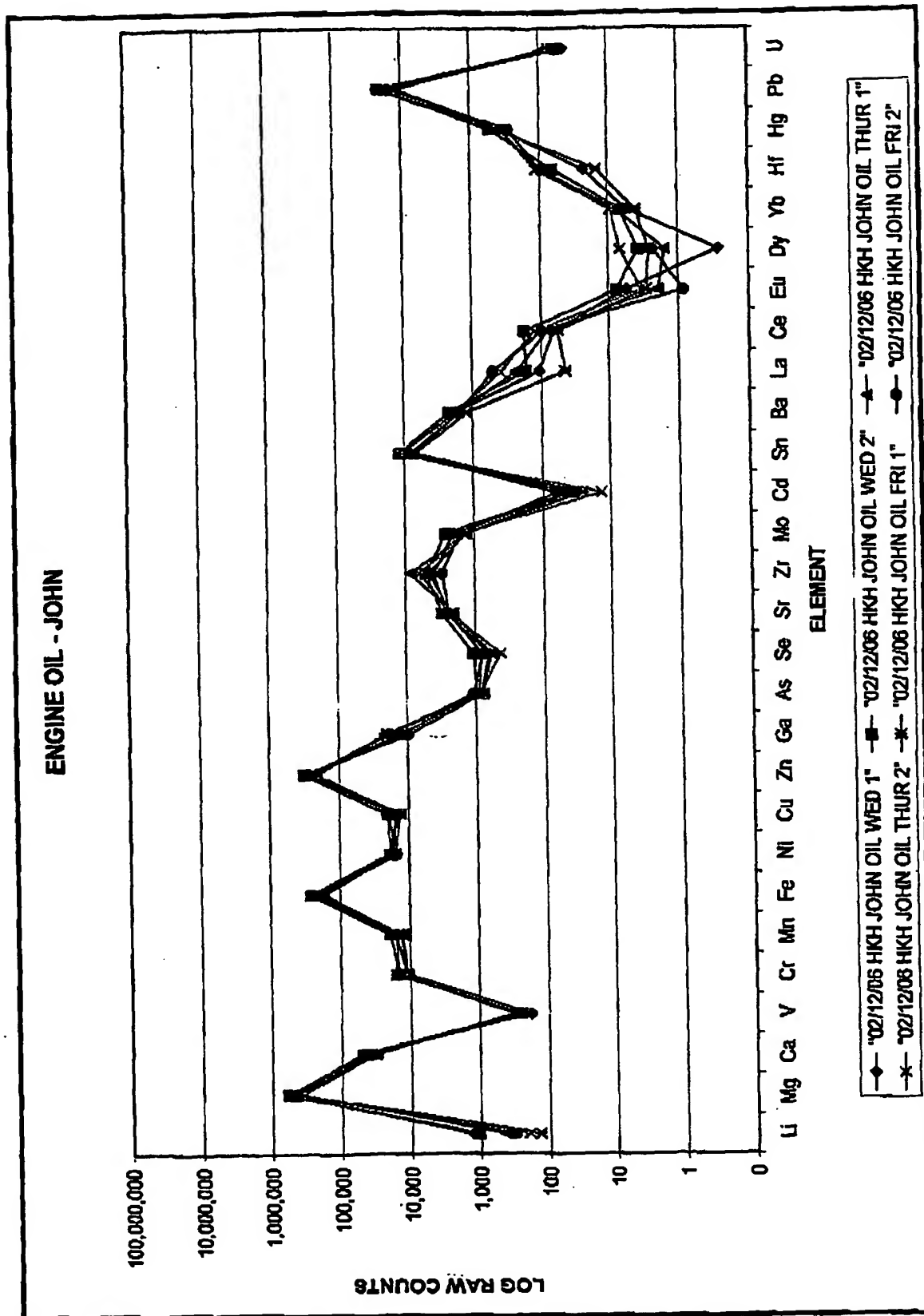

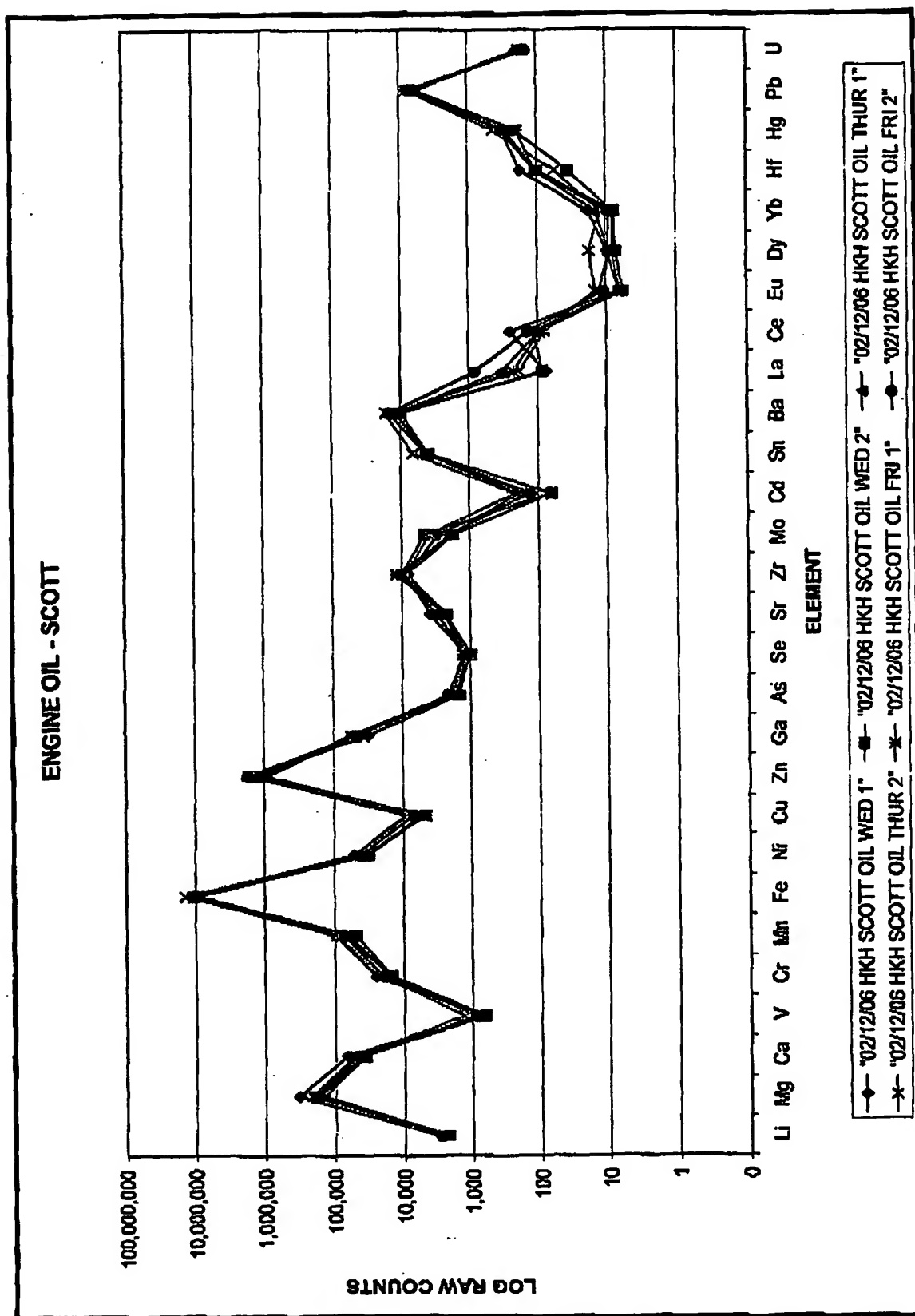

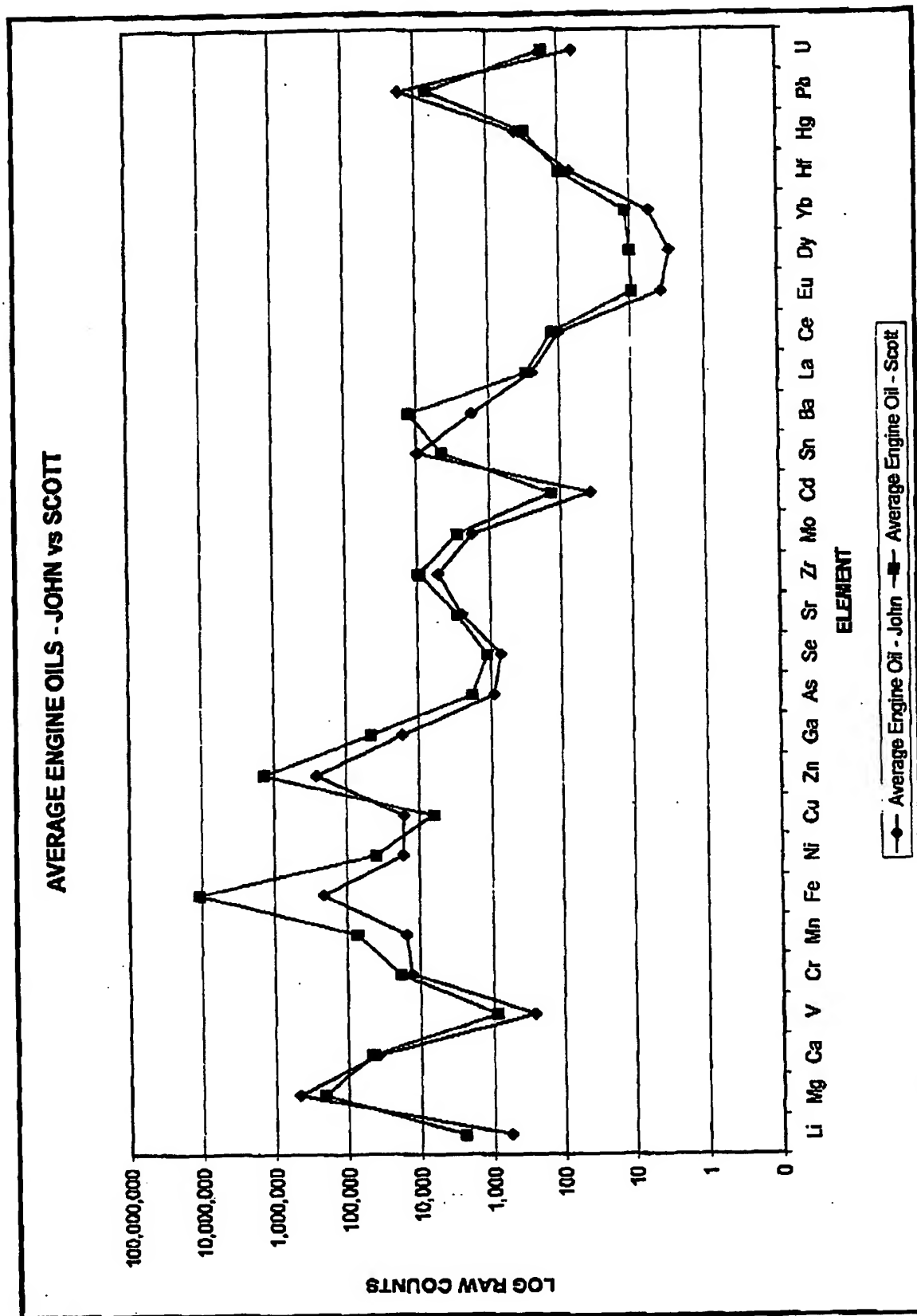

APPENDIX EXPERIMENT M1

Run	Normalized Data	7Li	9Be	51V	52Cr	55Mn	59Co	60Ni	65Cu	66Zn	68Ga	75As	82Se	85Rb	88Sr	89Y	90Zr
1	Blank TE 15/02/2003	8	0	182	281	42	25	111	23	18	20	18	4	20	21	1	33
2		8	1	184	281	41	24	112	23	18	20	17	4	19	21	1	28
3		8	1	150	283	42	24	110	24	19	20	17	4	18	21	1	26
4		8	0	140	288	42	24	112	24	18	20	18	5	18	22	1	24
5		8	0	132	286	42	24	110	23	19	20	17	4	17	21	1	23
	Mean	8.2	0.5	153.6	283.7	41.8	24.1	111.1	23.4	18.5	19.7	17.9	4.4	18.5	21.3	1.2	28.7
	Standard Deviation	0.2	0.0	19.9	3.0	0.4	0.3	0.8	0.5	0.3	0.2	1.0	0.1	1.3	0.2	0.0	3.7
	Coefficient of Variation	3.0	4.8	12.9	1.1	1.6	1.4	0.7	2.3	1.5	0.9	5.8	1.8	6.9	0.9	4.2	13.9
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
6	18-Feb-03	8	0	128	268	43	23	111	24	19	20	16	4	16	21	1	21
7		8	0	122	264	42	23	110	23	19	20	16	5	16	22	1	20
8		8	1	117	266	42	23	110	24	18	20	18	5	15	21	1	19
9		8	0	111	267	43	23	111	23	18	20	16	5	15	21	1	18
10		7	0	108	269	42	23	110	23	19	20	15	5	15	21	1	18
	Mean	7.8	0.4	118.7	268.5	42.2	23.1	110.8	23.4	18.4	20.1	15.5	4.6	15.6	21.3	1.1	18.1
	Standard Deviation	0.2	0.0	7.5	1.6	0.5	0.3	0.8	0.5	0.5	0.1	0.3	0.1	0.5	0.4	0.1	1.4
	Coefficient of Variation	2.8	10.4	6.4	0.8	1.3	1.4	0.7	2.0	2.4	0.6	2.3	3.2	3.3	1.7	6.7	7.1
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
1	0.1 ppm 15/02/2003	24	4	137	377	408	103	113	50	332	80	24	8	105	172	159	39
2		24	5	137	377	405	99	114	49	329	82	24	8	108	174	159	37
3		24	5	135	379	401	101	114	49	322	82	24	6	108	174	157	37
4		24	5	137	377	408	101	115	51	332	82	24	8	102	172	158	38
5		24	6	137	378	401	103	114	51	328	82	24	6	106	171	158	40
	Mean	23.7	4.8	138.4	377.7	403.9	101.4	114.2	49.8	328.7	81.8	23.8	5.8	104.8	172.8	157.8	37.8
	Standard Deviation	0.2	0.1	1.1	1.0	2.6	1.7	1.0	0.9	4.2	1.1	0.2	0.1	1.7	1.3	1.3	1.8
	Coefficient of Variation	0.7	2.7	0.8	0.3	0.7	1.7	0.9	1.8	1.3	1.2	1.0	1.9	1.7	0.7	0.8	4.3
	Count Limit 3 sigma	0.02	0.08	0.02	0.01	0.02	0.05	0.03	0.05	0.04	0.04	0.03	0.08	0.05	0.02	0.02	0.13
6	18-Feb-03	24	5	136	380	403	101	115	50	332	80	24	8	103	175	159	40
7		24	5	135	377	409	100	115	49	330	82	23	8	101	173	153	41
8		23	4	134	371	403	99	115	73	329	91	23	6	104	174	157	43
9		24	5	134	373	404	101	115	49	338	90	24	6	108	173	156	43
10		23	5	134	373	400	98	114	48	327	81	24	6	102	171	160	44
	Mean	23.7	4.5	134.7	374.8	403.7	100.1	114.3	53.7	328.7	80.9	23.5	5.8	103.1	173.2	158.4	42.2
	Standard Deviation	0.5	0.1	1.2	3.8	2.8	1.0	0.9	10.8	2.3	1.0	0.5	0.2	1.8	1.5	2.5	1.9
	Coefficient of Variation	2.0	1.4	0.9	0.9	0.7	1.0	0.8	20.1	0.7	1.0	2.8	3.4	1.8	0.8	1.6	4.2
	Count Limit 3 sigma	0.08	0.04	0.03	0.03	0.02	0.05	0.02	0.50	0.02	0.03	0.08	0.10	0.05	0.03	0.05	0.13
1	0.2 ppm 15/02/2003	38	9	211	444	585	178	131	69	268	164	34	7	183	382	397	80
2		38	8	211	432	553	173	130	68	263	183	33	7	195	387	312	81

APPENDIX EXPERIMENT M1

Run	Normalised Data	88Rb	90Mg	111Cd	120Sn	121Sb	128Te	138Ba	140Ce	141Pr	146Nd	153Eu	157Gd	158Tb	163Dy	168Tb
Blank TE 15/02/2003																
1		77	5	0	7	1	1	807	1	0	0	1	0	1	0	0
2		62	6	0	7	1	1	822	1	0	0	1	0	0	0	0
3		53	5	0	6	1	1	815	1	1	0	1	0	1	0	1
4		47	6	0	6	1	1	811	1	1	0	1	0	1	0	1
5		41	4	0	6	1	1	799	0	1	0	1	0	0	0	0
Mean		56.7	4.7	0.2	6.3	1.2	0.8	810.9	0.5	0.7	0.2	0.7	0.3	0.5	0.2	0.5
Standard Deviation		14.0	0.5	0.0	0.8	0.0	0.1	8.5	0.1	0.0	0.0	0.1	0.0	0.1	0.0	0.1
Coefficient of Variation		25.2	10.5	18.1	10.2	1.6	18.8	1.1	13.9	7.2	19.2	8.2	10.2	10.2	21.8	18.6
Count Limit 3 sigma		N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
16-Feb-03																
6		37	4	0	5	1	1	816	1	1	0	1	0	0	0	0
7		33	4	0	5	1	1	821	0	1	0	1	0	0	0	0
8		51	4	0	5	1	1	822	0	1	0	1	0	0	0	0
9		28	3	0	5	1	0	828	0	0	0	1	0	0	0	0
10		28	4	0	5	1	0	827	0	1	0	1	0	0	0	0
Mean		31.4	3.6	0.2	5.1	1.0	0.5	824.2	0.4	0.5	0.1	0.8	0.2	0.3	0.1	0.4
Standard Deviation		9.7	0.2	0.0	0.4	0.1	0.0	5.9	0.1	0.1	0.0	0.1	0.0	0.1	0.0	0.1
Coefficient of Variation		11.9	0.8	24.5	7.1	0.8	9.4	0.7	23.7	15.0	20.2	12.2	10.3	28.6	27.0	20.9
Count Limit 3 sigma		N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
0.1ppm 15/02/2003																
1		115	41	17	72	62	8	528	188	238	42	144	45	276	71	274
2		115	38	16	71	61	8	523	188	238	43	142	45	271	70	279
3		114	38	16	72	61	8	527	183	233	43	142	45	269	68	277
4		118	40	17	73	61	8	574	187	229	43	140	48	270	71	280
5		118	40	16	72	61	8	571	188	238	43	140	45	271	70	274
Mean		116.1	39.8	16.4	72.1	61.5	8.0	552.4	188.0	235.1	42.7	141.4	45.4	271.4	70.8	278.9
Standard Deviation		1.8	0.7	0.4	0.5	0.5	0.3	21.1	2.2	3.5	0.2	1.3	0.2	2.6	1.3	2.8
Coefficient of Variation		1.5	1.8	2.4	0.7	0.8	3.3	3.8	1.2	1.9	0.4	1.3	0.4	1.0	1.9	0.9
Count Limit 3 sigma		0.05	0.05	0.07	0.02	0.03	0.10	0.11	0.04	0.05	0.01	0.04	0.01	0.03	0.05	0.03
16-Feb-03																
6		118	39	16	72	61	8	578	188	237	44	144	45	267	70	275
7		112	40	17	73	61	8	584	188	233	44	140	44	269	70	272
8		114	40	16	72	60	8	573	185	237	42	143	44	268	67	275
9		114	38	16	72	60	8	571	187	231	43	141	46	268	68	278
10		114	48	16	73	62	8	586	184	230	42	142	44	268	68	268
Mean		114.3	38.4	16.2	72.8	60.8	7.8	574.0	188.0	233.7	42.9	141.8	44.8	268.2	68.9	273.5
Standard Deviation		2.8	0.4	0.4	0.7	1.1	0.1	6.5	1.5	3.2	1.0	1.8	0.5	0.7	1.2	3.7
Coefficient of Variation		1.8	1.0	2.3	0.9	1.9	1.8	1.1	0.8	1.4	2.3	1.2	1.1	0.3	1.7	1.3
Count Limit 3 sigma		0.05	0.03	0.07	0.03	0.06	0.05	0.03	0.02	0.04	0.07	0.03	0.03	0.01	0.05	0.04
0.2ppm 15/02/2003																
1		219	72	32	135	107	15	404	360	469	84	281	91	540	135	549
2		215	70	31	134	106	16	405	371	456	83	281	89	525	138	542

APPENDIX EXPERIMENT M1

Run	Normalized Data	168Ez	169Tm	172Tb	175Lu	178Hf	181Ta	182W	205Tl	208Pb	209Bi	232Th	238U
1	Blank TE 15/02/2003	0	1	0	1	49	13	49	3	19	10	33	1
2		0	1	0	1	41	11	43	3	19	8	25	1
3		0	1	0	1	38	9	40	2	19	7	21	1
4		0	1	0	1	34	10	40	2	21	8	18	1
5		0	1	0	1	29	8	34	2	19	5	16	1
	Mean	0.2	0.7	0.2	0.8	37.7	10.4	41.5	2.4	18.2	7.3	22.8	0.8
	Standard Deviation	0.0	0.1	0.0	0.1	7.7	1.6	5.4	0.4	0.7	2.0	8.8	0.1
	Coefficient of Variation	11.2	15.6	28.1	12.4	20.4	15.8	13.1	15.7	3.7	27.4	29.4	17.5
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
16-Feb-03													
6		0	1	0	1	27	9	32	2	18	5	14	1
7		0	1	0	1	25	7	30	2	19	4	13	1
8		0	1	0	0	23	9	30	1	18	4	12	1
9		0	0	0	0	21	7	28	2	19	4	11	1
10		0	0	0	0	21	7	27	2	18	3	11	0
	Mean	0.1	0.5	0.1	0.4	23.5	7.9	28.0	1.8	19.0	4.0	12.4	0.5
	Standard Deviation	0.0	0.1	0.0	0.1	2.7	1.0	2.3	0.1	0.4	0.5	1.3	0.0
	Coefficient of Variation	33.1	18.9	21.9	28.2	11.4	12.8	8.1	7.9	2.0	13.9	10.8	5.8
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
0.1ppm 15/02/2003													
1		94	291	88	289	43	232	33	186	620	198	209	232
2		93	288	66	294	45	232	35	188	625	188	208	221
3		94	288	67	297	47	236	33	188	622	192	214	239
4		94	290	65	290	48	236	32	194	622	188	212	239
5		92	288	68	298	46	233	33	180	621	193	214	235
	Mean	93.4	290.8	65.9	294.9	45.8	234.0	33.8	188.1	622.1	196.7	211.5	235.2
	Standard Deviation	0.9	4.2	0.9	2.9	1.9	2.7	1.2	2.1	1.7	3.1	2.9	3.5
	Coefficient of Variation	0.9	1.4	1.4	1.0	4.1	1.2	3.6	1.1	0.3	1.6	1.4	1.5
	Count Limit 3 sigma	0.03	0.04	0.04	0.03	0.12	0.03	0.11	0.03	0.01	0.05	0.04	0.04
18-Feb-03													
6		91	288	68	286	48	231	32	188	631	194	218	236
7		93	292	67	293	50	230	32	185	623	193	213	235
8		92	288	64	289	50	228	32	183	623	187	220	227
9		93	287	64	291	52	231	31	185	627	199	216	232
10		92	282	64	291	51	229	30	182	608	191	215	230
	Mean	92.2	288.9	66.1	292.3	50.4	229.7	31.4	184.2	622.3	194.7	218.5	232.3
	Standard Deviation	1.2	3.5	1.5	2.6	1.4	1.5	0.7	1.7	0.8	3.1	2.8	3.6
	Coefficient of Variation	1.3	1.2	2.3	0.9	2.8	0.7	2.1	0.9	1.4	1.6	1.3	1.6
	Count Limit 3 sigma	0.04	0.04	0.07	0.03	0.09	0.02	0.06	0.03	0.04	0.05	0.04	0.05
0.2ppm 15/02/2003													
1		180	571	129	560	100	431	60	370	946	394	422	486
2		187	568	128	581	102	433	64	397	946	393	426	459

APPENDIX EXPERIMENT M1

Run	Normalized Data	7Li	9Be	51V	52Cr	55Mn	58Co	60Ni	65Cu	68Zn	69Ga	75As	82Se	85Rb	86Sr	89Y	90Zr
3		38	8	212	438	561	173	130	89	205	183	33	7	195	281	305	84
4		38	8	209	437	547	177	132	70	208	161	34	7	190	285	310	85
5		39	8	208	420	555	176	130	72	203	162	33	7	194	283	305	85
	Mean	38.2	8.4	209.6	434.1	554.8	175.4	130.4	69.8	204.9	182.4	33.6	6.8	189.6	285.3	308.4	82.0
	Standard Deviation	0.3	0.4	2.2	8.1	6.7	2.2	0.9	1.5	1.3	1.3	0.4	0.3	2.2	3.6	2.6	2.4
	Coefficient of Variation	0.9	4.2	1.1	2.1	1.2	1.2	0.7	2.2	0.7	0.8	1.2	4.1	1.2	1.3	0.8	2.9
	Count Limit 3 sigma	0.03	0.13	0.03	0.08	0.04	0.04	0.02	0.07	0.02	0.02	0.04	0.12	0.03	0.04	0.03	0.03
18-Feb-03																	
6		38	8	209	405	552	171	130	69	203	181	33	7	193	285	308	88
7		39	8	206	415	548	178	130	67	201	182	33	7	192	277	312	81
8		39	8	209	410	550	174	132	68	198	182	33	7	191	279	302	92
9		38	8	208	404	558	170	130	69	207	165	33	7	188	282	305	85
10		37	8	207	409	559	176	132	68	203	163	33	7	188	278	303	85
	Mean	38.3	8.4	208.1	408.3	553.1	173.2	131.1	68.2	203.4	182.8	32.9	6.9	188.8	280.4	308.0	81.9
	Standard Deviation	0.6	0.3	0.9	4.4	4.3	2.6	1.1	0.8	4.0	1.5	0.3	0.1	2.4	3.3	3.8	3.7
	Coefficient of Variation	1.8	3.3	0.4	1.1	0.8	1.5	0.8	1.2	2.0	0.9	1.0	2.1	1.3	1.2	1.3	4.0
	Count Limit 3 sigma	0.05	0.10	0.01	0.03	0.02	0.04	0.02	0.03	0.08	0.03	0.03	0.08	0.04	0.03	0.04	0.12
1ppm 15/02/2003																	
1		154	39	768	873	1038	737	258	204	280	730	104	18	876	1240	1415	488
2		155	38	771	858	1029	739	254	203	281	738	103	18	887	1250	1430	509
3		155	39	757	850	1040	718	253	202	275	717	103	15	853	1226	1413	517
4		151	39	754	848	1039	732	253	201	277	732	104	18	872	1241	1418	531
5		154	39	759	867	1028	730	253	202	276	719	104	18	859	1243	1421	507
	Mean	153.8	38.7	759.5	859.4	1034.4	728.8	254.1	202.8	278.0	728.9	103.7	15.6	868.8	1240.5	1419.5	524.8
	Standard Deviation	1.7	0.6	8.5	10.9	8.3	8.0	2.4	1.2	2.4	8.4	0.7	0.2	8.3	8.8	6.5	34.5
	Coefficient of Variation	1.1	1.5	0.9	1.3	0.8	1.1	0.9	0.6	0.9	1.2	0.8	1.4	1.8	0.7	0.5	6.6
	Count Limit 3 sigma	0.03	0.04	0.03	0.04	0.02	0.03	0.03	0.02	0.03	0.03	0.02	0.04	0.03	0.02	0.01	0.20
18-Feb-03																	
6		156	39	768	856	1007	727	258	206	274	728	104	15	871	1248	1418	500
7		165	39	759	870	1016	725	254	200	282	721	104	16	872	1225	1397	583
8		150	38	763	866	1022	718	258	201	285	704	102	18	868	1231	1414	575
9		158	38	754	850	1025	720	250	203	278	722	104	18	877	1223	1400	577
10		155	38	782	863	1041	727	257	201	272	724	102	16	851	1258	1432	589
	Mean	154.2	38.4	759.4	851.6	1022.3	728.1	254.5	202.0	273.9	718.7	103.3	15.6	871.9	1237.4	1412.2	578.7
	Standard Deviation	2.5	0.8	4.0	8.4	12.8	4.8	2.9	2.4	8.1	8.1	1.2	0.4	8.5	15.1	14.2	14.1
	Coefficient of Variation	1.9	1.8	0.5	1.0	1.5	0.7	1.1	1.2	2.2	1.3	1.2	2.7	1.0	1.2	1.0	2.4
	Count Limit 3 sigma	0.05	0.05	0.02	0.03	0.04	0.02	0.03	0.04	0.07	0.04	0.04	0.08	0.03	0.04	0.03	0.07
15/02/2003																	
1		757	182	3569	5163	5035	3400	860	818	823	3965	462	69	4341	6338	7240	3083
2		745	188	3559	5124	5065	3388	861	901	815	3857	464	80	4408	6238	7247	3051
3		744	185	3558	5128	4988	3318	870	921	828	3872	469	58	4327	6239	7147	3097
4		758	189	3604	5134	4915	3484	853	907	817	3972	458	59	4318	6208	7174	3055
5		750	188	3568	5134	4955	3508	858	925	822	3817	463	58	4333	6383	7209	3141

APPENDIX EXPERIMENT M1

Run	Normalized Data	931b	934b	111Gd	120Sn	121Sb	126Te	138Ba	139La	140Ce	141Pr	148Nd	153Eu	157Gd	159Tb	163Dy	168Ho
3		214	72	32	134	109	18	413	364	358	456	85	279	89	531	135	548
4		212	71	32	135	108	15	408	365	358	450	83	278	89	532	136	541
5		214	68	32	135	108	16	404	367	358	453	83	278	89	532	136	547
	Mean	214.9	70.6	31.7	134.5	107.8	15.3	405.9	365.3	358.2	456.7	83.4	278.1	88.8	530.3	135.6	545.3
	Standard Deviation	2.5	1.3	0.7	1.0	1.0	0.4	3.7	4.0	2.2	7.3	0.8	1.9	0.8	1.2	0.8	3.8
	Coefficient of Variation	1.1	1.8	2.1	0.6	1.0	2.7	0.9	1.1	0.6	1.6	0.9	0.7	1.0	1.2	0.9	0.7
	Count Limit 3 sigma	0.03	0.05	0.08	0.02	0.03	0.08	0.03	0.03	0.02	0.06	0.03	0.02	0.03	0.04	0.03	0.02
18-Feb-03																	
6		212	71	31	133	109	15	409	367	352	455	85	274	88	522	136	542
7		214	68	31	136	108	14	404	358	358	458	85	271	87	518	136	527
8		217	68	31	134	107	15	410	364	358	449	82	276	88	522	136	538
9		212	69	31	134	106	15	421	369	353	457	84	273	88	531	137	537
10		212	70	31	135	107	18	424	365	353	458	83	278	88	516	137	531
	Mean	213.4	69.5	31.2	134.5	107.5	15.1	413.6	362.6	353.9	455.2	83.8	274.0	88.9	521.6	138.4	534.8
	Standard Deviation	2.0	0.7	0.3	1.2	1.0	0.5	8.3	4.1	3.3	3.7	1.5	2.4	1.0	8.1	0.5	6.0
	Coefficient of Variation	0.9	1.0	1.0	0.9	1.0	3.0	2.0	1.1	0.9	0.8	1.7	0.9	1.1	1.2	0.4	1.1
	Count Limit 3 sigma	0.03	0.03	0.03	0.03	0.03	0.08	0.08	0.03	0.03	0.02	0.05	0.03	0.03	0.04	0.01	0.03
19pm 15/07/2003																	
1		948	333	148	595	436	70	1387	1721	1681	2127	392	1329	408	2503	838	2580
2		982	327	147	604	443	68	1404	1721	1689	2147	390	1293	418	2520	842	2805
3		928	332	142	600	433	69	1385	1704	1630	2129	385	1307	413	2484	828	2578
4		967	325	148	607	440	70	1430	1682	1688	2171	388	1301	412	2474	840	2555
5		950	332	144	592	437	70	1350	1688	1629	2113	387	1288	411	2458	849	2573
	Mean	948.9	328.7	145.3	598.5	437.8	69.5	1403.2	1701.1	1657.5	2132.2	390.2	1300.9	412.4	2487.4	839.0	2580.0
	Standard Deviation	12.7	5.3	2.3	6.2	3.4	0.8	18.0	23.4	28.0	22.3	6.1	10.1	3.9	24.9	7.6	15.1
	Coefficient of Variation	1.3	1.6	1.6	1.0	0.8	1.2	1.1	1.4	1.7	1.0	1.3	0.8	1.0	1.0	1.2	0.9
	Count Limit 3 sigma	0.04	0.03	0.05	0.03	0.02	0.04	0.03	0.04	0.06	0.03	0.04	0.02	0.03	0.03	0.04	0.02
18-Feb-03																	
6		935	330	142	593	430	69	1421	1688	1651	2168	390	1305	410	2484	838	2587
7		951	326	142	588	439	69	1400	1681	1647	2124	388	1295	411	2468	837	2582
8		951	324	143	588	430	69	1398	1725	1670	2188	391	1312	411	2484	843	2551
9		955	328	147	600	439	68	1424	1684	1645	2147	397	1326	421	2489	843	2518
10		942	328	147	600	435	71	1417	1701	1688	2180	389	1328	414	2512	844	2557
	Mean	946.8	328.9	144.2	598.6	434.6	69.4	1410.2	1695.4	1658.2	2163.0	391.1	1312.3	413.7	2488.5	841.1	2553.0
	Standard Deviation	8.2	2.4	2.6	1.7	4.1	1.0	14.8	18.5	11.8	27.7	3.3	13.1	4.5	15.6	3.2	23.5
	Coefficient of Variation	0.9	0.7	1.8	0.3	0.9	1.4	1.1	1.1	0.7	1.3	0.8	1.0	1.1	0.6	0.5	0.9
	Count Limit 3 sigma	0.03	0.02	0.05	0.01	0.03	0.04	0.03	0.03	0.02	0.04	0.02	0.03	0.03	0.02	0.01	0.03
19pm 15/07/2003																	
1		4787	1561	712	3011	2175	344	7131	8768	8594	11080	1939	6821	2088	13083	3245	13341
2		4821	1589	721	3008	2185	344	7162	8754	8287	11135	1974	6897	2081	12888	3258	14281
3		4810	1591	710	2984	2223	338	7198	8687	8509	10918	1948	6727	2095	12895	3218	13474
4		4758	1580	708	2940	2143	338	7041	8685	8477	11165	1968	6745	2102	13112	3221	13424
5		4720	1577	710	2964	2182	332	7312	8884	8539	11088	1934	6888	2059	12828	3243	13383

APPENDIX EXPERIMENT M1

Run	Normalized Date	166Ez	169Tm	1727b	175Ju	178H#	181Tb	182W	206T1	208Pb	209Bi	232Th	238U
3		182	554	129	578	105	427	62	382	833	382	433	450
4		181	585	130	575	106	429	70	370	827	384	431	450
5		181	558	128	568	105	423	62	359	818	383	430	455
	Mean	182.3	568.7	129.9	578.1	105.5	428.7	63.6	365.4	834.3	387.2	428.0	455.9
	Standard Deviation	2.8	12.1	0.7	5.2	2.5	3.7	3.9	5.1	11.5	6.0	4.2	6.8
	Coefficient of Variation	1.5	2.1	0.5	0.9	2.4	0.9	6.2	1.4	1.4	1.5	1.0	1.5
	Count Limit 3 sigma	0.06	0.08	0.02	0.03	0.07	0.03	0.19	0.04	0.04	0.05	0.03	0.04
16-Feb-03													
6		183	568	127	561	106	428	69	364	822	388	424	454
7		179	580	126	570	115	425	61	359	816	387	432	457
8		178	561	128	567	112	424	61	368	824	382	430	458
9		180	584	128	570	113	440	63	368	820	378	428	454
10		177	583	130	572	117	432	62	385	841	390	431	444
	Mean	178.7	582.8	128.3	568.2	112.3	430.0	63.1	382.6	826.0	385.4	428.2	452.9
	Standard Deviation	2.9	2.7	1.5	4.1	3.7	50.5	3.1	8.2	9.3	5.1	3.1	5.3
	Coefficient of Variation	1.1	0.5	1.2	0.7	3.3	11.2	5.0	1.7	1.1	1.3	0.7	1.2
	Count Limit 3 sigma	0.03	0.01	0.04	0.02	0.10	0.34	0.15	0.06	0.03	0.04	0.02	0.03
19pm 15/02/2003													
1		853	2720	805	2735	611	2283	303	1736	1188	1808	2080	2210
2		853	2722	813	2742	619	2326	305	1744	1178	1830	2082	2207
3		850	2688	815	2725	648	2328	308	1838	1178	1816	2112	2145
4		850	2727	818	2758	658	2315	441	1886	1191	1821	2061	2184
5		868	2704	813	2714	674	2312	404	1716	1183	1784	2088	2189
	Mean	858.3	2712.5	811.2	2735.0	641.1	2312.9	351.8	1718.7	1183.3	1811.3	2082.7	2188.8
	Standard Deviation	6.3	15.8	3.7	16.8	25.9	18.2	88.0	25.0	5.3	17.8	21.8	25.8
	Coefficient of Variation	0.7	0.6	0.6	0.6	4.0	0.8	18.8	1.5	0.4	1.0	1.0	1.2
	Count Limit 3 sigma	0.02	0.02	0.02	0.02	0.12	0.02	0.58	0.04	0.01	0.03	0.03	0.04
16-Feb-03													
6		855	2699	811	2738	687	2284	308	1783	1208	1789	2083	2183
7		850	2681	807	2724	674	2287	305	1728	1174	1839	2076	2184
8		855	2725	807	2710	683	2271	300	1734	1172	1778	2088	2150
9		847	2677	808	2717	686	2300	345	1711	1189	1782	2078	2158
10		852	2684	802	2735	679	2283	304	1720	1178	1838	2074	2159
	Mean	854.8	2693.2	807.0	2734.8	677.7	2287.0	312.3	1731.1	1180.1	1804.8	2073.1	2168.0
	Standard Deviation	6.8	18.6	3.2	12.0	7.0	10.9	18.8	18.7	18.0	31.3	13.1	18.8
	Coefficient of Variation	0.8	0.7	0.5	0.4	1.0	0.5	6.0	1.1	1.4	1.7	0.6	0.9
	Count Limit 3 sigma	0.02	0.02	0.02	0.01	0.03	0.01	0.18	0.03	0.04	0.05	0.02	0.03
5pm 15/02/2003													
1		4352	14247	3083	14951	3580	11584	1572	8887	6821	8318	10830	11280
2		4338	14147	3060	14839	3608	11557	1504	8858	6806	8284	10826	11251
3		4378	14099	3148	14440	3723	11788	1608	8915	6820	8338	10805	11388
4		4327	14571	2884	14728	3689	11433	1559	8789	6882	8334	10776	11294
5		4319	14782	3125	15081	4051	11388	1573	8774	6886	8288	10889	10857

APPENDIX EXPERIMENT M1

Run	Normalized Data	7Li	8Be	51V	52Cr	56Fe	58Co	60Ni	64Cu	68Zn	69Ga	76As	82Se	83Rb	88Sr	89Y	90Zr
1	Mean	750.4	188.1	3571.8	3136.7	4853.1	3535.7	859.8	914.3	820.9	3556.7	463.0	53.9	4344.9	8277.4	7203.4	3087.8
2	Standard Deviation	5.9	2.5	18.6	15.3	80.3	58.8	8.3	10.1	4.9	23.1	4.7	0.8	38.5	68.4	42.8	37.2
3	Coefficient of Variation	0.8	1.3	0.5	0.5	1.2	1.6	0.7	1.1	0.6	0.6	1.0	1.4	0.8	1.1	0.8	1.2
4	Count Limit 3 sigma	0.02	0.04	0.02	0.01	0.04	0.05	0.02	0.03	0.02	0.02	0.03	0.04	0.03	0.03	0.02	0.04
5	18-Feb-03																
6	Mean	757	189	3557	3058	5028	3581	858	915	829	3553	464	59	4375	8383	7286	3170
7	Standard Deviation	754	191	3594	3181	5030	3597	858	927	828	3558	462	60	4388	8398	7276	3210
8	Coefficient of Variation	749	185	3605	3183	5003	3523	850	904	824	3553	457	58	4275	8272	7208	3182
9	Count Limit 3 sigma	752	191	3583	3167	4971	3481	845	903	811	3558	463	59	4302	8204	7266	3147
10	Mean	748	188	3580	3180	4908	3497	859	909	812	3522	460	59	4318	8189	7116	3077
11	Standard Deviation	731.4	188.4	3578.8	3153.9	4988.9	3535.8	851.2	916.4	820.7	3528.3	461.1	59.2	4320.7	8285.0	7231.1	3181.1
12	Coefficient of Variation	4.3	2.8	28.3	32.1	51.1	50.7	10.3	10.1	8.6	35.9	2.8	0.4	47.8	87.5	74.2	48.5
13	Count Limit 3 sigma	0.8	1.5	0.8	1.0	1.0	1.4	1.2	1.1	1.1	1.0	0.8	0.7	1.1	1.1	1.0	1.5
14	18-Feb-03																
15	Mean	1531	372	7229	6163	10888	7201	1804	1832	1332	7371	913	111	8804	12637	15704	8640
16	Standard Deviation	1524	374	7177	6120	11088	7218	1821	1845	1342	7259	914	109	8801	12644	15875	7263
17	Coefficient of Variation	1502	370	7257	6100	11047	7084	1810	1841	1332	7348	913	112	8831	12653	15740	8418
18	Count Limit 3 sigma	1514	355	7167	5891	10848	7092	1806	1888	1329	7299	888	109	8811	12698	15882	8500
19	Mean	1549	371	7202	5877	11031	7077	1822	1819	1332	7421	903	110	8829	12680	15757	8469
20	Standard Deviation	1524.1	370.4	7206.8	6070.3	11020.9	7130.3	1806.5	1841.3	1333.8	7321.7	903.3	110.2	8835.3	12371.1	15891.6	8537.7
21	Coefficient of Variation	17.9	3.3	37.0	82.0	63.8	73.0	10.5	10.6	8.1	85.7	7.3	1.5	78.9	134.9	133.5	352.2
22	Count Limit 3 sigma	1.2	0.9	0.5	1.4	0.5	1.8	0.7	1.0	0.4	1.2	0.8	1.3	0.9	1.0	0.9	5.3
23	18-Feb-03																
24	Mean	1488	378	7188	6051	11055	7054	1807	1821	1313	7401	881	109	8802	12780	15749	7038
25	Standard Deviation	1525	373	7245	5970	10973	7122	1802	1809	1318	7310	890	110	8748	12870	15888	7083
26	Coefficient of Variation	1493	375	7249	6108	11027	6988	1800	1784	1322	7310	889	109	8735	12888	15561	7102
27	Count Limit 3 sigma	1542	389	7187	6131	11074	7109	1813	1823	1301	7285	882	111	8729	12888	15537	8415
28	Mean	1535	389	7284	6139	10718	7150	1807	1842	1325	7343	888	110	8821	12787	15844	6381
29	Standard Deviation	1518.3	372.8	7222.5	6078.8	10889.4	7084.3	1807.8	1817.7	1315.8	7328.9	892.0	109.8	8788.5	12728.3	15823.7	6801.8
30	Coefficient of Variation	22.5	3.4	42.4	71.2	65.4	65.0	17.1	17.8	8.3	44.5	4.1	0.9	41.8	60.5	70.5	384.8
31	Count Limit 3 sigma	1.5	0.9	0.6	1.2	1.5	0.9	1.1	1.0	0.7	0.6	0.5	0.8	0.5	0.5	0.5	5.4
32	18-Feb-03																
33	Mean	878	141	2800	4180	6308	129	180	481	3033	11700	758	15	142867	5478	88252	104528
34	Standard Deviation	861	140	2800	4113	63217	137	182	885	3051	11732	758	14	140280	5503	85103	103077
35	Coefficient of Variation	873	140	2461	4125	61853	142	185	888	3007	11390	768	15	138428	5379	83325	103209
36	Count Limit 3 sigma	865	140	2413	4088	63183	147	189	877	2988	11482	763	15	140351	5388	82819	102867
37	Mean	887	139	2779	4176	61808	151	187	880	3031	11680	784	16	141545	5334	84342	101882
38	Standard Deviation	883.4	140.1	2530.7	4132.4	62658.9	141.2	183.5	880.3	3024.1	11590.8	788.1	14.8	140744.1	5404.0	84388.2	103047.8
39	Coefficient of Variation	6.1	0.7	172.3	35.6	687.2	8.6	2.8	11.8	21.0	157.5	9.8	0.5	1877.8	82.2	1978.4	981.0
40	Count Limit 3 sigma	0.7	0.5	6.8	0.9	1.1	6.1	1.4	1.3	0.7	1.4	1.5	3.2	1.2	1.5	1.5	1.0

APPENDIX EXPERIMENT M1

Run	Normalized Data	834No	834Mo	111Cd	120Sn	121Sb	126Te	138Ba	139La	140Ce	141Pr	148Nd	152Eu	157Gd	159Tb	163Dy	165Ho
1	Mean	4700.7	1575.6	710.8	2960.9	2183.7	337.4	7188.3	8757.7	8478.0	11079.4	1890.3	6886.8	2086.7	12861.0	3238.5	13580.8
2	Standard Deviation	42.1	14.8	7.7	32.0	28.0	7.0	98.1	79.2	114.5	88.2	21.6	47.0	18.9	128.0	18.4	394.8
3	Coefficient of Variation	0.9	0.9	1.1	1.1	1.3	2.1	1.4	0.9	1.4	0.8	1.1	0.7	0.8	1.0	0.5	2.9
4	Count Limit 3 sigma	0.03	0.03	0.03	0.03	0.04	0.06	0.04	0.03	0.04	0.03	0.03	0.02	0.02	0.03	0.02	0.09
5	18-Feb-03																
6	Mean	4812	1570	685	3033	2182	338	7247	8796	8432	11049	1878	6855	2082	13154	3274	13431
7	Standard Deviation	4785	1587	714	3085	2188	344	7178	8911	8622	11088	1888	6867	2108	12878	3281	13461
8	Coefficient of Variation	4789	1586	725	2978	2176	335	7109	8731	8628	11091	1888	6848	2059	13008	3241	13512
9	Count Limit 3 sigma	0.03	0.03	0.03	0.03	0.04	0.06	0.04	0.03	0.04	0.03	0.03	0.02	0.02	0.03	0.02	0.09
10	Mean	4784	1582	718	3033	2179	337	7284	8847	8463	11119	1889	6852	2087	13286	3278	13412
11	Standard Deviation	4788.5	1583.8	718.5	3028.6	2180.8	337.8	7184.3	8785.1	8522.1	11088.5	1878.1	6881.8	2088.4	13188.3	3272.6	13440.9
12	Coefficient of Variation	23.9	8.7	13.4	35.7	8.8	3.7	80.8	97.7	94.9	61.4	18.6	63.8	21.2	310.2	18.8	47.7
13	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
14	10ppm 15/02/2003																
15	Mean	8570	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
16	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
17	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
18	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
19	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
20	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
21	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
22	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
23	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
24	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
25	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
26	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
27	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
28	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
29	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
30	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
31	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
32	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
33	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
34	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
35	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
36	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
37	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
38	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
39	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
40	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
41	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
42	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
43	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
44	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
45	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
46	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
47	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
48	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
49	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
50	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
51	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
52	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
53	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
54	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
55	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
56	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
57	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
58	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
59	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
60	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
61	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
62	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
63	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
64	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
65	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
66	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
67	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
68	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412	7315	28382
69	Coefficient of Variation	9853	3182	1439	6042	4388	650	14633	19085	19091	25060	4492	13841	4223	28259	7279	28307
70	Count Limit 3 sigma	0.02	0.02	0.02	0.04	0.01	0.03	0.03	0.03	0.03	0.01	0.03	0.03	0.03	0.07	0.02	0.01
71	Mean	8571	3173	1385	6112	4431	680	15127	18559	19183	24335	4488	13748	4268	27807	7221	28835
72	Standard Deviation	9709	3219	1445	6180	4444	685	14820	19154	19185	24384	4483	13854	4350	28412</		

APPENDIX EXPERIMENT M1

Run	Normalized Date	166E1	166Tm	172Tb	175Lw	178H	181Ta	182W	205T1	208Pb	209Bi	227Th	238U
Mean		4354.9	14377.0	3061.9	14802.1	3730.6	11645.7	1583.0	8934.3	5963.4	9294.1	10821.1	11292.3
Standard Deviation		23.5	320.6	59.8	238.1	188.8	149.8	21.8	62.0	51.7	49.5	105.4	159.6
Coefficient of Variation		0.5	2.3	1.9	1.6	5.1	1.3	1.4	0.7	0.9	0.5	1.0	1.4
Count Limit 3 sigma		0.02	0.07	0.08	0.05	0.15	0.04	0.04	0.02	0.03	0.02	0.03	0.04
16-Feb-03													
6		4355	14151	3085	14746	3673	11621	1800	8808	5900	8651	10789	11474
7		4328	14252	3068	15103	4179	11809	1828	8825	6029	9278	10879	11438
8		4418	14630	3043	14809	3753	11501	1821	8872	5911	8253	10877	11210
9		4386	14754	3303	14888	3749	11563	1800	9055	6040	9339	11043	11518
10		4357	14859	3129	14938	3756	11803	1818	8857	5882	9297	10878	11260
Mean		4370.5	14549.3	3125.2	14971.3	3822.1	11817.4	1813.2	8903.1	5948.5	9343.8	10913.5	11540.0
Standard Deviation		36.2	340.3	104.5	173.8	202.8	117.0	12.9	88.8	80.7	120.1	98.3	113.4
Coefficient of Variation		0.8	2.3	3.3	1.2	5.3	1.0	0.8	1.0	1.4	1.3	0.9	1.0
Count Limit 3 sigma		0.02	0.07	0.16	0.04	0.16	0.03	0.02	0.03	0.04	0.04	0.03	0.03
10ppm 15/02/2003													
1		9721	30208	6885	30387	8454	24217	3889	18134	13381	18532	22131	23473
2		9446	29326	6825	30522	7735	24228	3816	18248	13755	19103	22222	23314
3		9634	29785	6863	30154	7610	24203	3744	18315	13578	18888	22578	23644
4		9520	29272	6722	30687	7694	24214	3768	18120	13931	19067	22488	23546
5		9429	28888	6781	30240	8370	24249	3784	18154	13910	19150	22714	23718
Mean		9549.8	29654.8	6777.2	30420.0	7972.8	24221.6	3762.2	18193.3	13568.9	18830.0	22428.0	23537.7
Standard Deviation		125.2	378.9	88.8	235.6	404.9	174	44.2	82.9	138.1	282.1	244.0	267.0
Coefficient of Variation		1.3	1.3	1.3	0.8	5.1	0.7	1.2	0.5	1.0	1.5	1.1	1.1
Count Limit 3 sigma		0.04	0.04	0.04	0.02	0.15	0.00	0.04	0.01	0.03	0.04	0.03	0.03
16-Feb-03													
6		9403	28895	6782	30331	7734	24003	3849	18281	13762	18840	22563	23825
7		9704	29891	6579	30388	8389	23770	3885	18282	13405	18774	22185	23980
8		9430	29078	6734	30151	8389	23802	3705	18093	13508	18591	22048	23306
9		9407	29084	6653	30041	8284	23889	3685	18295	13238	18571	22209	23385
10		9888	30071	6735	30611	8373	23909	3679	18498	13585	18581	22481	23784
Mean		9522.2	30043.4	6683.5	30284.3	8231.8	23894.5	3688.5	18283.8	13465.0	18887.4	22334.9	23388.0
Standard Deviation		148.6	48.4	67.7	188.0	281.4	103.1	27.2	132.8	183.9	118.7	208.6	247.0
Coefficient of Variation		1.6	0.2	1.0	0.6	3.4	0.4	0.7	0.7	1.4	0.6	0.9	1.1
Count Limit 3 sigma		0.05	0.00	0.03	0.02	0.10	0.01	0.02	0.02	0.04	0.02	0.03	0.03
SARM 1 15/02/2003													
1		8015	2888	4425	2813	5823	7200	570	747	22058	303	58265	21244
2		8028	2884	4425	2859	5821	7221	585	748	22048	279	58687	21419
3		5825	2827	4422	2844	6338	7288	554	757	21912	283	58824	21907
4		5885	2854	4434	2888	6229	7183	583	771	22272	251	58784	21844
5		5978	2814	4388	2859	6118	7267	582	754	21238	256	59188	21438
Mean		5874.1	28507	4420.9	2844.6	5833.7	7227.3	582.6	755.2	21894.7	270.7	58807.4	21450.8
Standard Deviation		61.3	32.1	13.7	21.4	208.5	48.8	8.0	8.6	431.6	21.8	393.3	244.4
Coefficient of Variation		0.9	1.1	0.3	0.8	3.6	0.7	1.1	1.3	2.0	8.1	0.8	1.1

APPENDIX EXPERIMENT M1

Run	Normalized Data	71U	886	51V	52C	56In	56Co	80Ni	86Cu	86Zn	68Ga	75As	82Se	85Rb	88Sr	88Y	90Zr
	Count Limit 3 sigma	0.02	0.01	0.20	0.03	0.03	0.16	0.04	0.04	0.02	0.04	0.04	0.09	0.04	0.05	0.04	0.03
18-Feb-03																	
6		871	159	2353	4131	82582	158	187	801	3072	11882	778	14	13202	5383	93272	102872
7		872	141	2335	4113	82005	158	184	880	3010	12153	763	14	138167	5420	93883	101857
8		872	142	2347	4171	83173	163	184	884	3043	11658	782	15	142107	5444	93807	103917
9		871	140	2338	4138	82500	187	163	885	3045	11855	778	15	141194	5436	94601	104928
10		883	144	2335	4307	82290	167	182	880	3043	11823	788	15	138891	5452	92323	102557
	Mean	871.0	141.2	2342.0	4171.9	82453.9	162.2	184.1	880.0	3042.7	11788.5	777.1	14.8	140110.2	5428.9	93887.5	102212.6
	Standard Deviation	1.8	1.8	6.0	78.2	435.1	5.1	1.7	8.4	21.9	222.8	9.1	0.4	1564.4	23.3	1355.9	1189.2
	Coefficient of Variation	0.2	1.3	0.3	1.9	0.7	3.1	0.9	0.9	0.7	1.8	1.2	2.5	1.1	0.4	1.4	1.2
	Count Limit 3 sigma	0.01	0.04	0.01	0.06	0.02	0.09	0.03	0.03	0.02	0.06	0.03	0.07	0.03	0.01	0.04	0.03
SARUM 3 14/02/2003																	
1		2716	459	27808	3459	2842853	788	290	980	21148	23316	331	8	81810	2808768	16978	3873433
2		2729	465	27590	3512	2818864	789	286	981	20959	22818	325	8	83744	2786570	16948	3854318
3		2764	484	28082	3532	2588870	813	284	1003	21463	22807	322	8	82043	2800893	16810	3908431
4		2778	470	28088	3530	2820828	815	285	1005	20859	23428	318	6	82151	2808053	16982	3933845
5		2781	472	27868	3547	2818720	820	283	1004	21430	23407	315	6	82479	2820968	17404	3857583
	Mean	2749.6	465.1	27827.5	3527.9	2817845.0	808.5	283.9	984.8	21191.8	23344.8	322.4	6.0	82445.4	2801208.2	17024.2	3905752.8
	Standard Deviation	25.8	5.1	203.6	25.3	18831.2	10.8	2.4	13.0	240.3	364.7	6.3	0.2	784.8	19183.8	233.7	42333.1
	Coefficient of Variation	0.9	1.1	0.7	0.7	0.7	1.3	0.8	1.3	1.1	1.6	2.0	2.8	0.9	0.7	1.3	1.1
	Count Limit 3 sigma	0.03	0.03	0.02	0.02	0.02	0.04	0.02	0.04	0.03	0.05	0.06	0.08	0.03	0.02	0.04	0.03
18-Feb-03																	
6		2769	483	28183	3329	2831153	801	287	986	21286	23286	311	6	82643	2827025	17244	3938447
7		2768	472	27850	3543	2834834	823	282	988	21540	23118	318	6	82747	2845483	17271	3924167
8		2787	473	28538	3483	2838253	820	280	1003	21548	23537	307	6	82821	2782828	16987	3988818
9		2827	477	28801	3390	2815682	825	282	989	21380	23304	308	6	82248	2770810	17094	3889781
10		2758	477	28733	3489	2823253	817	288	1011	21508	23388	302	6	83087	2827573	17330	3903391
	Mean	2781.3	473.5	28445.0	3525.2	2837811.1	817.0	281.7	988.2	21452.7	23424.4	307.4	6.2	82871.2	2812888.0	17180.9	3910478.8
	Standard Deviation	27.8	3.8	309.6	40.9	13078.5	8.7	3.1	10.9	113.7	280.2	3.4	0.3	304.1	30316.1	140.0	19247.6
	Coefficient of Variation	1.0	0.8	1.3	1.2	0.5	1.2	1.1	1.0	0.5	1.2	1.1	4.6	0.4	1.1	0.9	0.5
	Count Limit 3 sigma	0.03	0.02	0.04	0.03	0.01	0.04	0.03	0.03	0.02	0.04	0.03	0.14	0.01	0.03	0.03	0.01
SARUM 46 15/02/2003																	
1		986	17	61357	144481	4086008	21859	8216	44421	325432	5257	35385	13	9143	22031	8825	20942
2		995	17	61478	144517	4044171	21881	8636	43410	323332	5060	35180	13	9189	21484	8386	20448
3		977	16	60780	142017	4041842	21690	8892	42885	315842	5002	34887	13	9087	21670	9658	20202
4		981	16	60887	139845	4045884	21747	8883	42943	322296	4888	35142	12	8888	21305	9648	18880
5		1001	18	60818	141077	4077948	21425	8870	43093	325767	5064	35401	12	8891	21648	9734	19875
	Mean	987.9	16.4	61081.4	142287.5	4087880.6	21640.4	8375.4	43350.8	327533.8	5048.0	35197.4	12.4	9035.8	21527.3	8182.1	20271.4
	Standard Deviation	9.8	0.7	312.0	2288.2	22885.3	124.8	500.8	889.8	4010.4	140.7	208.9	0.2	147.3	288.0	630.2	444.0
	Coefficient of Variation	1.0	4.1	0.5	1.6	0.6	0.6	5.9	1.8	1.2	2.8	0.6	1.8	1.6	1.2	6.9	2.2
	Count Limit 3 sigma	0.03	0.12	0.03	0.05	0.02	0.02	0.18	0.05	0.04	0.08	0.02	0.05	0.05	0.04	0.21	0.07
18-Feb-03																	

APPENDIX EXPERIMENT M1

Run	Normalized Data	Count Limit 3 sigma	83Nb	89Mo	111Cd	120Sn	121Sb	126Te	138Ba	139La	140Ce	141Pr	148Nd	153Eu	157Gd	159Tb	163Dy	165Ho
6			28945	441	23	1279	185	1	80420	107747	103835	26870	15809	225	3489	3632	6128	5408
7			28753	442	24	1201	185	1	77820	104333	188028	28217	15807	230	3512	3837	6040	5421
8			30159	447	24	1212	185	1	78162	105605	188790	28083	16178	228	3502	3834	6135	5387
9			29900	438	24	1188	184	1	78333	105823	189171	28202	16258	224	3498	3833	6135	5474
10			30142	441	24	1201	186	1	78604	106357	182158	28082	16198	227	3510	3838	6132	5445
Mean			28882.2	441.9	23.9	1217.9	185.0	0.9	78547.9	105913.4	190400.0	28432.9	16107.7	227.2	3504.2	3888.4	6113.4	5432.8
Standard Deviation			335.3	3.4	0.6	34.4	0.9	0.0	1032.9	1255.6	2829.5	394.6	148.0	2.5	6.3	30.0	41.0	28.8
Coefficient of Variation			1.1	0.8	2.3	2.8	0.5	3.0	1.3	1.2	1.3	1.5	0.9	1.1	0.2	0.8	0.7	0.5
Count Limit 3 sigma			0.03	0.02	0.07	0.08	0.01	0.09	0.04	0.04	0.04	0.04	0.03	0.03	0.01	0.02	0.02	0.02
SARM 3 15/02/2019																		
1			358068	207	648	2139	27	10	274441	203870	258146	24857	8321	715	1348	723	904	827
2			374578	218	654	2158	27	8	274949	204313	257300	25180	10025	708	1356	725	918	832
3			387379	215	628	2155	28	10	271782	204782	258150	24856	10588	717	1370	728	907	830
4			395187	215	682	2242	28	9	271775	204231	255886	25280	11020	714	1342	721	917	827
5			398428	218	658	2210	28	9	271878	207002	258483	25284	10887	718	1374	738	925	838
Mean			384309.9	213.7	647.2	2180.9	27.6	9.4	272860.8	204418.7	257397.4	25277.3	10488.1	714.2	1357.7	728.8	913.8	830.5
Standard Deviation			16483.1	3.8	12.0	43.4	0.7	0.3	1594.3	1757.3	1481.7	284.0	874.8	3.2	14.4	7.3	8.4	3.7
Coefficient of Variation			4.1	1.7	1.9	2.0	2.5	3.8	0.8	0.9	0.6	1.1	8.4	0.4	1.0	1.0	0.9	0.4
Count Limit 3 sigma			0.12	0.06	0.08	0.06	0.06	0.11	0.02	0.03	0.02	0.03	0.19	0.01	0.03	0.03	0.03	0.01
16-Feb-03																		
6			378980	218	638	2158	28	9	277272	205041	254633	24708	10789	728	1385	736	938	831
7			370898	211	653	2155	28	8	275893	202689	255067	25113	10016	716	1378	752	918	828
8			398830	211	642	2158	27	9	269719	204582	255844	25050	10838	723	1385	738	923	840
9			384344	208	633	2138	28	9	274209	204408	253812	25104	10322	714	1377	744	931	827
10			382847	214	635	2157	27	9	271684	202943	281481	25316	10879	718	1361	740	937	838
Mean			387889.8	212.3	640.1	2152.5	27.4	8.8	273898.4	203808.8	258121.2	25078.3	10888.8	720.4	1378.8	741.8	928.8	832.5
Standard Deviation			5428.4	2.9	8.1	8.4	0.4	0.2	3851.3	1074.6	3875.4	185.3	80.1	6.0	12.7	6.3	8.2	5.8
Coefficient of Variation			1.5	1.4	1.3	0.4	1.8	1.7	1.1	0.5	1.2	0.7	0.7	0.8	0.9	0.9	0.8	0.7
Count Limit 3 sigma			0.04	0.04	0.04	0.04	0.05	0.05	0.03	0.02	0.04	0.02	0.02	0.03	0.03	0.03	0.03	0.02
SARM 48 15/02/2019																		
1			3517	117	3417	1841	253189	3	117294	15682	54588	4226	2881	458	676	584	781	575
2			3318	114	3421	1780	247832	3	118338	13838	54903	4188	2830	449	688	572	782	574
3			3245	113	3451	1805	250065	3	115431	13947	54782	4142	2815	435	680	578	788	583
4			3193	113	3366	1777	247112	2	118788	13775	54613	4149	2858	438	680	575	788	578
5			3005	114	3478	2342	249327	3	114801	13863	55134	4088	2885	450	673	581	778	583
Mean			3256.5	114.1	3422.3	1808.2	248911.1	2.9	116148.4	14188.9	54788.9	4159.8	2973.2	448.1	675.0	577.9	777.1	573.7
Standard Deviation			188.7	1.7	38.2	243.4	1378.5	0.1	677.9	855.1	223.8	53.8	47.3	9.2	5.8	4.8	5.3	4.3
Coefficient of Variation			5.7	1.4	1.1	12.7	0.6	4.9	0.8	6.0	0.4	1.3	1.6	2.1	0.8	0.8	0.7	0.7
Count Limit 3 sigma			0.17	0.04	0.03	0.38	0.02	0.15	0.03	0.18	0.01	0.04	0.06	0.06	0.02	0.02	0.02	0.02
16-Feb-03																		

APPENDIX EXPERIMENT M1

Run	Normalized Data	1685Z	1801m	1727b	1730u	1781H	1811b	1824W	2051I	2087b	2088b	237Th	238U
	Count Limit 3 sigma	0.03	0.03	0.01	0.02	0.12	0.02	0.03	0.04	0.06	0.24	0.02	0.03
15-Feb-03													
6	5938	2828	4412	2880	5231	7228	570	764	21730	242	58818	21493	
7	5932	2835	4448	2791	5222	7147	587	749	21549	253	57886	21183	
8	5865	2805	4384	2829	5207	7175	584	757	21348	323	58343	21422	
9	5835	2788	4334	2808	5178	7205	587	741	21804	390	58539	21247	
10	6059	2882	4344	2829	5148	7324	582	751	22280	373	58687	21702	
	Mean	5897.2	4380.4	2827.3	5187.2	7218.0	586.2	750.6	21743.9	318.2	58574.5	21411.6	
	Standard Deviation	50.0	27.9	48.3	34.6	33.7	3.0	6.1	382.9	87.8	887.5	282.6	
	Coefficient of Variation	0.8	1.0	1.1	1.2	0.6	0.9	0.5	0.8	1.8	21.4	1.2	1.0
	Count Limit 3 sigma	0.03	0.03	0.03	0.04	0.02	0.03	0.02	0.02	0.05	0.64	0.04	0.03
SARM 3 15022003													
1	993	489	821	538	98200	18377	1780	265	25387	714	68745	18839	
2	887	488	837	584	94146	18288	1875	260	26384	718	67780	18572	
3	1013	501	829	604	86120	19737	1861	262	25890	701	68357	18877	
4	1005	488	838	608	88344	19777	1879	254	25288	700	69001	18708	
5	1010	501	844	597	88356	18583	1871	256	25850	710	68851	18821	
	Mean	1001.7	488.8	833.3	588.5	18532.1	1813.2	258.0	25677.3	708.3	68788.8	18743.7	
	Standard Deviation	11.2	5.0	8.7	6.2	844.5	1389.1	86.2	3.0	381.2	7.5	703.2	137.5
	Coefficient of Variation	1.1	1.0	1.0	0.9	0.9	7.3	4.5	1.2	1.5	1.1	1.0	0.7
	Count Limit 3 sigma	0.03	0.03	0.03	0.03	0.03	0.22	0.13	0.04	0.05	0.03	0.03	0.02
16-Feb-03													
6	1015	514	840	604	88288	18183	1888	282	25883	712	68841	18713	
7	1009	485	823	605	98408	17404	1858	253	25717	778	68080	18888	
8	1021	500	837	608	98118	17819	1884	258	26133	855	68879	18880	
9	1015	502	825	588	94511	17880	1848	251	26384	839	68032	18887	
10	1009	501	838	606	97290	17358	1878	253	25878	794	68808	18883	
	Mean	1008.5	482.5	832.5	604.3	17850.3	1800.4	258.3	25878.8	785.3	68104.0	18884.2	
	Standard Deviation	6.5	7.0	7.9	4.1	1052.5	328.2	51.5	4.5	288.4	58.5	757.3	151.3
	Coefficient of Variation	0.8	1.4	0.9	0.7	1.1	1.9	2.7	1.8	1.1	7.1	1.1	0.8
	Count Limit 3 sigma	0.02	0.04	0.03	0.02	0.03	0.08	0.06	0.06	0.03	0.21	0.03	0.02
SARM 48 15022003													
1	538	217	319	200	614	407	550	220	7904988	9039	8488	1286	
2	544	218	314	203	681	412	578	222	7902718	8842	8494	1328	
3	528	216	313	200	585	408	575	220	8182828	9000	9173	1312	
4	540	218	307	201	588	408	582	221	8032868	9018	9407	1303	
5	535	218	315	208	582	408	570	222	8080575	8988	9428	1308	
	Mean	538.8	217.6	313.7	202.1	596.9	408.9	221.1	7984838.8	9023.3	8983.6	1308.8	
	Standard Deviation	6.0	1.1	4.4	2.4	12.4	4.4	4.7	1.0	84758.3	558	128.2	11.3
	Coefficient of Variation	1.1	0.5	1.4	1.2	2.1	1.1	0.8	0.5	1.1	0.6	1.4	0.9
	Count Limit 3 sigma	0.03	0.02	0.04	0.04	0.08	0.03	0.02	0.01	0.03	0.02	0.04	0.03
16-Feb-03													

APPENDIX EXPERIMENT M1

Run	Normalized Data	7Li	8Be	51V	52Cr	53Mn	56Co	60Ni	65Cu	66Zn	68Ga	75Se	82Se	86Sr	87Sr	90Zr
6		982	16	80435	139560	3994885	21239	8820	42848	318870	4953	35180	13	8782	21401	8931
7		1001	18	59838	138183	4015559	21423	10048	42441	321762	4886	34591	12	8830	21521	8973
8		991	16	60878	140588	4050870	21393	8825	42322	317054	4918	34337	12	8778	21245	8860
9		1000	17	60413	142639	4059888	21471	8918	42027	317286	4901	34489	12	8888	21201	8958
10		1008	18	60254	138839	4026499	21272	8730	42548	312219	4807	34334	12	8781	21577	8733
	Mean	988.5	16.3	60385.3	140079.5	4023520.4	21395.5	9028.9	42037.4	317574.2	4883.1	34844.1	12.0	8825.8	21388.9	8838.7
	Standard Deviation	6.8	0.2	3727	1807.7	25395.9	98.6	580.4	285.6	3556.4	54.3	378.9	0.4	81.5	188.2	82.4
	Coefficient of Variation	0.7	1.3	6.6	1.1	0.8	0.5	6.4	0.7	1.1	1.1	0.8	0.4	0.9	0.8	0.8
	Count Limit 3 sigma	0.02	0.04	0.02	0.03	0.02	0.01	0.19	0.02	0.03	0.03	0.03	0.08	0.03	0.02	0.03
Spinn check 15/02/2003																
1		794	200	3787	3098	4828	3519	888	902	871	3558	482	59	4341	6113	6885
2		808	200	3688	2883	4832	3455	883	825	880	3584	485	68	4331	6181	6860
3		824	201	3683	3130	4914	3901	888	928	851	3573	484	57	4412	6148	7009
4		808	202	3682	3067	4887	3491	886	937	858	3528	479	58	4347	6114	6882
5		802	199	3624	3080	4888	3392	884	946	842	3524	478	59	4293	6088	6879
	Mean	807.0	200.4	3688.1	3071.1	4880.9	3470.1	889.3	927.5	868.4	3548.8	484.2	57.5	4344.8	6130.8	6980.2
	Standard Deviation	11.2	1.1	50.8	54.7	58.4	49.0	8.9	7.9	10.7	22.6	5.9	1.0	43.0	33.7	18.5
	Coefficient of Variation	1.4	0.6	1.4	1.8	1.2	1.4	1.2	0.8	1.2	0.8	1.2	1.7	1.0	0.8	0.3
	Count Limit 3 sigma	0.04	0.02	0.04	0.05	0.04	0.04	0.03	0.03	0.04	0.02	0.04	0.05	0.03	0.02	0.02
18-Feb-03																
6		802	195	3621	3081	4822	3469	850	912	834	3543	472	57	4229	6087	6850
7		788	197	3576	2884	4839	3410	848	908	845	3430	486	57	4206	5983	6782
8		810	197	3593	3003	4783	3444	842	918	840	3488	485	58	4227	5987	6884
9		788	197	3554	2873	4794	3388	880	901	830	3535	488	56	4264	6021	6882
10		777	183	3544	2889	4755	3384	839	907	828	3448	458	55	4193	6023	6814
	Mean	793.0	195.7	3573.8	3010.0	4820.7	3417.1	845.4	908.7	834.9	3485.0	483.1	56.1	4228.5	6018.7	6830.3
	Standard Deviation	12.7	1.9	20.1	48.4	64.0	28.4	4.8	6.7	6.4	51.1	6.2	1.0	34.8	50.3	41.7
	Coefficient of Variation	1.9	1.0	0.8	1.5	1.3	0.8	0.8	0.8	0.8	1.5	1.3	1.8	0.8	0.8	0.6
	Count Limit 3 sigma	0.05	0.03	0.03	0.03	0.04	0.02	0.02	0.02	0.02	0.04	0.04	0.05	0.02	0.03	0.02
Blank TE 15/02/2003																
1		8	1	88	287	45	23	108	25	28	21	16	4	14	25	1
2		8	0	86	270	45	23	111	25	30	21	15	5	14	25	1
3		7	0	84	269	44	23	114	25	30	21	14	4	13	25	1
4		8	0	82	271	44	23	111	26	30	22	15	5	14	26	1
5		7	0	80	270	45	23	112	26	30	21	15	5	13	26	1
	Mean	7.5	0.4	84.1	289.4	44.4	22.7	111.6	25.4	28.4	21.3	14.9	4.5	13.6	25.6	1.0
	Standard Deviation	0.2	0.1	3.3	1.8	0.9	0.2	1.7	0.6	0.8	0.3	0.5	0.2	0.3	0.3	0.1
	Coefficient of Variation	2.4	13.1	3.5	0.6	1.3	1.0	1.5	2.4	2.8	1.4	3.3	4.4	2.2	1.1	8.0
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
18-Feb-03																
6		7	1	88	273	43	22	111	27	29	21	14	5	13	25	1
7		8	1	86	272	44	23	112	26	29	21	15	5	13	26	1
8		7	0	83	283	44	23	113	25	28	20	14	5	12	24	1

APPENDIX EXPERIMENT M1

Run	Normalized Data	83Nb	88Mo	111Cd	120Sn	121Sb	126Te	138Ba	139La	140Ce	141Pr	146Nd	153Eu	157Gd	159Tb	163Dy	165Ho
6		3222	111	3347	2378	241889	3	114377	13745	53228	4147	2688	433	677	568	771	571
7		3252	113	3377	1743	243262	3	112085	13369	54024	4082	2846	437	650	562	768	573
8		3131	112	3383	1739	243333	2	113482	13528	58153	4051	2818	436	669	571	767	579
9		3070	110	3365	1753	248370	2	114849	13933	54393	4087	2845	434	663	578	760	569
10		3003	110	3334	2344	246422	3	114205	13460	55172	4068	2859	442	673	570	770	567
	Mean	3135.4	111.1	3361.3	1981.4	244897.3	2.5	113807.6	13560.4	54081.8	4067.0	2851.1	436.8	668.2	568.3	770.7	571.8
	Standard Deviation	104.7	1.3	20.6	337.7	2631.7	0.2	1076.0	188.7	838.2	38.5	24.5	3.5	7.0	4.8	5.4	4.8
	Coefficient of Variation	3.3	1.2	0.6	17.0	1.1	6.2	0.9	1.4	1.5	0.9	0.9	0.8	1.0	0.6	0.7	0.8
	Count Limit 3 sigma	0.10	0.04	0.02	0.51	0.03	0.16	0.03	0.04	0.05	0.03	0.03	0.02	0.03	0.03	0.02	-0.02
50pm check 15/02/2003																	
1		4658	1523	680	2859	2125	318	8568	8489	8109	10089	1905	6340	1888	12422	3142	12778
2		5316	1528	682	2867	2117	322	8600	8398	8241	10532	1889	6409	2008	12430	3155	12874
3		4658	1527	681	2867	2107	327	8650	8464	8284	10585	1852	6344	1869	12354	3140	12910
4		4780	1501	649	2889	2109	322	8619	8457	8111	10582	1921	6408	2005	12687	3196	13012
5		4658	1515	688	2854	2078	324	8465	8278	8118	10593	1859	6400	2040	12742	3156	13207
	Mean	4813.0	1518.9	685.4	2875.3	2107.3	322.5	8582.4	8416.9	8172.4	10540.2	1882.8	6379.9	2071.8	12528.2	3148.8	12876.5
	Standard Deviation	284.5	11.2	6.4	21.9	18.1	3.3	73.6	85.7	83.9	97.4	25.7	34.7	28.3	177.3	8.1	182.8
	Coefficient of Variation	5.9	0.7	0.8	0.8	0.8	1.0	0.8	1.0	1.0	0.9	1.4	0.5	1.3	1.4	0.3	1.3
	Count Limit 3 sigma	0.18	0.02	0.02	0.02	0.03	0.03	0.03	0.03	0.03	0.03	0.04	0.02	0.04	0.04	0.01	0.04
16-Feb-03																	
6		5075	1525	675	2831	2091	319	8486	8215	8137	10478	1854	6384	1877	12500	3155	13122
7		5334	1488	681	2837	2097	324	8418	8184	8203	10539	1870	6384	1867	12353	3119	12980
8		4682	1503	688	2884	2059	322	8403	8284	8091	10583	1838	6318	1989	12340	3068	12745
9		4580	1485	678	2797	2054	321	8400	8344	8132	10284	1853	6386	1872	12261	3114	12943
10		4704	1481	672	2855	2065	319	8349	8290	8032	10450	1845	6351	1988	12405	3078	12707
	Mean	4935.1	1498.4	678.7	2830.9	2075.3	318.7	8410.8	8263.3	8088.9	10488.3	1852.0	6342.5	1878.2	12371.7	3102.8	12781.5
	Standard Deviation	385.9	17.1	6.5	20.8	18.3	4.4	48.2	63.8	73.2	78.0	11.6	38.4	3.3	88.3	28.0	222.4
	Coefficient of Variation	7.8	1.1	1.0	0.7	0.9	1.4	0.6	0.8	0.9	1.2	0.6	0.6	0.5	0.7	0.9	1.7
	Count Limit 3 sigma	0.23	0.03	0.03	0.02	0.03	0.04	0.02	0.02	0.03	0.04	0.02	0.02	0.01	0.02	0.03	0.05
Blank TE 15/02/2003																	
1		21	3	0	5	1	0	655	0	0	0	0	1	0	0	0	0
2		19	3	0	5	1	1	850	0	0	0	0	1	0	0	0	0
3		18	3	0	5	1	0	852	0	0	0	0	0	0	0	0	0
4		18	3	0	5	1	0	888	0	0	0	0	0	0	0	0	0
5		18	3	0	6	1	0	801	0	0	0	0	0	0	0	0	0
	Mean	18.6	3.1	0.2	5.3	0.6	0.4	870.7	0.3	0.4	0.1	0.1	0.5	0.2	0.1	0.0	0.1
	Standard Deviation	1.8	0.1	0.0	0.2	0.1	0.1	25.5	0.1	0.0	0.0	0.0	0.1	0.0	0.0	0.0	0.0
	Coefficient of Variation	9.8	2.8	19.7	4.5	10.9	12.7	2.9	28.3	9.8	23.3	27.2	10.9	24.8	19.3	28.1	31.3
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
16-Feb-03																	
6		16	3	0	6	1	0	877	0	0	0	0	0	0	0	0	0
7		15	3	0	6	1	0	875	0	0	0	0	0	0	0	0	0
8		15	3	0	6	1	0	880	0	0	0	0	1	0	0	0	0

APPENDIX EXPERIMENT M1

Run	Normalized Data	168E3	169Tm	172Tb	179Lu	178Hf	181Ta	182W	205Tl	208Pb	209Bi	232Th	238U
6		534	215	307	199	689	403	581	221	8147881	9229	9541	1324
7		534	215	314	203	687	405	581	219	8083497	9219	9498	1322
8		532	214	308	200	576	399	578	219	8128488	9219	9288	1303
9		532	220	308	202	588	391	582	219	8100315	9160	9464	1318
10		528	212	308	200	1082	392	1034	218	8104505	9182	9555	1319
	Mean	532.0	214.8	308.7	200.9	682.8	398.0	683.8	218.9	8079832.8	9201.7	9491.0	1318.7
	Standard Deviation	2.8	3.4	2.8	1.9	212.8	6.3	207.3	1.8	61330.2	28.2	138.8	8.2
	Coefficient of Variation	0.5	1.6	0.9	0.9	31.2	1.6	31.2	0.8	0.8	0.3	1.4	0.6
	Count Limit 3 sigma	0.01	0.05	0.03	0.03	0.94	0.05	0.94	0.02	0.02	0.01	0.04	0.02
Sppm check 15/02/2003													
1		4245	13801	2989	13831	3500	11085	1631	8585	5829	9157	10719	11380
2		4211	13529	3025	14728	3911	11377	1558	8844	5921	9170	10737	11238
3		4211	13714	2989	14429	3904	11505	1548	8827	5888	9387	10970	11348
4		4239	13788	3010	13988	3577	11478	1535	8744	5918	9168	10687	11415
5		4225	13783	2987	13940	3604	11229	1534	8731	5884	9208	10905	11284
	Mean	4226.4	13876.8	3010.1	14145.4	3705.3	11398.6	1545.1	8888.2	5883.8	9213.8	10798.8	11333.0
	Standard Deviation	15.7	176.4	17.1	410.8	282.2	173.3	11.5	66.8	38.6	88.1	130.8	71.7
	Coefficient of Variation	0.4	0.6	0.6	2.9	5.5	1.5	0.7	0.8	0.7	1.0	1.2	0.6
	Count Limit 3 sigma	0.01	0.02	0.02	0.08	0.18	0.05	0.02	0.02	0.02	0.03	0.04	0.02
16-Feb-03													
6		4196	13731	2980	14079	3585	11138	1542	8639	5893	9110	10813	11305
7		4118	13568	2952	13889	3666	11041	1538	8618	5749	9044	10640	11150
8		4154	13594	2958	13920	3689	11188	1541	8558	5807	9008	10531	11100
9		4125	13320	2988	13859	3538	11159	1504	8524	5732	8899	10504	11045
10		4173	13981	2993	13885	3583	10857	1488	8428	5798	8949	10428	10880
	Mean	4153.1	13392.8	2952.8	13808.9	3587.8	11098.9	1522.3	8633.3	5795.8	9071.8	10543.3	11088.0
	Standard Deviation	32.8	283.0	40.3	213.6	20.5	85.8	25.6	83.5	63.1	82.1	85.5	151.8
	Coefficient of Variation	0.8	2.1	1.4	1.5	0.9	0.9	1.7	1.0	1.1	0.9	0.8	1.4
	Count Limit 3 sigma	0.02	0.08	0.04	0.05	0.02	0.03	0.05	0.03	0.03	0.03	0.02	0.04
Blank TE 15/02/2003													
1		0	0	0	0	19	8	22	1	18	2	8	0
2		0	0	0	0	18	5	21	1	15	2	8	0
3		0	0	0	0	17	5	20	2	15	2	7	0
4		0	0	0	0	18	5	18	1	16	2	7	0
5		0	0	0	0	18	5	18	1	16	2	8	0
	Mean	0.0	0.2	0.0	0.2	18.9	5.3	19.9	1.4	15.8	2.1	7.2	0.3
	Standard Deviation	0.0	0.0	0.0	0.0	1.5	0.2	1.6	0.1	0.3	0.1	0.7	0.0
	Coefficient of Variation	81.1	12.9	20.8	22.9	8.8	4.2	7.9	7.5	2.0	7.0	8.9	14.9
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
16-Feb-03													
6		0	0	0	0	15	5	17	1	18	2	8	0
7		0	0	0	0	14	5	18	1	16	2	8	0
8		0	0	0	0	14	5	17	1	16	2	8	0

APPENDIX EXPERIMENT M1

Run	Normalized Data	7Li	9Be	51V	52Cr	53Mn	58Co	60Ni	62Zn	66Zn	75As	82Se	86Rb	88Sr	89Y	90Zr
9		8	0	82	271	44	28	110	28	28	14	5	13	25	1	12
10		8	0	80	273	44	23	113	28	28	14	4	13	25	1	13
	Mean	7.8	0.5	83.6	274.4	43.9	22.8	111.7	28.4	28.4	14.3	4.8	12.7	25.1	1.0	12.9
	Standard Deviation	0.2	0.0	3.4	4.7	0.4	0.3	1.2	0.8	0.7	0.5	0.2	0.3	0.4	0.1	0.4
	Coefficient of Variation	2.5	9.7	4.1	1.7	0.9	1.3	1.1	3.1	2.5	3.3	1.7	2.1	1.8	5.5	3.5
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
SARM 1 15/02/2003																
1		876	141	2800	4160	60018	129	180	3033	3033	768	15	142657	5476	86252	104528
2		861	140	2800	4113	63217	137	182	3051	3051	768	14	140280	5603	85108	103077
3		873	140	2881	4125	61858	142	185	3007	3007	768	15	138428	5378	85325	103207
4		885	140	2413	4088	63195	147	189	2988	2988	763	15	140351	5328	82819	102587
5		887	139	2378	4178	61938	151	187	3031	3031	784	15	141545	5334	84342	101882
	Mean	868.4	140.1	2530.7	4132.4	62658.9	141.2	188.5	3024.1	3024.1	768.1	14.8	140714.1	5404.0	84368.2	103047.5
	Standard Deviation	8.1	0.7	172.3	35.8	697.2	8.8	2.6	21.0	21.0	9.8	0.5	1077.8	82.2	1378.4	981.0
	Coefficient of Variation	0.7	0.5	8.8	0.9	1.1	6.1	1.4	0.7	0.7	1.3	3.2	1.2	1.5	1.5	1.0
	Count Limit 3 sigma	8.02	0.01	0.20	0.03	0.03	0.18	0.04	0.02	0.02	0.04	0.08	0.04	0.05	0.04	0.03
16-Feb-03																
6		871	139	2363	4131	62332	158	187	3072	3072	778	14	138202	5393	83272	102872
7		872	141	2336	4113	62005	168	184	3010	3010	763	14	139187	5420	83983	101857
8		872	142	2347	4171	63173	163	184	3043	3043	782	15	142707	5444	86907	103817
9		871	140	2339	4138	62500	167	183	3045	3045	778	15	141184	5438	84801	104928
10		868	144	2333	4307	62280	187	182	3043	3043	788	15	138881	5432	82328	102587
	Mean	871.0	141.2	2342.0	4171.8	62463.8	162.2	184.1	3042.7	3042.7	777.1	14.8	140110.2	5428.9	83887.5	103212.8
	Standard Deviation	1.8	1.8	8.0	78.2	435.1	5.1	1.7	21.9	21.9	9.1	0.4	1564.4	23.3	1563.9	1189.2
	Coefficient of Variation	0.2	1.3	0.3	1.9	0.7	3.1	0.9	0.7	0.7	1.2	2.5	1.1	0.4	1.4	1.2
	Count Limit 3 sigma	0.01	0.04	0.01	0.08	0.02	0.09	0.03	0.02	0.02	0.03	0.07	0.03	0.01	0.04	0.03
Average SARM 1																
		870	141	2438	4152	62561	152	188	3033	3033	773	15	140412	5418	84183	103130
	SARM Certified Value	12.80	7.75	2.00	12.00	154.68	0.36	8.00	90.00	90.00	18.80	0.01	325.00	10.00	143.00	300.00
	Counts per ppm	72	18	1218	348	404	421	23	61	61	40	1232	432	542	659	344
Concentrations in CRMs																
Based on SARM 1																
SARM 3 15/02/2003																
	Repeat	38	28	23	10	6481	2	13	349	349	8	<1	191	5172	26	11562
	Repeat	38	28	23	10	6331	2	13	354	354	8	<1	191	5193	26	11575
SARM 46 15/02/2003																
	Repeat	14	1	50	411	10071	51	403	5316	5316	12	<1	21	40	14	59
	Repeat	14	1	50	405	8861	51	388	5285	5285	11	<1	20	39	15	58
SARM 3 Cert Val																
		Li	Be	V	Cr	Mn	Co	Ni	Zn	Cu	As	Se	Rb	Sr	Y	Zr
	SARM 46 Cert Val	48.00	29.5	81	10	5903	2.44	2.20	305	13	1.92	0.01	190	4600	22	11020
				186	593		94	122	6200	583			18	28		95

APPENDIX EXPERIMENT M1

Run	Normalised Data	93Nb	89Mo	111Cd	123Sn	121Sb	126Te	138Ba	139La	140Ce	141Pr	146Nd	153Eu	157Gd	159Tb	163Dy	165Ho
9		15	3	0	5	1	0	881	0	0	0	0	1	0	0	0	0
10		15	3	0	6	1	0	882	0	0	0	0	0	0	0	0	0
	Mean	13.1	3.0	0.2	5.7	0.6	0.4	888.8	0.3	0.4	0.1	0.1	0.5	0.2	0.2	0.0	0.1
	Standard Deviation	0.4	0.1	0.0	0.2	0.0	0.1	8.3	0.0	0.0	0.0	0.0	0.1	0.0	0.0	0.0	0.0
	Coefficient of Variation	2.9	3.4	17.0	3.8	4.7	12.7	1.0	11.2	4.6	14.8	18.0	14.3	15.7	28.7	34.2	23.6
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
SARM 1 15/02/2003																	
1		30012	441	24	1219	186	1	80029	108943	194833	27029	16762	223	3500	3718	6097	5485
2		30183	458	24	1431	188	1	78824	108604	180505	26568	16142	231	3483	3718	6130	5442
3		28909	437	24	1204	188	1	78517	108331	189588	26808	16241	228	3483	3802	6025	5383
4		28565	445	23	1186	185	1	80247	108221	191387	26580	16372	228	3448	3883	6165	5800
5		28353	442	25	1183	184	1	79453	107173	182228	26406	16166	220	3494	3887	6134	5389
	Mean	28822.5	444.2	23.9	1245.1	185.7	0.9	79778.0	107134.4	191688.2	26344.2	16338.8	229.9	3477.6	3881.6	6110.3	5440.9
	Standard Deviation	347.1	7.3	0.5	104.5	1.4	8.1	888.6	1070.5	2093.2	224.4	253.9	2.1	21.8	47.6	53.3	57.8
	Coefficient of Variation	1.2	1.7	2.2	8.4	0.7	10.6	1.1	1.0	1.0	0.9	1.6	0.9	0.6	1.3	0.9	1.1
	Count Limit 3 sigma	0.03	0.05	0.07	0.25	0.02	0.32	0.03	0.03	0.03	0.03	0.05	0.03	0.02	0.04	0.03	0.03
16-Feb-03																	
6		28343	441	23	1278	185	1	80420	107747	181905	26670	15988	225	3499	3832	6128	5428
7		28753	442	24	1201	185	1	77920	104533	188028	26217	15807	230	3512	3887	6040	5421
8		30159	447	24	1212	185	1	78162	105505	188710	26803	16176	229	3502	3894	6165	5387
9		28908	438	24	1188	184	1	78633	105823	189711	26202	16258	224	3488	3883	6135	5474
10		30142	441	24	1201	185	1	78804	105337	182158	26382	16188	227	3510	3838	6182	5445
	Mean	28829.2	441.9	23.9	1217.9	185.0	0.9	78947.9	105913.4	190408.0	26432.9	16107.7	227.2	3504.2	3888.4	6113.4	5432.6
	Standard Deviation	335.3	9.4	0.8	34.4	0.9	0.0	1882.0	1256.6	2520.5	394.8	148.0	2.5	6.3	30.0	41.0	28.8
	Coefficient of Variation	1.1	0.8	2.3	2.8	0.5	3.0	1.3	1.2	1.3	1.5	0.9	1.1	0.2	0.8	0.7	0.5
	Count Limit 3 sigma	0.03	0.02	0.07	0.08	0.01	0.09	0.04	0.04	0.04	0.04	0.03	0.03	0.01	0.02	0.02	0.02
Average SARM 1																	
	SARM Certified Value	53.00	2.84	0.11	3.30	1.19	0.01	128.00	198.00	185.00	19.50	72.00	0.35	14.00	3.00	17.00	3.60
	Counts per ppm	583	156	211	373	158	129	881	977	980	1381	225	653	249	1225	380	1510
Concentrations in CRM's																	
Based on SARM 1																	
SARM 3 15/02/2003																	
	Repeat	675	1	3	6	<1	<1	413	208	253	18	47	1	5	1	3	1
	Repeat	653	1	3	6	<1	<1	414	208	261	18	48	1	6	1	3	1
SARM 48 15/02/2003																	
	Repeat	6	1	16	5	158B	<1	178	14	56	3	13	1	3	<1	2	<1
	Repeat	6	1	16	5	1571	<1	172	14	56	3	13	1	3	<1	2	<1
SARM 3 Cert Val																	
	SARM 48 Cert Val	9810	1.21	0.91	7.40	0.13	0.01	650	250	240	16	48	1.20	3.60	0.70	3.10	0.80
		28															

APPENDIX EXPERIMENT M1

Run	Normalized Data	165Er	169Tm	172Yb	175Lu	178Hf	181Ta	182W	205Tl	208Pb	208Bi	232Th	238U
8		0	0	0	0	14	5	16	1	15	2	8	0
10		0	0	0	0	14	5	17	1	15	2	6	0
	Mean	0.1	0.2	0.1	0.2	14.1	4.8	16.9	1.3	15.4	2.0	6.0	0.3
	Standard Deviation	0.0	0.0	0.0	0.0	0.8	0.2	0.6	0.1	0.3	0.1	0.2	0.0
	Coefficient of Variation	23.5	12.4	36.8	13.9	4.3	5.5	3.6	7.8	1.7	4.9	4.1	12.9
	Count Limit 3 sigma	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
SARM 1 15022003													
1		6015	2898	4425	2813	5825	7200	570	747	22056	305	58245	21264
2		6028	2884	4425	2858	5821	7221	565	748	22046	279	58957	21419
3		5925	2827	4422	2844	5828	7286	554	757	21512	283	58824	21307
4		5985	2854	4434	2868	5829	7183	543	771	22272	251	59704	21844
5		5916	2814	4398	2839	5188	7287	582	754	21238	258	59188	21439
	Mean	5974.1	2850.7	4420.9	2844.8	5383.7	7227.3	562.8	755.2	21824.7	270.7	59307.4	21450.8
	Standard Deviation	61.3	32.1	13.7	21.4	283.5	48.8	8.0	8.6	431.6	21.8	388.3	234.4
	Coefficient of Variation	0.9	1.1	0.3	0.8	3.9	0.7	1.1	1.5	2.0	8.1	0.8	1.1
	Count Limit 3 sigma	0.03	0.03	0.01	0.02	0.12	0.02	0.03	0.04	0.06	0.24	0.02	0.03
18-Feb-03													
6		5936	2829	4412	2880	5831	7228	570	764	21750	242	58618	21493
7		5982	2835	4448	2781	5722	7147	567	749	21549	253	57688	21193
8		5955	2805	4394	2829	5837	7175	564	757	21346	323	58343	21422
9		6035	2789	4334	2888	5178	7205	587	741	21804	390	58539	21247
10		6059	2882	4344	2829	5149	7234	662	751	22280	373	58887	21702
	Mean	5987.2	2824.1	4380.4	2827.3	5197.2	7246.0	588.2	750.8	21743.9	318.2	58574.5	21411.8
	Standard Deviation	50.0	27.9	48.3	33.6	33.7	67.7	3.0	6.1	362.9	87.8	697.5	203.6
	Coefficient of Variation	0.8	1.0	1.1	1.2	0.6	0.9	0.5	0.8	1.6	21.4	1.2	1.0
	Count Limit 3 sigma	0.03	0.03	0.03	0.04	0.02	0.03	0.02	0.02	0.05	0.94	0.04	0.03
Average SARM 1													
	SARM1 Certified Value	8888	2837	4401	2806	5280	7222	594	753	21784	293	58081	21431
	Counts per ppm	10.50	2.00	14.20	2.00	12.40	4.90	1.45	0.83	48.80	0.26	51.00	15.00
		570	1419	310	1418	428	1474	389	810	545	1067	1189	1428
Concentrations in CRU's													
Based on SARM 1													
SARM 3 15022003													
2		<1	3	3	<1	224	13	5	<1	47	1	59	14
Repeat		2	<1	3	<1	225	12	5	<1	48	1	60	14
SARM 48 15022003													
1		<1	1	1	<1	1	<1	1	<1	14880	8	8	1
Repeat		1	<1	1	<1	2	<1	2	<1	14834	9	8	1
SARM 3 Cert Yd													
	Er	Tm	Yb	Lu	Hf	Ta	W	Pb	Bi	Th	U		
SARM 48 Cert Yd	2.80	3.00	0.40	291.00	25.20	0.28	0.33	43	0.47	68	14		

APPENDIX EXPERIMENT M1

Run	Normalized Data	7U	9Ba	51V	52Cr	55Mn	56Co	60Ni	65Cu	66Zn	68Ga	75As	82Se	85Rb	88Sr	89Y	90Zr
	Samples diluted 250x prior to analysis																
	Calculated																
	Detection Limit Data																
	Based on standards																
	Concs in ppb	7U	9Ba	51V	52Cr	55Mn	56Co	60Ni	65Cu	66Zn	68Ga	75As	82Se	85Rb	88Sr	89Y	90Zr
		10	24	13	28	6	9	30	15	8	6	13	55	8	6	6	218

APPENDIX EXPERIMENT M1

Run	Normalized Data	92Nb	98Nb	111Cd	120Sn	121Sb	128Te	138Ba	140Ce	141Pr	148Nd	153Eu	157Gd	159Tb	163Dy	165Ho
	Samples diluted 250x prior to analysis															
	Calculated															
	Detection Limit Data Based on standards.															
		92Nb	98Nb	111Cd	120Sn	121Sb	128Te	138Ba	140Ce	141Pr	148Nd	153Eu	157Gd	159Tb	163Dy	165Ho
	Conc in ppb	8	9	9	6	6	18	9	5	7	8	5	6	7	4	5

APPENDIX EXPERIMENT M1

Run	Normalized Data	168Er	169Tm	172Yb	176Lu	178Hf	181Ta	182W	205Tl	208Pb	208Bi	232Th	238U
	Samples diluted 250x prior to analysis												
	Calibrated												
	Detection Limit Data												
	Based on standards:-												
	Concs in ppb	168Er 8	169Tm 8	172Yb 5	176Lu 5	178Hf 24	181Ta 38	182W 75	205Tl 9	208Pb 8	208Bi 9	232Th 5	238U 8

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.